



OCT 12, 2023

OPEN ACCESS



DOI:

dx.doi.org/10.17504/protocols.io.kqdg3xx5zg25/v1

Protocol Citation: Sabina Marciano, Roberta Marongiu 2023. Fecal Water Content Assay. **protocols.io**

<https://dx.doi.org/10.17504/protocols.io.kqdg3xx5zg25/v1>

MANUSCRIPT CITATION:

<https://pubs.rsc.org/en/content/articlelanding/2017/FO/C6F001559J>

License: This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](#), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

Protocol status: Working
We use this protocol and it's working


Fecal Water Content Assay

Sabina

Marciano¹, Roberta Marongiu¹

¹Weill Cornell Medicine

ASAP Collaborative Research Network

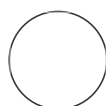

Dylan J Barbera

Weill Cornell Medicine

ABSTRACT

This test is used to check for signs of constipation

Created: Aug 14, 2023

Last Modified: Oct 12, 2023

PROTOCOL integer ID:
86468

Keywords: ASAPCRN

- 1 Label, weight and record the weight of each Eppendorf tube
- 2 Collect 2 pellets from each mouse and weigh.
- 3 Dry the pellets in the oven at 60 degrees C for 24 hrs
- 4 Weigh the pellets again after drying
- 5 Calculate fecal water content according to the equation:

(fecal weight before drying – fecal weight after drying)/feces weight before drying × 100.