



NOV 23, 2023

OPEN ACCESS



DOI:
dx.doi.org/10.17504/protocols.io.6qpvr3222vmk/v1

Protocol Citation: Marina Lorente Picón, Núria Peñuelas, Ariadna Laguna, Miquel Vila 2023. Tail Suspension Test.
protocols.io
<https://dx.doi.org/10.17504/protocols.io.6qpvr3222vmk/v1>

License: This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](#), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

Protocol status: Working
We use this protocol and it's working

Created: Nov 23, 2023

🌐 Tail Suspension Test


Marina Lorente
Picón¹,
Miquel Vila¹

Núria
Peñuelas¹,

Ariadna Laguna¹,

¹Vall d'Hebron Research Institute

Vilalab Public

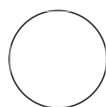

Miquel Vila

ABSTRACT

Tail Suspension Test for mice

- 1** Suspend animals by their tails with tape in a suspension bar.
- 2** To avoid the tail climbing behavior, pass a 2cm methacrylate tube through the tail before suspending the animal.
- 3** Quantify the escape-oriented behaviors (i.e. fore and hind limbs movement) during six minutes.
- 4** Calculate the total immobilization time (s) as the sum of all the time the animal was not performing any escape-oriented behavior.