




DEC 19, 2023

Motor behavioral evaluation of mice

Tae-Un Han¹

¹National Institute of Health

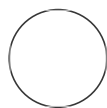

Tae-Un Han

ABSTRACT

This is a general protocol to test motor coordination of mice used models of neurological disorder

OPEN  ACCESS



Protocol Citation: Tae-Un Han 2023. Motor behavioral evaluation of mice.
protocols.io
<https://protocols.io/view/motor-behavioral-evaluation-of-mice-c6iwzcf>

License: This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](#), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

Protocol status: Working
We use this protocol and it's working

Created: Dec 18, 2023

Last Modified: Dec 19, 2023

PROTOCOL integer ID:
92470

Keywords: ASAPCRN,
Rotarod, Beam walk, Motor
behavior

Funders
Acknowledgement:

ASAP
Grant ID: ASAP-000458

Rotarod test

- 1 On first day, mice are acclimated to standing on a non-rotating rod (Rotamex 5, Columbus Instrument), 11m for three 00:01:00 trials with 00:10:00 interval.
- 2 Three hours later, mice are acclimated to walking on the rod over three 00:01:30 trials at a speed of 4 rpm, 1m 30s
- 3 On next day, latency of mice staying in the rotarod is recorded with accelerating speed of 4-40 rpm increasing by 1 rpm every 8.3 seconds for a maximum of 00:05:00 per trial. 5m
- 4 Repeat experiment for three trials daily for two days

Beam walk test

- 5 Mice are trained to cross beams of 1~1.5m length of round beam. 24mm and 18mm diameter of beams are used. 5m
- 6 Beam-crossing training is done for three times for two days
- 7 During the experimental day, its motor coordination will be measured by determining the number of slips, the latency required to cross the beam as the mouse traverses the beam.

- 8 For actual experiment, mice crosses both beams twice per trial for three trials with 5m rest periods in between trials