



JAN 11, 2024

OPEN  ACCESS



Protocol Citation: Anastasi
a Filimontseva, YuHong Fu,
Glenda Halliday 2024.
Neuromelanin Processing and
Image Acquisition.

protocols.io
<https://protocols.io/view/neuromelanin-processing-and-image-acquisition-c7eszjee>

License: This is an open
access protocol distributed
under the terms of
the [Creative Commons
Attribution License](#), which
permits unrestricted use,
distribution, and reproduction
in any medium, provided the
original author and source
are credited

Protocol status: Working
We use this protocol and it's
working

Created: Jan 11, 2024

Last Modified: Jan 11, 2024

PROTOCOL integer ID:
93362


Neuromelanin Processing and Image Acquisition

YuHong

Anastasia Filimontseva^{1,2}, Fu^{1,2}, Glenda Halliday^{1,2}

¹Brain and Mind Centre & Faculty of Medicine and Health School of
Medical Sciences, The University of Sydney, Sydney, NSW 2050, Australia;

²Aligning Science Across Parkinson's (ASAP) Collaborative Research
Network, Chevy Chase, MD, 20815

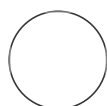

Anastasia Filimontseva

ABSTRACT

Protocol for preparing post-mortem tissue for intracellular NM quantification.

Tissue Processing

- 1 Formalin-fixed paraffin-embedded tissue blocks from the pontine and midbrain regions were sectioned at 6µm, collected onto adhesive microscope slides and allowed to dry in the oven at 37°C for 48 hours.
- 2 Slides were incubated in the oven at 60°C for 1 hour to melt the paraffin.
- 3 To remove the paraffin, sections were submerged in HistoChoice Clearing Agent (Sigma-Aldrich, H2779) for 2 x 7 minutes, followed by rehydration in decreasing ethanol concentrations (100% ethanol for 2 x 3 minutes, 95% ethanol for 3 minutes, 70% ethanol for 3 minutes) and distilled H2O for 3 minutes.
- 4 Slides were submerged in 1xPBS and coverslipped using Fluorescence Mounting Medium (DAKO, S302380-2).

Image Acquisition

- 5 Sections were scanned using 20x objective (NA=0.8) with pre-set focusing and exposure parameters for optimal NM signal quality with an automated Slide Scanner (SLIDEVIEW VS200, Tokyo, Japan).
- 6 Identical parameters were applied to each scanned section to ensure consistency in capturing NM.
- 7 Once the scan was completed, sections were submerged in 1xPBS to remove coverslips and continue with IF processing.