




May 13, 2021

Ultra-Rapid Sequencing (LAMP)

Jack Wadden¹¹University of Michigan - Ann Arbor

1 Works for me

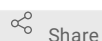

Share

This protocol is published without a DOI.

Jack Wadden
University of Michigan - Ann Arbor

ABSTRACT

This protocol accompanies the paper "Ultra-Rapid Somatic Variant Detection via Real-Time Threshold Sequencing." This protocol was followed to initiate a sequencing run that resulted in a somatic variant call from known tumor tissue in less than 30 minutes. The protocol outlines DNA extraction, LAMP, library preparation for Oxford Nanopore Sequencing, and sequencer preparation and loading.

PROTOCOL CITATION

Jack Wadden 2021. Ultra-Rapid Sequencing (LAMP). [protocols.io](https://protocols.io/view/ultra-rapid-sequencing-lamp-btvmnn46)
<https://protocols.io/view/ultra-rapid-sequencing-lamp-btvmnn46>

LICENSE

This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](#), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

CREATED

Apr 01, 2021

LAST MODIFIED

May 13, 2021

PROTOCOL INTEGER ID

48781

DISCLAIMER:

DISCLAIMER – FOR INFORMATIONAL PURPOSES ONLY; USE AT YOUR OWN RISK

The protocol content here is for informational purposes only and does not constitute legal, medical, clinical, or safety advice, or otherwise; content added to protocols.io is not peer reviewed and may not have undergone a formal approval of any kind. Information presented in this protocol should not substitute for independent professional judgment, advice, diagnosis, or treatment. Any action you take or refrain from taking using or relying upon the information presented here is strictly at your own risk. You agree that neither the Company nor any of the authors, contributors, administrators, or anyone else associated with protocols.io, can be held responsible for your use of the information contained in or linked to this protocol or any of our Sites/Apps and Services.

Preparation

- 1 Prepare the LAMP reaction tube by mixing the following components in a **0.2 mL** PCR tube

WarmStart LAMP Kit (DNA and RNA) - 500 rxns **New England**

- **12.5 µl** **Biolabs Catalog #E1700L**
- **2.5 µl** **[M]10 Micromolar (µM)** LAMP primer mix
- **9 µl** **Nuclease-free Water** **Contributed by users**

2 Prepare the DNA extraction tube by mixing the following components in a **0.2 mL** PCR tube

 **Epicentre QuickExtract™ DNA Extraction**

-  **100 µl** of **Solution Epicentre Catalog #QE09050**

3 Prepare ONT sequencing buffer

 **SQB (Sequencing Buffer) Oxford Nanopore**

-  **34 µl** **Technologies Catalog #SQK-RAD004**

SQB

 **LB (Loading Beads) Oxford Nanopore**

-  **25.5 µl** **Technologies Catalog #SQK-RAD004**

-  **4.5 µl** **nuclease free water Contributed by users**


4 Prepare ONT fragmentation mix

 **FRA (Fragmentation Mix) Oxford Nanopore**

-  **2.5 µl** **Technologies Catalog #SQK-RAD004**


-  **5.6 µl** **nuclease free water Contributed by users**

5 Prepare ONT loading buffer

 **FLT (Flush Tether) Oxford Nanopore**

- Add  **30 µl** **Technologies Catalog #EXP-FLP002**

into tube of

 **FB (Flush Buffer) Oxford Nanopore**

Technologies Catalog #EXP-FLP002

- Mix and spin down

6 Prime flow cell with  **800 µl** of loading buffer and let incubate for **00:05:00**

5m

7 Start sequencing run in MinKNOW software and pause run immediately after flow cell QC check

8 Pre-heat thermocycler initiate, then pause, combined Ultra-Rapid Program

19m

DNA Extraction

 **65 °C** for **00:01:00**

 **98 °C** for **00:02:00**

hold at  **20 °C**

LAMP

 **65 °C** for **00:14:00**

Fragmentation

 **30 °C** for **00:01:00**

 **80 °C** for **00:01:00**

hold at  **20 °C**


9 When acquired, place  20 µg tumor tissue into DNA extraction tube

10 Vortex on high for  00:00:30

30s


11 Place in pre-heated thermocycler and unpause combined program

12 Let thermocycler run through the DNA extraction protocol

13 Either let thermocycler draw DNA down to room temperature or quench in ice for  00:00:10 . Pause thermocycler program to pre-heat for LAMP. Briefly vortex on high and spin down.

10s


Target Amplification

14 Add  1 µl extracted DNA to LAMP reaction tube, being careful to avoid leftover tissue or other debris in tube.

15 Place in thermocycler and unpause combined protocol

16 Let thermocycler run through the LAMP protocol

26m

17 Either let thermocycler draw DNA down to room temperature or quench in ice water for  00:00:10 . Pause thermocycler program to pre-heat for ONT rapid library preparation tagmentation.

10s

Library Preparation

18 Add  1.9 µl LAMP product to tagmentation mix, flick to mix, and spin down


19 Place in thermocycler and unpause combined protocol

20 Let thermocycler run through the tagmentation protocol

21 Either let thermocycler draw DNA down to room temperature or quench in ice water for  00:00:10

10s

22

 [RAP \(Rapid Adapter\)](#) **Oxford Nanopore**

Add  **1 µl** **Technologies Catalog #SQK-RAD004**

to

tagmentation tube, flick to mix, and spin down

23

Incubate sequencing library for  **00:02:00** at room temperature, flicking to mix occasionally and spinning down ^{2m}

23.1

During the final 1m of incubation, re-prime the ONT flow-cell with SpotOn port open, using  **200 µl**

 [FB \(Flush Buffer\)](#) **Oxford Nanopore**

of **Technologies Catalog #EXP-FLP002**

mix

Sequencing

24

Add the entire 11 µl library to the sequencing mix tube

25

Pipette entire  **75 µl** sequencing mix into the MinION Spot On port, by slowly squeezing droplets onto the ramp/port.

26

Unpause sequencing run in MinKNOW software