




DEC 15, 2023

Ea(52-68) uptake assays

Rebecca
Wallings¹

¹University of Florida

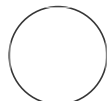

Rebecca Wallings

ABSTRACT

Ea(52-68) uptake assays on pMacs

OPEN  ACCESS



DOI:

dx.doi.org/10.17504/protocols.io.14egn3r3pl5d/v1

Protocol Citation: Rebecca Wallings 2023. Ea(52-68) uptake assays. **protocols.io** <https://dx.doi.org/10.17504/protocols.io.14egn3r3pl5d/v1>

License: This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](#), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

Protocol status: Working
We use this protocol and it's working

Created: Dec 14, 2023

Last Modified: Dec 15, 2023

PROTOCOL integer ID:
92327

Funders

Acknowledgement:

ASAP

Grant ID: ASAP-020527

- 1 Reconstitute MHC-II Ea chain (Ea) (52–68) peptide (Anaspec) in sterile distilled H₂O to a final concentrate of 1mg/mL.

Once peritoneal macrophages are adhered to plates, add 5µg per well of 24-well plate and incubate for 18-hours

Harvest cells as per flow cytometry protocol and stain with eBioY-Ae antibody (ThermoFisher) to assess Ea52-68 peptide bound to I-Ab