


POSTNATAL
VENTRAL
MIDBRAIN
DOPAMINE
NEURONAL
CULTURE
PROTOCOL

May 18, 2021

Postnatal ventral midbrain dopamine neuronal culture protocols

Dave Sulzer¹, Ellen Kanter¹¹Columbia University

1 Works for me

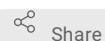

Share

dx.doi.org/10.17504/protocols.io.bpaamiaie

Yaqian Xu

ABSTRACT

This protocol details the specific steps and media solutions required to dissect and triturate the postnatal ventral midbrain dopamine neuronal cultures on mice.

EXTERNAL LINK

<https://s3.amazonaws.com/protocols-files/public/b5b2fc87b46ef008096859ffe7b0cfa338a5b35c7c4bc7c5b758e95c7ca4816a/cprvbg4xf.pdf>

ATTACHMENTS

[cprvbg4xf.pdf](#)

DOI

dx.doi.org/10.17504/protocols.io.bpaamiaie

EXTERNAL LINK

<https://s3.amazonaws.com/protocols-files/public/b5b2fc87b46ef008096859ffe7b0cfa338a5b35c7c4bc7c5b758e95c7ca4816a/cprvbg4xf.pdf>

PROTOCOL CITATION

Dave Sulzer, Ellen Kanter 2021. Postnatal ventral midbrain dopamine neuronal culture protocols.
protocols.io
<https://dx.doi.org/10.17504/protocols.io.bpaamiaie>

KEYWORDS

mice, neuronal culture, postnatal, midbrain, dopamine

LICENSE

————— This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](#), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

CREATED

Oct 31, 2020

LAST MODIFIED

May 18, 2021

OWNERSHIP HISTORY

Oct 31, 2020 Megan Freund

Nov 24, 2020 Yaqian Xu

Acknowledgments

Thanks to Steve Rayport, who initiated this project, and Serge Przedborski for TH immunolabel and cryosectioning of the figure in the dissection section, and to Robert Burke for collaboration on growth factor studies in this system. Thanks also to the technicians in our lab who advanced these methods, particularly Viviana Davila, Johanna Bogulavsky, Gerald Behr, Anthonia Hananiya, Kester Phillips, Candace Makeda Moore, Carl St. Remy, Megan Antonelli, from our lab, and Geetha Rajendran of Steve Rayport's lab. Thanks to Christian Bjerggaard, who visited our lab to learn these techniques in April 2004, and revised and added to an earlier addition of this booklet.

Introduction

This manual is designed to instruct in a detailed manner the steps our group uses to produce postnatally-derived ventral midbrain cultures. The approach used here has been used in numerous studies from our and Stephen Rayport's labs (Sulzer and Rayport, 1990; Rayport et al., 1992; Sulzer et al., 1993; Cubells et al., 1994; Rayport and Sulzer, 1995; Pothos et al., 1996; Przedborski et al., 1996; Sulzer et al., 1996; Fon et al., 1997; Mena et al., 1997a; Mena et al., 1997b; Burke et al., 1998; Pothos et al., 1998; Sulzer et al., 1998; Kholodilov et al., 1999; Pothos et al., 2000; Sulzer et al., 2000). In addition, these cultures have been used in studies from Louis Trudeau's lab (Bourque and Trudeau, 2000; Michel and Trudeau, 2000; Congar et al., 2002), Susan Amara's lab (Prasad and Amara, 2001; Ingram et al., 2002) and Mark Cookson's groups (Petrucelli et al., 2002).

Previously, dissociated culture studies of dopamine neurons use embryonically-derived cultures that produce a fraction of dopamine neurons of 1% or less. The postnatally derived cultures developed by Stephen Rayport and myself when I was a postdoc in his lab, have several advantages, the most striking being that the fraction of dopamine neurons is between 20-70%. The fraction is chiefly depending upon the area dissected.

Other advantages are that the neurons are more mature than embryonically-derived neurons. Not only do many more display tyrosine hydroxylase (TH), but the dopamine uptake transporter appears to be more highly expressed in the postnatally-derived neurons. Under optimal conditions, these cultures have survived for as long as 6 months, although most survive about 2 months. We suspect that the death is mostly due to increased osmolarity of the medium, and could probably be controlled by weekly partial exchange with conditioned media, testing the osmolarity of the medium over time with an osmometer and compensating with the addition of water, using a dish lid that allows air transfer but not evaporation. The latter approach has recently been introduced by Steve Potter, but has not been attempted to date for the postnatal dopamine neurons.

The techniques shown here have also been adapted by our group to produce postnatal cultures of other brain areas, such as the dorsal striatum. However, we suspect that the dopamine neurons are particularly challenging to culture, mostly because dopamine easily oxidizes to produce fairly potent neurotoxins. An odd choice for a neurotransmitter; perhaps we will eventually discover why the property of auto-oxidation to a reactive substance provides a selective advantage. Alternatively, maybe catecholamines were chiefly useful as skin and fur pigment, and simply adapted to produce additional neurotransmitters!

Overview of Techniques


First, we outline how we prepare the glass coverslip petri dishes in our lab. While we have used various approaches, including multiwell plates, to good effect, we generally culture on glass coverslips that are glued to petri dishes, so that the neurons produce small, dense cultures. This appears to keep them healthier, and to require no feeding with new media. The glass allows optimal microscopy with DIC/Nomarski imaging and fluorescent imaging. The central glass well allows all of the neurons to be approached with electrodes. We also make neuronal cultures on free-floating coverslips, so they can be removed from the dish for microscopy later. An alternate approach we also outline is the use of Aclar coverslips, which are well adapted for processing for electron microscopy.

We have found it necessary to culture these neurons on glia monolayers. We have also attempted the "Banker" technique, in which glia are grown close by on an inverted coverslip. However, to date, the dopamine neurons die

within 3 days unless they actually touch the glia. Thus, we outline the techniques to produce glia monolayers followed by a protocol how to dissociate the neurons and plate them onto the glia.

We use similar techniques for culturing of other types of neurons (hippocampal, cortical, striatal, etc.).

Figure indicating dissection of an TH immunolabel in mouse sections

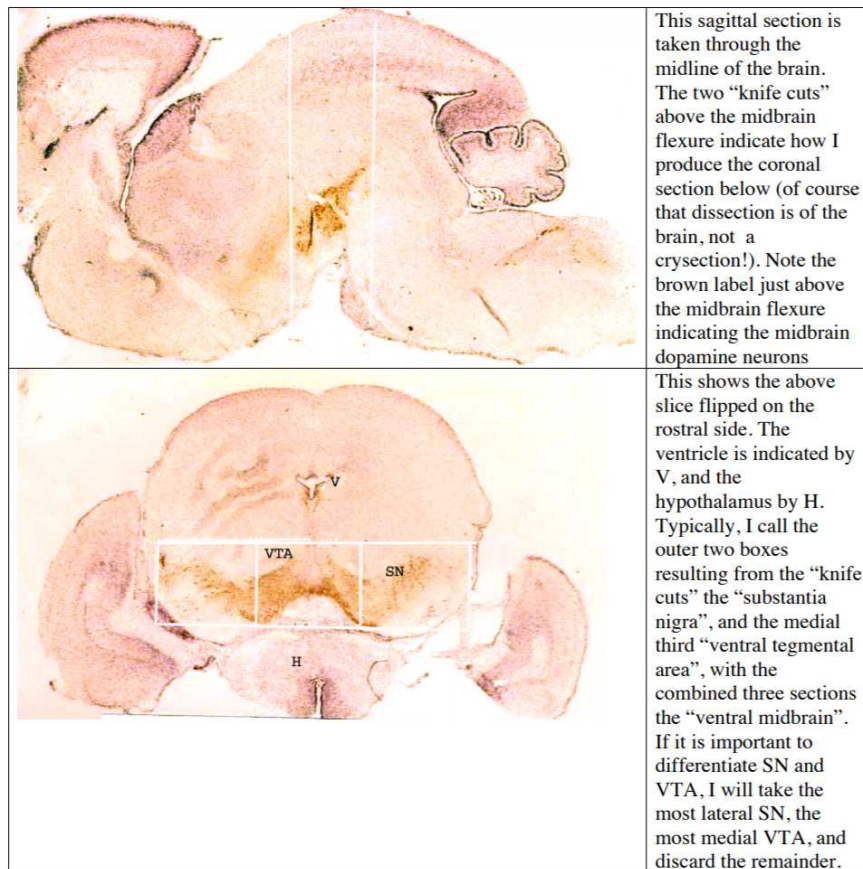

Special Instructions for Mouse Neurons

Some types of mouse tails (for ex, VMATs) need to be saved for genotyping. Prepare a numbered 1.5ml tube for each mouse. After the mouse has been decapitated, remove the tail using a fresh, sterile razor blade and place tail in numbered tube. Put all tubes on ice during the dissection. After dissection, store tubes in -80°C freezer. Alternatively, number the pups with the marker on their back, cut the tails and genotype the pups before the cultures, and no later than 1- day old. This way the pups with the same genotype can be cultured together rather than individually.

MATERIALS TEXT

Glass Coverslip Prep Materials

- Coverslips
- Beaker, 200 ml
- 95% EtOH
- Forceps
- Poly-D-Lysine, MW 70,000-150,000
- Water, “nanopure” or tissue culture grade
- Glass dish

Dissection Materials

- Rat pups, P1-P3
- Ketamine (100mg/mL)
- Dissecting scope & light
- Laminin, 70 µg aliquot
- 3 cc Syringe

- Dissecting instruments
- Aluminum foil pouch for bodies
- Ice chips in bucket
- Aluminum foil square for heads
- Chilled, sterile PBS
- 70% EtOH
- Scissors
- 18 gauge needles
- sterile disposable filters
- sterile microstir bars
- Petri dish [Falcon 100x15] w/ Sylgard square stuck on

Neuronal Preparation on Coverslips with Glia Monolayer Materials

- Magnetic Stirrer w/ 37°C water bath
- Carbogen, which is a mixture of 5% CO₂ and 95% O₂
- Papain solution
- 18 Gauge needles
- Micro stir bars (Fisher #14-511-68)
- 0.22 µm syringe filters
- Dissociation vials, Nunc specimen vials
- Dissecting instruments:
 - 2 fine forceps and tooth forceps
 - fine scissor
 - larger scissors
 - weighing scoop
 - scalpel with size 11 blade
 - hemostat
- 70% EtOH in a 50 ml tube
- Sterile slide rings
- Cold, sterile PBS
- Sterile, conditioned prewarmed SF1C

Microglial Culture Materials

- Rat or Mouse pups, P1-P3
- Ketamine (100mg/mL)
- 3 cc Syringe
- Dissecting instruments
- Aluminum foil pouch for bodies
- Ice chips in bucket
- Aluminum foil square for heads
- Chilled, sterile PBS
- 70% EtOH
- Scissors
- 18 gauge needles
- sterile disposable filters
- sterile microstir bars
- Dissecting scope & light
- Petri dish [Falcon 100x15] w/ Sylgard circle stuck on

M10C-G (for glia)

200 ml (can be made the day before):

Component	Amount
MEM	180 ml
Calf serum fresh, not heat-inactivated	20 ml
Glucose 45%	1.5 ml
Pen-strep	0.24 ml
Insulin (25 mg/ml)	40 µl, dissolve in 20 mM HCl
Glutamine (200 mM)	0.5 ml

1. Use a bent neck tissue culture flask to measure and mix solutions.
2. Filter using a 0.22 µm filter unit, label, and store at 4°C till needed.
3. Please note that L-Glutamine stock (200mM) is for most common culture media recipes is 100x (final concentration in media 2mM). We use less – final concentration 0.5mM.

SF1C (for neurons)

200 mL (can be made the day before):

Ingredient	Amount	Notes
BSA	0.50 g	fraction V. 0.25% final conc.
MEM liquid	94.0 ml	Sigma
DME liquid	80.0 ml	Sigma
F-12 liquid	20.0 ml	Sigma
Glucose 45% liquid	1.5 ml	Sigma solution
Glutamine 200 mM	0.50 ml	Aliquoted Sigma Solution
DiPorzio Conc.	2.0 ml	see below
Calf Serum	2.0 ml	Gibco, heat- inactivated. 1% final conc.
Liquid Catalase	0.10 ml	Sigma solution
Kynurenic acid 0.5 M	200 ul	0.5 M kynurenic acid is 946 mg in 10 ml. Dissolve in 1N NaOH
HCl 5N	50 ul	-

1. Combine ingredients in a bent necked flask.
2. Shake till BSA goes into solution.
3. Filter using a 0.22 µm filter unit, label and refrigerate
4. Note – neurons can be cultured in serum-free media. Make sure to condition the serum-free SF1C over astrocytes before using this media on neurons.

Laminin

NEED: Laminin 70 µg aliquot (enough for 50-70 dishes). Stock solution – mouse Laminin, 0.5mg/ml, from Chemicon, Sterile, cat # CC095. One 140ul-aliquot = 70ug.

15 ml centrifuge tube
7 ml cold sterile MEM

1. Dilute 1 aliquot laminin (70 µg) in 7 ml MEM. Final concentration 10ug/ml.
2. Add at least 100 µl per round well, enough to cover the surface, 1 hour at room temp in the hood before plating cells (can be done overnight at +4°C).

Glia Derived Neurotropic Factor (GDNF)

Human Recombinant GDNF – Chemicon, cat # GF030, 10ug.

1. Add 4.8ml of H₂O (tissue culture grade, Sigma # W3500) to 10ug of GDNF
2. Aliquot (76.9µl per tube – enough for 8 cell dishes) and freeze at -20°C
3. Add 725µl SF1C media to a tube with 77-µl aliquot and add 100µl of this solution to each cell dish (per 2ml media). This procedure results in a final concentration of 10ng/ml media.

5-Fluorodeoxyuridine (FDU)-solution

Make 1000X Stock:

<i>Ingredient</i>	<i>Amount</i>	<i>Final Conc.</i>
Uridine	247 mg	16.5 mg/ml
5-FDU	100 mg (bottle)	6.7 mg/ml
Tissue Culture Water	15 ml	-

1. Make a little over 15 ml of 16.5 mg/ml uridine
2. Add 15 ml uridine solution to 100 mg bottle of FDU to make 6.7 mg/ml solution FDU.
3. Divide FDU stock into 200µl aliquots and freeze at -20°C.
4. **(Adding FDU solution to cultures):** Add 1.8 ml of MEM into the 200 µl aliquot of FDU stock. Add 20 µl of diluted FDU solution to each dish (20µl per 2ml media).

Calf Serum

Use Heat Inactivated calf serum for neurons (SF1C), and Not-Heat Inactivated calf serum for glia (M10C-G).

1. Thaw serum in the refrigerator overnight.
2. Next morning, finish thawing in room temperature water bath.
3. Put serum on counter, raise bath temperature to exactly 56°C. Then place serum in bath for 30 min exactly, mixing every 10 min by turning bottle.

Papain Solution for GLIA

- prepare *ON the day of culturing*

10 ml solution:

Ingredient	Amount	Final Conc.
Cysteine Water	7.87 ml	1 mM cysteine
Papain	varies	20 units/ml
H&B Conc.	2 ml	-
0.5% Phenol red	20 µl	0.001%
Carbogen	95% O ₂ + 5% CO ₂	-

Papain Solution for NEURONS

-prepare *ON the day of culturing*

20 ml solution

Ingredient	Amount	Final Conc.
Cysteine Water	15.75 ml	1 mM cysteine
Papain	varies	20 units/ml
H&B Conc.	4 ml	-
HCL 5N	15 ul	-
0.5% Phenol red	40 ul	0.001%
Kynurenate 0.5M	20 ul	0.5 mM
Carbogen	95% O2 + 5% CO2	-

Notes

- Dissolve papain in cysteine water first by mixing well. Then add remaining ingredients. Otherwise, papain will not dissolve and will be removed by filtering.
- Make the Papain solution just before dissection, sterilize by filtering through 0.22u-syringe filter into NUNC dissociation vial, and effuse with Carbogen at least 30 minutes before you are ready to dissect (judge by the change of color, from pink-red to red-orange).
- Check the H&B concentrate before adding to make sure there is no crystal formation in solution. If there is, and you don't have time to make more, put solution in warm water bath for approximately 5 minutes or until crystals disappear.
- The Papain solution should be a light magenta color before effusion with Carbogen. Check periodically to make sure the magenta color is converted into a light salmon color (optimal pH for dissociation).
- For **Neurons, YOU MUST ADD** Kynurenate (kynurenic acid) and depending on the color, if it's too purple then add small quantities (<2ul) of 5N HCL at a time. **Do not add kynurenate or HCL to the papain solution for glia!**
- Use syringe filter to **STERILIZE** the Papain solution and put 2-10 ml of solution in each Nunc dissociation vial. Each vial could now hold a maximum of 7-10 animals' VM tissue. For Glia culture, you should only put 1 rat cortex in each tube with 10ml due to the excessive amount of tissue.
- Always make sure to label the vials with the **type** of cells and the **time** after the tissue put into the vials.
- The dissociation process should range from 7 to 20-25 minutes at the most. Keep an eye on the cells periodically. Do not let the stir bar spin too quickly and adjust the stirring to slow down the dissociation process, if necessary.

DiPorzio Concentrate

Single substance stocks to make up the DiPorzio Conc. Stock (see below)

Need:			Combine:		Aliquot:			
Additive	Solvent	Tube	Amount	ml	ml/tube	Conc	Amount	#Aliq.
Insulin	20mM HCl	plastic	250 mg / bottle	10	1	25 mg/ml	25 mg	10
Transferrin	Hank's		500 mg / bottle	5	1	100 mg/ml	100 mg	5
SOD	Hank's		70 mg / bottle	14	1	5 mg/ml	5 mg	14
Putrescine	Hank's		50 mg	2.5	0.12	20 mg/ml	2.4 mg	21
Na2SeO3	Hank's		0.104 mg	10	0.5	10.4 ug/ml	5.2 ug	20
T3	10mM NaOH		2 mg	10	0.1	0.2 mg/ml	20 ug	100
Progesterone	100% EtOH	glass	12.5 mg	10	0.05	1.25 mg/ml	62.5 ug	1
Cortisol	100% EtOH	glass	20 mg	10	0.02	2.0 mg/ml	40 ug	1

Notes:

SOD=superoxide dismutase, T3=3,3',5-triiodo-L-thyronine (Na salt) or Liothyronine; Stock solutions should be kept frozen at -80°C.

1. 20mM HCl = 41.5µl conc HCl (12N)/ 25ml water
2. Na2SeO3 - Make 1mg/ml stock and add 104 µl to 10 ml.
3. Keep stocks at -80°C

Combine the following to make up the DiPorzio Conc. Stock (100X) - 10 ml

Additive	Amount (ml)	Final Conc. in Stock	Final Conc. in SF1C Buffer		ml for two batches
Progesterone	0.05	6.25 ug/ml	62.5 ng/ml	200 nM	0.1
Cortisol	0.02	4 ug/ml	40 ng/ml	125 nM	0.04
Hank's BSS	6.21	-	-	-	12.42
Insulin	1	2.5 mg/ml	25 ug/ml	-	2
Na2SeO3	0.5	0.52 ug/ml	5.2 ng/ml	30 nM	1
T3	0.1	2 ug/ml	20 ng/ml	30 nM	0.2
SOD	1	0.5 mg/ml	5 ug/ml	-	2
Putrescine	0.12	0.24 mg/ml	2.4 ug/ml	15uM	0.24
Transferrin	1	10 mg/ml	100 ug/ml	-	2

Use a 15 ml polypropylene test tube (translucent) and perform quasi-sterility.

1. Add the progesterone and the cortisol.
2. Use an aspirator pipet and vacuum to speed evaporation of EtOH: Use a 5ml pipet broken in half and place half-way down into the tube being very careful not to aspirate EtOH liquid. Hold pipet securely in place with Kimwipes and then bring the tip down to the 500µl mark located on the side of the tube.
3. Add the subsequent aliquots in order
4. After the addition of the insulin which makes the solution cloudy, add 20 µl of 1M NaOH to neutralize the pH. Solution should go from yellow to pink and immediately clear. ALSO, after the addition of transferrin, immediately add 20µl of 1N NaOH, to neutralize the pH and prevent the formation of precipitates.
5. Draw up 10 ml into a serological pipet and divide into 5 aliquots, 2ml each (one batch). Store at -20°C.

H&B Concentrate

(100 ml of 5X)

Ingredient	M.W.	Powder/50 ml H2O	Conc. (M)	Combine (ml)	Final Conc. (mM)
NaCL	58.44	11.70 g	4	14.5	116
KCl	74.56	3.728 g	1	2.7	5.4
NaHCO3	84.01	4.2 g	1	13	26
NaH2PO4·H2O	137.99	6.90 g	1	1	2
MgSO4	120.38	6.02 g	1	0.5	1
EDTA	292	Sigma 5%	0.17	1.46	0.5
Glucose	180	Sigma 45%	2.5	5	25
Tissue culture H2O	-	-	-	61.84	-

1. Combine amounts from stock solutions.
2. Filter sterilize (0.22 µm), divide into 4 ml aliquots and freeze at -20°C.
3. Filter sterilize (0.22 µm), divide into 4 ml aliquots and freeze at -20°C.

Cysteine Water

157.5 ml of 1.25X

1. Make stock of 0.5 M CaCl2 and stock of 0.02 M cysteine in water.
2. Filter sterilize into bottle and keep cold at +4°C.
3. Combine 146.9ml water, 10 ml 0.02M cysteine stock and 0.6 ml of 0.5 M CaCl2.

Components						
Ingredient	M.W.	Powder (mg)	H2O (ml)	Stock Conc. (mM)	Combine (ml)	Final Conc. (mM)
CaCl2	111	555	10	500	0.60	1.9
Cysteine	121.2	24	10	19.8	10	1.25
TC water	-	-	-	-	146.9	-

CATALOG NUMBERS

L-Glutamine – 200mM, Sigma # G2150

Pen-Strep - Sigma # P0781

Fetal Bovine Serum, heat-inactivated - Gibco # 16140-063

DMEM - Sigma # D5546

SIGMA CATALOG #'S/size		
A4503	BSA (albumin)	50 g
C2505	Cortisol	20 mg
C3155	Catalase (liquid)	50 mg
C7352	L-Cysteine (base)	25 g
D5546	DME	500 ml
F0503	5-fluorodeoxyuridine	
G8769	Glucose, 45% w/v	100 ml
G2150	L-Glutamine, 200 mM	20 ml
I5500	Insulin	250 mg
K3375	Kynurenic Acid	5 g
M2279	MEM	1000 ml
N4888	Nutrient Mixture F-12 (Ham's)	500 ml
D1408	PBS 10x	500 ml
P0781	Pen/Strep	20 ml
P0290	Phenol Red	100 ml
P0899	Poly-D-Lysine, MW70,000-150,000	10 mg
P0130	Progesterone	5 g
P7505	Putrescine	50 mg
S1382	Na ₂ SeO ₃	10 g
S2515	Superoxide dismutase (SOD)	70 mg/300,000 units
T2752	T3	2 mg
T1428	Transferrin	500 mg X2
W3500	Water, Tissue Culture Grade	1000 ml
U3003	Uridine	

OTHER CATALOG #'S		
10501-10	Aclar film	Ted Pella, Inc.
CC095	Laminin (mouse)	Chemicon
12-565-206	Dissociation vial	Nunc (Fisher)
6705-R12	Microscope slide rings	Thomas Scientific
PAP	Papain	Worthington
14170-120	Hank's Balanced Salt Solution with Phenol Red	Gibco
26140-087	Calf Serum	Gibco
16140-063	Calf Serum, heat inactivated	Gibco
184	Sylgard	Dow Corning Corporation
14-511-68	Micro stir bars	Fisher
351006	Petri dish 50x9mm	Falcon
63-3111	Cover slips #2	Carolina Biological Supply
GF030	GDNF	Chemicon
12-565-224	NUNC Bioassay trays	Fisher
N0 XX	Portable Light Duty Puncher	Roper Whitney

SAFETY WARNINGS

See the Safety Data Sheet (SDS) for any hazards or safety warnings.

BEFORE STARTING

Three Days Before Culturing

1. Make sure that there is an ample supply of glass dishes, Sylgard circles, Techtips and both types of plastic pipette tips (for 1000uL and 200uL).
2. Figure out the number of dishes and types of neurons required. Calculate the number of pups needed. Pups needed vary by dissection technique, but we have found 1 rat cortex is enough material to easily plate ~110-130 dishes, and one mouse VM enough to plate 2-3 dishes.
3. Clean and sterilize dissection tools (This should be done right after each dissection!)
 - a. First wipe clean with 70% ethanol.
 - b. Wipe again with 100% ethanol to speed drying process.
 - c. Let dry!
 - d. Leave in the oven at 200°C for at least 1 hour.
4. Clean and sterilize the plastic rings and magnetic stir-bars:
 - a. Wipe off any debris from rings and stir-bars.
 - b. Let soak in 1N HCL for 1 hour.
 - c. Rinse at least 20 times with de-ionized water and let dry.
 - d. Place rings into 100 ml beakers and individually distribute the stir-bars into small tubes (the 1.5 ml snap-top tubes are fine) and autoclave.
5. Once-a-month, autoclave the flask and glass tube (NOT THE RUBBER STOPPER) that is used for bubbling the carbogen through water for papain dissociation.

On the Day Before Culturing

1. Prepare appropriate media as needed (SF1C, M10C-G, PBS).
2. *For neuronal culturing*, take glial dishes and wash them twice with 2ml of cold MEM. Replace MEM with 2 ml of SF1C. Alternatively, neurobasal medium A/B27 can also be used (Table 16.9). While A/B27 is easier to make, we have had slightly better success with SF1C in certain toxicity assays (unpublished observations).

3. Place sterile ring in each dish to the side of the glass coverslip well, using forceps with sterile tips. (see dish and media prep. for neurons). For free-floating coverslips, place the ring in the middle of a coverslip.
4. Leave dishes with SF1C in the incubator overnight. This allows the media to condition over the glial cells. Alternatively, use preconditioned SF1C media on the day of the culture.

TO CONDITION MEDIA:

- have T-225 flask with astrocytes monolayer prepared at least 10 days prior;
- replace media in the flask with freshly made 200ml of M10CG media. Incubate for 1-2 days. Collect conditioned media and filter to sterilize.

To condition neuronal media, add 200ml of SF1C media to the flask with astrocytes and incubate overnight to 2 days. Collect conditioned media and filter-sterilize.

Aclar Dishes for Electron Microscope Preparations

5m

5m

- 1 Cut Aclar sheets into squares ~ 1.5 cm and drop into Nanopure dd H₂O for ⌚ 00:05:00 .

- 2 Swirl the squares in a beaker (if they don't all immediately float to the top).

- 3 Rinse twice in Nanopure dd H₂O for 5 min:

- 3.1 Rinse 1/2: rinse in Nanopure dd H₂O for ⌚ 00:05:00 .

5m

- 3.2 Rinse 2/2: rinse in Nanopure dd H₂O for ⌚ 00:05:00 .

5m

- 4 Expose to 70% ethanol in Nanopure dd H₂O for ⌚ 00:30:00 .

30m

- 5 Rinse twice in Nanopure dd H₂O for 5 min:

- 5.1 Rinse 1/2: rinse in Nanopure dd H₂O for ⌚ 00:05:00 .

5m

- 5.2 Rinse 2/2: rinse in Nanopure dd H₂O for ⌚ 00:05:00 .

5m


- 6 Expose to **1 Molarity (M) HCl** in Nanopure dd H₂O for ⌚ 00:30:00 .

30m

7 Rinse again several times in water.

8 Maintain in 70% ethanol, covered.



9 Sylgard the squares to the dishes after the holes are prepared:

9.1 Wearing gloves, prepare about  **5 mL Sylgard** in a large weighing boat by combining 10 parts resin with 1 part catalyst (by weight) and stirring extensively with flame-sealed Pasteur pipette. Pour contents into new weighing boat to avoid use of resin without catalyst.

9.2 Dip an empty 15 ml Falcon tube or glass scintillation vial upside-down into the Sylgard, using it as a rubber stamp to make a circle of Sylgard around the hole, on the outside of the dish.





9.3 Place the Aclar sheets on the circle of Sylgard, covering the hole.

9.4 Allow Sylgard to spread (about  **00:15:00**). 15m

9.5 To cure the Sylgard, place dishes, upside-down, at  **37 °C**  **Overnight** , or at room temperature 15m
for 2 to 3 days (cured Sylgard is no longer sticky to the touch).

9.6 Store dishes in zip lock bags: initial and date.

10 Coat with PDL:

10.1 Add  **200 µl 40 µg/ml PDL solution** of the to each dish. Expose for  **02:00:00** - 6h
 **04:00:00** at  **Room temperature** in a culture hood without UV light.

10.2  5m

Rinse 3 times for 5 min each in tissue culture water:

Rinse 1/3: rinse in tissue culture water for  **00:05:00** .

10.3 Rinse 2/3: rinse in tissue culture water for 🕒 00:05:00 . 5m

10.4 Rinse 3/3: rinse in tissue culture water for 🕒 00:05:00 . 5m

10.5 Remove water and expose to UV for 🕒 01:00:00 to sterilize. 1h

Note that these dishes must be coated with laminin before use, just as we do with glass coverslips (see glia dissection protocol)

Left over sylgard from dish preparation can be used to make sylgard circles.

Glass Coverslip Preparation 1m

11 Drop 100 coverslips (15 x 15 mm) one-at-a-time into a 200 ml beaker containing 📦 100 mL 95% EtOH . (Leave approx. 🕒 00:01:00 at least) 1m

12 Prepare a solution of poly-d-lysine in "nanopure" or tissue culture grade water. The final concentration should be 📊 40 microgram per milliliter (µg/mL) and 250mL are adequate for 100 coverslips. Place the poly-Dlysine solution in a glass dish of enough area to give 3/4"/2 cm height to the poly-D lysine solution. (the walls of the dish should be at least 1.5"/4 cm, the area no less than 40"/100 cm square)

13 Throw (pour) away excess EtOH. Place ethanol soaked cover slips on to a paper or glass container.

14 Individually flame dry the coverslips (one pass through flame); use forceps.

15 After putting the slip through the flame, allow it to cool in the air in forceps for approximately 🕒 00:00:10 while moving back and forth (cracking or sizzling sounds mean it did not cool enough). 10s

16 Put the coverslip into the poly-D-lysine solution. Repeat for all slips.

17 Let sit for 🕒 01:00:00 , then pour off the poly-D-lysine, then let slips dry (🕒 Overnight in the hood with the blower on). 2h

18 Prepare dishes with 10 mm holes punched with a bench punch for wells where neurons and glia will grow. Sylgard the coverslips to the dishes after the holes are prepared.

18.1 Wearing gloves, prepare about **5 mL Sylgard** in a large weighing boat by combining 10 parts resin with 1 part catalyst (by weight) and stirring extensively with flame-sealed Pasteur pipette. Pour contents into new weighing boat to avoid use of resin without catalyst.

18.2 Dip an empty 15 ml Falcon tube or glass scintillation vial upside-down into the Sylgard, using it as a rubber stamp to make a circle of Sylgard around the hole. Use caution when applying the Sylgard as too much will result in excess Sylgard covering the culturing surface of the coverslip whereas inadequate Sylgard will result in a leaky well.

18.3 Place the coverslips on the circle of Sylgard, covering the hole.

18.4 Press down on coverslip gently using forceps to eliminate bubbles.

18.5 Allow Sylgard to spread (about **00:15:00**). 15m

18.6 To cure the Sylgard, place dishes, upside-down, at **37 °C Overnight**, or at room temperature 15m for 2 to 3 days (cured Sylgard is no longer sticky to the touch). Once the dishes are cured they can be stored in zip-lock bags at **4 °C** until needed.

Note that these glass coverslips must be coated with laminin before use (see glia dissection protocol)

Left over sylgard from dish preparation can be used to make sylgard circles.

If these coverslips ever become unavailable, the Bellco #2 photoetched coverslips have been recommended

These are used for triturating cells

Obtain a box of 1000 ul pipette tips (unfiltered). Wear gloves.

- 20 Flame all tips in order to seal off hole (we recently started using tips without sealing the hole, making additional holes with needles as described below).
- 21 Using a needle, poke a hole as close to the tip as possible through one side and out the other. This will result in two holes. The purpose of piercing right through both sides is to prevent cells getting caught in the pipette. Rather, they will have an outlet hole. Make 3 types of tips, with Large (21G 1½), Medium (22G 1½) and Small (26G 1½) size holes.
- 22 Autoclave as usual.

Sylgard Circles

23

For use during dissection as a platform to do the cutting:

Weigh out the Elastomer Base and the Elastomer Curing Agent. Mix in the ratio 10:1

- 24 Mix well and pour into 12 micro well plate, about 3-4 mm in height.

- 25 Cure in oven at 🔥 60 °C for ⌚ 00:30:00 or at 🔥 Room temperature for 5-7 days.

30m

- 26 Cut circles out of wells using razor and individually wrap each circle in aluminum foil, autoclave to sterilize.

Left over sylgard from dish preparation can be used to make sylgard circles.

Glia Preparation on Coverslip Petri Dishes

27 LAMININ COATING

Before the dissection/digestion time, dilute 1 aliquot of laminin (140ul of 0.5mg/ml laminin = 70 µg) in


📄 7 mL sterile MEM . Final concentration [M] 10 microgram per milliliter (µg/mL) .

- 28 Add at least 📄 100 µl 10 µg/mL laminin per round well (enough to cover the surface) for minimum ⌚ 01:00:00^{1h} at 🔥 Room temperature (can be done O/N at +4°C) before plating cells and leave in flow bench. Aspirate off laminin and wash each well with MEM. Place into 🔥 37 °C incubator to warm dishes before plating the cells. 1 aliquot is enough for about 50-60 dishes.


29 DISSECTION


Obtain one or two P1-P3 rats.

- 30 Prepare the papain dissociation solution. Each vial should hold tissue from a maximum of 2 animals though it's better to digest 1 rat cortex per vial.
- 31 Sterilize the top of each Nunc specimen vial (25 ml) with ethanol. Punch two holes in the top with an 18 gauge needle, using a separate needle to bore holes. In one hole place a new needle with a 0.22 μ m filter attached. Leave the other hole empty for ventilation.
- 32 Add an HCl-cleaned, autoclaved micro-stir bar to the vial. Fill each tube with **5 mL** - **10 mL freshly made papain solution** (sterile filtered with 0.22 μ m filter). Place in temperature controlled **37 °C** water bath with magnetic stirrers.
- 33 Arrange dissection tools, microscope, dissection light, a bunch of transfer pipettes, sterile Petri dish with Sylgard circle, beaker, and 70% ethanol in a 50 ml tube in the sterile hood.
- 34 Prepare two buckets of ice to chill the PBS and for the pups.
- 35 Cut aluminum foil squares (about 15 x 15 cm).
- 36 Perfuse dissociation tube (papain tube) with a steady flow of humidified carbogen (95% oxygen, 5% CO₂) delivered through the filter/needle (flow rate adjusted to about 1 bubble per second).
- 37 When the papain has reached a pH of **pH7.2** - **pH7.4** in about **00:30:00** - **01:00:00** ^{1h 30m} (indicated by a red color; purple indicates that the papain solution is too alkaline and orange that it is too acidic), anesthetize pups and begin dissection.
- 38 Anesthetize pups with an intraperitoneal injection of ketamine (0.01 and 0.05 ml of 100 mg/ml for mouse and rat pups, respectively).
- 39 When pups begin to show sedation, put them on ice till hypothermic.
- 40 Rinse aluminum foil square, scissors and head with 70% EtOH.
- 41 Decapitate pup, allowing head to fall onto aluminum foil square and move to hood.

- 42 Using (toothed) forceps hold head by the eyes (using left hand if you are right handed).
- 43 Using the other hand and a curved or angled forceps, pinch the scalp just behind the eyes and pull back, tearing it off.
- 44 Next, use a micro-scissor to cut around the circumference of the skull, and gently peel the skull off with forceps, being careful in case any tissue still connects it to the rest of the skull.
- 45 Gently remove brain (using a small scoop shaped scalpel) and place into Petri dish with one SYLGARD® circle stuck onto the bottom and filled with ice cold sterile PBS (use enough PBS such that the brain lying on the SYLGARD® circle is submerged)
- 46 Remove the olfactory lobes and brain stem (with V-cut). Hemi-sect brain and then peel off meninges. Shell out cortex, taking care to remove the hippocampus. Remove any torn areas to avoid DNA release. Cut into 1 mm³ chunks.
- 47 **DIGESTION** 15m
Place chunks in papain as with neuronal prep, *without kynurenate* (1 cortex/vial). Chunks will be digested enough to triturate during the initial ⌚ 00:05:00 - ⌚ 00:10:00 . We judge by the appearance of the chunks (like cotton candy) when they are ready.
- 48 **TRITURATION AND PLATING**
Warm up 🧊 45 mL sterile M10C-G for rinsing and trituration.
- 49 From this point forward, everything should be done in a sterile laminar hood with sterile technique.
- 50 Remove chunks from enzyme vial using transfer pipette and place in 15 ml conical test tube.
- 51 
Rinse chunks 3 times using 🧊 2 mL M10C-G , allowing them to settle each time and discard supernatant:
- 51.1 Rinse 1/3: Rinse chunks using 🧊 2 mL M10C-G , allow them to settle and discard supernatant.
- 51.2 Rinse 2/3: Rinse chunks using 🧊 2 mL M10C-G , allow them to settle and discard supernatant.

51.3 Rinse 3/3: Rinse chunks using  **2 mL M10C-G** , allow them to settle and discard supernatant.

52 Using  **1 mL warm M10CG** triturate up to a maximum of 10 times (3-5 times should be sufficient) using a regular 1 ml pipette tip.

53 Remove supernatant to a separate tube, and replace the volume with  **1 mL fresh M10C-G** .

54 Repeat trituration using the tech tips with smaller holes each time until the chunks have all fallen apart.

55 

Pellet the glial cells at  **1000 rpm, 00:05:00** .

56 . Resuspend pellet in  **10 mL M10C-G** , and count the cells using hemocytometer

57 Dilute to plating density of 1,000,000-1,500,000 cells/ml.

58 Remove laminin from wells and plate about 80,000 cells per well. Plate by dripping the cell suspension from about 10 cm above the dish to help minimize number of surviving neurons in glial prep.

59 Tap the sides of the trays (not for island cultures) in order to evenly distribute cell suspension.

60 

3h 30m


Allow cells to settle for  **01:30:00** -  **02:00:00** in incubator.

For island preps allow 2-5 hours.

61 Chill M10C-G in meantime.


62 

RINSING AND FEEDING


For monolayers: squirt  **2 mL cold MEM** into each well, moving the pipette around in a circular motion to wash each well evenly.

63 

Aspirate off MEM and repeat wash:


63.1 Squirt  **2 mL cold MEM** into each well, moving the pipette around in a circular motion to wash each well evenly.

63.2 Aspirate off MEM.

64 For island cultures: One dish at a time, lift the cover and tilt the dish at a 45° angle with the hood countertop. Using a transfer pipette very gently, drip  **2 mL cold MEM** onto a spot on the dish just above the well so that it slides down the incline and over the well.

65 Aspirate off the MEM.

66 Feed with  **2 mL cold M10C-G** .

67 Add FDU to the dishes when cells are ~70% confluent (about 3-5 days). If neurons are still present at this point, remove the medium and cold-feed with  **2 mL M10C-G** before adding FDU.

68 

Astrocytes dishes should be at least one week old and at most 4 weeks old at the time the neurons are plated. Feed with cold M10C-G at one month if not already used.




Neuronal Preparation on Coverslips with Glia Monolayer

69 Please see "Before Starting" in Guidelines & Warnings section on Things to Prepare during the Week before Plating Glia and Neurons.

70 First, turn on the dissociation water bath and let warm to  **37 °C** .

- 71 Make the papain dissociation solution (see Media section). Each vial can hold tissue from 1 to maximum 10 mice preps.

- 72 Sterilize the top of each Nunc specimen vial with EtOH. Make two holes in the top with an 18 gauge needle, using a separate needle to bore holes. In one hole place a new needle with a 0.22 µm filter attached. Leave the other hole empty for ventilation.

- 73 Add an HCl-cleaned, autoclaved micro-stir bar to the vial. Fill each tube with  **5 mL** -
 **10 mL freshly made papain solution** (sterile filtered with 0.22 µm filter). Place in temperature controlled
 **37 °C** water bath with magnetic stirrers.

- 74 Arrange dissection tools, microscope, dissection light, a bunch of transfer pipettes, sterile Petri dish with sylgard circle, beaker, and 70% ethanol in a 50 ml tube in the sterile hood.





- 75 Get 2 buckets of ice to chill the PBS and to anesthetize the pups.

- 76 Rip aluminum foil squares (about 15 cm²).

- 77 Fill “insulin” syringe with 100 mg/ml ketamine solution.

- 78 Obtain as many P1-P3 rat or mouse pups as needed (animals which are one to three days old)

- 79 Superfuse dissociation tube (papain tube) with a steady flow of humidified carbogen delivered through the filter/needle (Flow rate adjusted to about 1 bubble per second).

- 80 When the papain has reached a pH of  **pH7.2** -  **pH7.4** in about  **00:30:00** -  **00:45:00** ^{1h 15m} (indicated by the red color, not too purple, not too orange), anesthetization and dissection of pups can begin.

- 81 **DISSECTION**
 Anesthetize pups with an intraperitoneal injection of ketamine (0.01 ml and 0.05 ml of 100 mg/ml ketamine for mouse and rat pups, respectively).

- 82 When pups begin to show sedation, put them into ice till hypothermic.

- 83 Rinse aluminum foil square, scissors and head with 70% EtOH.

84 Decapitate pup, allowing head to fall onto aluminum foil square and move to hood.

Some types of mouse tails (for ex, VMATs) need to be saved for genotyping. Prepare a numbered 1.5ml tube for each mouse. After the mouse has been decapitated, remove the tail using a fresh, sterile razor blade and place tail in numbered tube. Put all tubes on ice during the dissection. After dissection, store tubes in -80°C freezer. Alternatively, number the pups with the marker on their back, cut the tails and genotype the pups before the cultures, and no later than 1- day old. This way the pups with the same genotype can be cultured together rather than individually.

85 Gently remove brain into ice cold PBS-filled Petri dish as described for glia. With anterior of brain facing left on the SYLGARD® circle (dissection will be described for right handed people) use a fine forceps in the left hand to pin the two hemispheres of the brain to the SYLGARD®. Remove any remaining meninges.

86 Place the brain, dorsal side down, anterior facing left onto the SYLGARD®. Using forceps, pin the brain to the SYLGARD® to immobilize.

87 Make the initial cut through the entire brain caudal to the midbrain flexure (right line, Fig 1 in Guidelines section), and the second cut rostral to the flexure, including the caudal aspect of the hypothalamus.

88 Lay the resulting slice flat on the SYLGARD® with the ventral side facing right, and the dorsal side to the left and the caudal aspect up.

89 Remove any remaining cortex.

90 Place the forceps through the aqueduct and push into the SYLGARD® to immobilize the slice.

91 Cut the ventral edge of the slice along the top of the hypothalamus as indicated by the line in Figure in the Guidelines section.

92 Next cut approximately halfway between the new ventral edge of the slice and the aqueduct.

93 Cut into 1 mm³ segments and place all three segments into the vial with papain solution for ventral midbrain cultures. For SN cultures only use the two outer segments and for VTA only use the inner segment. Use both SN and VTA for combined VM cultures.


94 Set the magnetic stirrer on low. Chunks should swirl around in the vial gently.


- 95 Digest for approx. ⌚ 00:07:00 - ⌚ 00:20:00 or until chunks have fallen apart and look something like cotton candy.^{27m} Periodically examine the tubes to see if segments are sticking to the stir bars; if so, tap vials to loosen tissue.

96 

DISH AND MEDIA PREPARATION

On day before dissection, take glia dishes and wash them twice with  2 mL cold MEM :

96.1 Wash 1/2: wash with  2 mL cold MEM .



96.2 Wash 2/2: wash with  2 mL cold MEM .

97 Replace MEM with  2 mL SF1C .

98 

20m

Leave dishes with SF1C in the incubator ⌚ Overnight . Alternatively, do this step on the day of the culture, using preconditioned media.

99 After dissection, place tubes containing SF1C for rinsing and triturating in the incubator until ready to triturate so that the pH stays between  pH7.2 and  pH7.4 , and the media is warm.

The media must be at the correct pH (indicated by the red color) or the neurons may die.

100 TRITURATION

From this point forward, everything should be done in a sterile laminar hood with sterile technique.

101 Remove chunks from enzyme vial using transfer pipette and place in 15 ml conical test tube.

IMPORTANT: Take great care to prevent bubble formation because bubbles are harmful to neurons. If bubbles occur, aspirate them off.

102 

Rinse chunks 2-3 times using about **2 mL SF1C** warmed to **37 °C** . **After adding SF1C for each wash, flick the tube so that the cells swirl around and shake off the papain.** Allow tissue to settle each time and discard as much of the supernatant as possible without exposing cells to air.

102.1 Rinse 1/2: Rinse chunks 2-3 times using about **2 mL SF1C** warmed to **37 °C** . **After adding SF1C, flick the tube so that the cells swirl around and shake off the papain.** Allow tissue to settle and discard as much of the supernatant as possible without exposing cells to air.

102.2 Rinse 2/2: Rinse chunks 2-3 times using about **2 mL SF1C** warmed to **37 °C** . **After adding SF1C, flick the tube so that the cells swirl around and shake off the papain.** Allow tissue to settle and discard as much of the supernatant as possible without exposing cells to air.

102.3 

Optional Rinse: Rinse chunks 2-3 times using about **2 mL SF1C** warmed to **37 °C** . **After adding SF1C, flick the tube so that the cells swirl around and shake off the papain.** Allow tissue to settle and discard as much of the supernatant as possible without exposing cells to air.

103 Triturate in **2 mL room temperature SF1C** with large-bore tech-tips a couple of times (when done properly the medium becomes slightly opalescent with dissociated cells).

104 Repeat the trituration process *up to three times* with a medium-, then a small-bore tech-tip. The chunks should now be completely dissociated into individual cells.

105 

Centrifuge cells at **1000 rpm, 00:05:00** .

106 Resuspend the pellet in 0.2-0.5 ml (for mouse cultures) or 1-2 ml (for rat cultures) of media, mix and count (see below).

107 Dilute cell suspension until **50 µl** - **100 µl** of media contains 80,000 cells.

108 Add the appropriate volume to the well (into the slide ring) of the SF1C-containing dishes so that each well contains 80,000 cells (use a circular motion to spread the cells around). Again, the rings should ensure that the neurons settle into the well and do not float away or get washed away when moving the dishes.

109 **ATTACHMENT** 1h 30m
After **01:30:00** , cells should have attached to the surface of the well. Remove rings from each well using forceps, which have sterile pipette tips at the end.

110 Add **100 µl GDNF (200 ng/ml)** to each culture for a final concentration of 10 ng/ml GDNF (see media and

solutions section how to prepare GDNF).

NEXT DAY - Mitotic Inhibition

111 The next day inhibit the growth of non-neural cells with 5-fluorodeoxyuridine (FDU). Dilute 1000X stock 1:10 by adding **200 μ l stock** to **1.8 mL MEM**.

112 Add **20 μ l diluted FDU solution** to each dish (2 ml media).

113 Media should be changed only when necessary (change of color) – aspirate half the media and add fresh preconditioned (pre-warmed at **37 °C**).

Try to avoid this step if possible as neuronal cells are very sensitive to the media exchange.

114 Disturb cultures as little as possible, avoiding multiple observations (especially in the first 2 days after culturing).

115 Process cells as needed for the individual experiment (ideally, between 5 days to 2 weeks old).

If necessary after fixation/immunostaining, the coverslips can be taken off the Petri dish (cut between the glass and Sylgard with a scalpel – careful not to break the glass!!!) and mounted on a glass slide. Alternatively, cells can be cultured on free-floating coverslips, which can be later removed from dishes for the experiments.

Microglial Culture 2w 0d 6h 27m

116


Prepare astrocyte cultures in tissue culture flasks (225 cm² for rat and 75 cm² for mice cultures)

Obtain one or two rat or five mouse pups P1-P3.

117 Anesthetize pups with an intraperitoneal injection of ketamine.


118 When pups begin to show sedation, put them **On ice** till hypothermic.

119 Rinse aluminum foil square, scissors and head with 70% EtOH.

- 120 Decapitate pup, allowing head to fall onto aluminum foil square and move to hood.
- 121 Gently remove brain into ice cold PBS filled Petri dish (with one Sylgard circle – should be covered with cold PBS!).
- 122 Remove the olfactory lobes and brain stem (with v-cut).
- 123 Hemi-sect brain and then peel off meninges.
- 124 Shell out cortex.
- 125 Remove any torn areas to avoid DNA release.
- 126 Cut into 1 mm³ chunks.
- 127 **DIGESTION** 17m
Place chunks in papain as with neuronal prep, *without kynurenate* (1 rat cortex/vial, or 5 mouse cortexes/vial). Chunks will be digested enough to triturate during the initial ⌚ 00:07:00 to ⌚ 00:10:00 . We judge by the appearance of the chunks (like cotton candy) when they are ready
- 128 **TRITURATION AND PLATING**
Warm up sterile M10C-G to ⌚ 37 °C for rinsing and trituration.
- 129 Remove chunks from enzyme vial using transfer pipette and place in 15 ml conical test tube
- 130 
Rinse chunks 3 times using 🧴 2 mL warm M10C-G , allowing them to settle each time and discard supernatant.
- 130.1 Rinse 1/3: Rinse chunks using 🧴 2 mL M10C-G , allow them to settle and discard supernatant.

130.2 Rinse 2/3: Rinse chunks using  **2 mL M10C-G** , allow them to settle and discard supernatant.


130.3 Rinse 3/3: Rinse chunks using  **2 mL M10C-G** , allow them to settle and discard supernatant.

131 Using  **1 mL warm M10CG** triturate 5-10 times using just the regular pipette tip.

132 Repeat trituration using the tech-tip with smaller size hole, until the chunks have all fallen apart.

133 



Pellet the glial cells at  **1000 rpm, 00:05:00** .

134 Resuspend pellet in  **15 mL M10C-G** .


135 Transfer cell suspension to tissue culture flask. Use big T-225 flask for rat culture and T-75 flask for mouse culture. The yield from two rat brains is ideal for one T-225 flask. And the yield from 5 mouse brains is ideal for one T-75 flask.

136 

4h

Allow cells to settle for at least  **02:00:00** , better leave them  **Overnight** in incubator.

137 **RINSING AND FEEDING**

For astrocyte monolayer remove M10CG from flask and squirt ~  **20 mL cold MEM** onto the base of upright tissue culture flask.

138 Gently allow the media to swirl around the bottom of the flask to remove unattached cells.

139 

Aspirate off MEM and feed with cold M10C-G (30 ml per T-225 flask, or 15ml per T-75 flask), and leave ~2 weeks in the incubator.

DO NOT ADD FDU TO FLASKS!!!!

140 MICROGLIA ISOLATION

Once the cell layers are stable and you start seeing rounded microglial cells on top of the layers (usually 10-14 d after the cortical culture is plated) and floating in the medium. You should harvest the floating microglia at this time.


141 Shake the flask with a back-and-forth motion several times to help lift the microglia off.

142 Shake additional ~  **00:10:00** placing the flask on a shaker at moderate speed.

10m


143 Collect the fluid in two 15 ml conical tube and centrifuge at  **3000 rpm, 00:05:00** .

Note that the speed is higher than for regular cells centrifugation, which is normally 1,000rpm.

144 Remove supernate and GENTLY resuspend the cell pellet in  **1 mL warm conditioned M10CG** with regular 1ml pipette tip.

145 

Count cells on hemacytometer to determine appropriate volume for cells density @ 60,000 cells/well on neuronal cultures or 150, 000 cells/well 96 well plates for biochemical assays.

146 When plating on neuronal cultures place plastic ring around the well and plate cells in the ring. Remove ring after  **02:00:00** . 2h

PLEASE BE SURE TO REMOVE NEURONAL MEDIA (SFIC) AND REPLACE WITH 2-3 DAY ASTROGLIAL CONDITIONED MEDIA BEFORE PLATING MICROGLIA.