

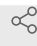


Sep 15, 2022

Immunoprecipitation protocol

Chuyu Chen¹¹Northwestern university

1 Works for me

 Sharedx.doi.org/10.17504/protocols.io.bp2l694zzlqe/v1 divya.darwinarulseeli

ABSTRACT

Immunoprecipitation procedure for brain tissue material

DOI

dx.doi.org/10.17504/protocols.io.bp2l694zzlqe/v1

PROTOCOL CITATION

Chuyu Chen 2022. Immunoprecipitation protocol. **protocols.io**
<https://protocols.io/view/immunoprecipitation-protocol-cgpqtmvw>



LICENSE

————— This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](https://creativecommons.org/licenses/by/4.0/), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

CREATED

Sep 15, 2022

LAST MODIFIED

Sep 15, 2022

PROTOCOL INTEGER ID

70096

- 1 Add 750ug tissue in 1ml IP buffer and antibody, incubate 4°C 2hrs on tube rotator
- 2 Prepare ProteinA/G beads 30ul/sample. Pool all beads in one tube to perform wash

- 3 Centrifuge 4°C 13000g 30sec, remove supernatant
- 4 Resuspend beads with 80ul PBS/sample. Centrifuge again at 13000g for 30sec.
- 5 Repeat wash 3 times.
- 6 Remove supernatant with fine tip
- 7 Resuspend beads with 30ul PBS/sample
- 8 Add 30ul beads to protein sample
- 9 Incubate 4°C overnight on tube rotator
- 10 Centrifuge 4°C 4000rpm 4min, remove supernatant
- 11 Wash with IP buffer (+ protease inhibitor) 800ul/sample
- 12 Repeat 3 times, total 4 wash
- 13 Remove supernatant with fine tip

14 Add 1x loading buffer 25ul and 5ul reducing agent. 75°C 5min

15 Centrifuge 4°C 4000rpm 4min, collect supernatant for western blotting