




JAN 16, 2024

Western Blot

Yiqin
Shen¹

¹University of California, San Diego

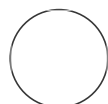

Yiqin Shen
University of California, San Diego

ABSTRACT

This is a protocol of Leo Parra-Rivas, PhD.

OPEN  ACCESS



Protocol Citation: Yiqin Shen 2024. Western Blot. protocols.io
<https://protocols.io/view/western-blot-c7fzzjp6>

License: This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](#), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

Protocol status: Working
We use this protocol and it's working

Created: Jan 11, 2024

Last Modified: Jan 16, 2024

PROTOCOL integer ID:
93401

Preparing Proteins

1 Dilute protein samples to $1.25 \mu\text{g}/\mu\text{L}$. For running a 10-well gel, each loading sample should be around $15 \mu\text{L}$ to $20 \mu\text{L}$.

2 Add 4x sample buffer to protein samples.

3 Denature protein at $95-100^\circ\text{C}$ for 00:05:00 to 00:10:00

15m

Running the Gel

4 Secure the gel (NuPAGE 4-12% Bis-Tris Gel; 1.5mm 10 wells) into the module. Fill the tank with run buffer. Load $15 \mu\text{L}$ of the ladder, and $20 \mu\text{L}$ of the samples (25ug/20ul, more if less concentrated).

Run at 70-80V for 01:15:00. Until the blue band reaches the end of the gel, run for 10 more min.

4.1 Ensure the strip on the bottom is removed. Remove the comb without breaking the wells.

4.2 Running Buffer: 1900ml H₂O, 100ml Running buffer 20X

4.3 You may add more buffer if the buffer level is low.
You may increase the voltage to 100V.

Transferring the Gel

5 Assemble the module according to instructions (membrane = 0.2ul pores; gel needs to be flipped). Ensure all bubbles are removed.
Add transfer buffer (ice cold) to the tank. Securely locate the module. Add ice around the tank before running.

Run at 25-30V for 01:25:00

5.1 transfer buffer: 1500ml H₂O, 400ml methanol, 100ml transfer buffer 20X

Staining

6 Fixation: rinse membrane with 1xPBS. Incubate in 0.4% PFA for 00:30:00 . 30m

7 Blocking: wash with 1xPBS, 3 times 00:10:00 . Block with 5-10% milk for 01:00:00 . 1h 10m

8 Primary antibody: Specific primary antibody diluted in 5% nonfat dried milk in TBST, 4 °C overnight 30m
with rocker.
Wash: Recover the primary antibody, and rinse with wash buffer. Wash for 00:30:00 with wash buffer

8.1 Wash buffer: 1800ml H₂O, 200ml 10X PBS, 2ml Tween-20

9 Secondary antibody: Secondary antibody 1:2000 in milk, 01:00:00 Room temperature . Add 1h
approximately 10 ml to each tray.
Wash: Rinse and wash with wash buffer. The membrane can be stored in the wash buffer

Imaging

10 Mix imaging reagent 1:1. Either apply reagent to membrane before imaging, or incubate the membrane with mixed reagent for 00:05:00 while rocking. 5m
Avoid bubbles when placing the membrane on the machine.
Single image, start with 0.1s exposure and gradually increase until desired exposure is achieved.
More reagent can be added if signal is reduced.

