




Jun 23, 2021

Modifications to the "HTAPP_TST- Nuclei isolation from frozen tissue V.2" Protocol

Xian Adiconis¹¹Broad Institute

1 Works for me

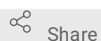

Share

This document is published without a DOI.

Human Cell Atlas Method Development Community

Xian Adiconis

ABSTRACT

A few modifications

DOCUMENT CITATION

Xian Adiconis 2021. Modifications to the "HTAPP_TST- Nuclei isolation from frozen tissue V.2" Protocol.
protocols.io
<https://protocols.io/view/modifications-to-the-34-htapp-tst-nuclei-isolation-bv2sn8ee>

LICENSE

This is an open access document distributed under the terms of the [Creative Commons Attribution License](#), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

CREATED

Jun 23, 2021

LAST MODIFIED

Jun 23, 2021

DOCUMENT INTEGER ID

50994

ABSTRACT

A few modifications

When working with larger amount tissues (i.e. combined kidney tissue from 3 mice), the following modifications are made to the original "HTAPP_TST- Nuclei isolation from frozen tissue V.2" protocol.

1) In step 2, used a razor blade (RNaseZapped and rinsed) instead of spring scissor to chop the tissues in a petri dish with 1 ml TST buffer. This is to increase the homogenization efficiency for the larger tissue quantity.

2) In step 5, the 3 ml TST buffer was first added to the petri dish to give the additional rinse then pass through the 40 µm filter.

3) After step 8, do another 2 washes by adding 3 ml 1x ST buffer and spinning at same condition each time.

4) Given there is certain buffer condition requirement for compatibility with the downstream applications (i.e. 10x SC RNA-seq assay), <https://kb.10xgenomics.com/hc/en-us/articles/115001937123-What-buffers-can-be-used-for-washing-and-cell-resuspension>, the final loading volume (if re-suspended in ST buffer) should be less than 6.17 µl, as this would give the final Mg concentration at max allowed 3mM (ST buffer contains 21 mM Mg and the H2O and sample volume for 10x V3.1 assay is 43.2 µl). If lower nuclei concentration (therefore higher volume is needed) is anticipated, the final nuclei suspension should be made with 1x PBS (containing 1% BSA, 0.2U/µl RNaseIn).