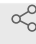




Oct 27, 2022

# Glycolysis stress test in organoids

[gustavo.parfitt](#)<sup>1</sup><sup>1</sup>Icahn School of Medicine1 *Works for me* Share

This protocol is published without a DOI.

**Ahfeldt Lab**

gustavo.parfitt

## ABSTRACT

Glycolysis stress test in organoids

## PROTOCOL CITATION

gustavo.parfitt 2022. Glycolysis stress test in organoids. **protocols.io**  
<https://protocols.io/view/glycolysis-stress-test-in-organoids-cig6ubze>



## LICENSE

————— This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](#), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

## CREATED


Oct 27, 2022

## LAST MODIFIED

Oct 27, 2022

## PROTOCOL INTEGER ID

71934

- 1  **Seahorse XFe96 Spheroid Microplates** **Agilent Technologies Catalog ##102905-100**


were coated with Laminin solution in L15+Bicarbonate: final concentration of 10 µg/mL .


- 2 Day 35 midbrain organoids were placed in the plates and media was changed every other day.

3 On the day of the assay

 [Agilent Seahorse XF Base Medium](#) **Agilent**

**Technologies Catalog #103334-100**

supplemented with B27 + N2 was warmed at  **37 °C** and the pH as adjusted to 7.4. And 175 ul was added to the organoid plates.

4 The plate containing the organoids were warmed at  **37 °C** in a non-CO2 incubator for 30 min

5 Using the

 [Seahorse XFp Glycolysis Stress Test](#) **Agilent**

**Technologies Catalog #103017-100**

kit the cartridges were loaded with loaded with Glucose (10mM), Oligomycin (1uM) and 2-DG (50mM).

6 The run was done in

Seahorse XFe96 Analyzer

Agilent BLE2100298 

7 After the run the organoids were lysed in RIPA buffer and total protein was quantified using BCA kit.

8 The total protein was used to normalize the ECAR and OCR.