



VERSION 1

MAY 13, 2023

OPEN  ACCESS

DOI:

dx.doi.org/10.17504/protocols.io.8epv59mr6g1b/v1

Protocol Citation: sandra.crameri 2023. NCEM Drop - Inactivation of sample on grid (TM-014). **protocols.io** <https://dx.doi.org/10.17504/protocols.io.8epv59mr6g1b/v1>
Version created by [sandra.crameri](#)

License: This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](#), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

Protocol status: Working
We use this protocol and it's working

Created: May 13, 2023


Last Modified: May 13, 2023

PROTOCOL integer ID:
81819

NCEM Drop - Inactivation of sample on grid (TM-014) V.1

sandra.crameri¹

¹CSIRO

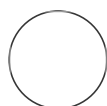






sandra.crameri

ABSTRACT

NCEM Drop - inactivation of samples on grid

GUIDELINES

This procedure is used for inactivation of 3Z agents, useful if the virus has low titre.

- 1 Adsorb  10 μ L sample to grid  00:10:00 10m
- 2 Immerse grid in drop of  2.5 % volume 10m
 Glutaraldehyde aqueous solution 25% EM Grade ampoules ProSciTech Catalog #C002 in
 0.1 Molarity (M) Sorenson's phosphate buffer  00:10:00
- 3 Wash  0.1 Molarity (M) Sorenson's phosphate buffer
- 4 Using filter paper wick excess buffer from grid - do not dry out, leave grid surface wet.
- 5 Stain  nano-W Contributed by users Catalog ##2018-5ML  00:01:00 1m
- 6 Using filter paper gradually wick stain dry from grid