

Version 2

Nov 21, 2020


PBMC- 01a - Isolation of Human PBMC from Buffy Coat

V.2

Marco Cosentino¹, Elisa Storelli¹, Alessandra Luini¹, Massimiliano LM Legnaro¹, Emanuela Rasini¹, Marco Ferrari¹, Franca Marino¹

¹Center for Research in Medical Pharmacology, University of Insubria (Varese, Italy)

1 Works for me dx.doi.org/10.17504/protocols.io.bpxjmpkn

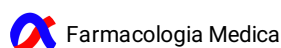

ABSTRACT

List of published work using this protocol

- Kustrimovic, N., Comi, C., Magistrelli, L., Rasini, E., Legnaro, M., Bombelli, R., Aleksic, I., Blandini, F., Minafra, B., Riboldazzi, G., Sturchio, A., Mauri, M., Bono, G., Marino, F., & Cosentino, M. (2018). Parkinson's disease patients have a complex phenotypic and functional Th1 bias: cross-sectional studies of CD4+ Th1/Th2/T17 and Treg in drug-naïve and drug-treated patients. *Journal of neuroinflammation*, 15(1), 205. <https://doi.org/10.1186/s12974-018-1248-8>
- Kustrimovic, N., Rasini, E., Legnaro, M., Bombelli, R., Aleksic, I., Blandini, F., Comi, C., Mauri, M., Minafra, B., Riboldazzi, G., Sanchez-Guajardo, V., Marino, F., & Cosentino, M. (2016). Dopaminergic Receptors on CD4+ T Naive and Memory Lymphocytes Correlate with Motor Impairment in Patients with Parkinson's Disease. *Scientific reports*, 6, 33738. <https://doi.org/10.1038/srep33738>
- Cosentino M., Ferrari M., Kustrimovic N., Rasini E., Marino F. (2015). Influence of dopamine receptor gene polymorphisms on circulating T lymphocytes: A pilot study in healthy subjects. *Human immunology*, 76, 10, 747-752. <https://doi.org/10.1016/j.humimm.2015.09.032>

DOI

dx.doi.org/10.17504/protocols.io.bpxjmpkn

PROTOCOL CITATION

Marco Cosentino, Elisa Storelli, Alessandra Luini, Massimiliano LM Legnaro, Emanuela Rasini, Marco Ferrari, Franca Marino 2020. PBMC- 01a - Isolation of Human PBMC from Buffy Coat. **protocols.io**
<https://dx.doi.org/10.17504/protocols.io.bpxjmpkn>
Version created by Farmacologia Medica

KEYWORDS

PBMC, Buffy Coat, Neuroimmune-Pharmacology, Parkinson's Disease, Cell isolation, Primary cell culture

LICENSE

— This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](https://creativecommons.org/licenses/by/4.0/), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

CREATED

Nov 21, 2020

LAST MODIFIED

Nov 21, 2020

PROTOCOL INTEGER ID

44747

MATERIALS TEXT

MATERIALS

 Fetal bovine serum

(FBS) BioWest Catalog #S181B-500

 Ficoll Paque PLUS Ge

Healthcare Catalog #17144003-500 ml

 RPMI

1640 EuroClone Catalog #ECM 0495L- 500 ml

 Trypan Blue Solution 0.4% Thermo Fisher

Scientific Catalog #15250061

Instrumentation required:

- Laminar flow hood
- Optical Microscope (manual cell count)

ABSTRACT

List of published work using this protocol

- Kustrimovic, N., Comi, C., Magistrelli, L., Rasini, E., Legnaro, M., Bombelli, R., Aleksic, I., Blandini, F., Minafra, B., Riboldazzi, G., Sturchio, A., Mauri, M., Bono, G., Marino, F., & Cosentino, M. (2018). Parkinson's disease patients have a complex phenotypic and functional Th1 bias: cross-sectional studies of CD4+ Th1/Th2/T17 and Treg in drug-naïve and drug-treated patients. *Journal of neuroinflammation*, 15(1), 205. <https://doi.org/10.1186/s12974-018-1248-8>
- Kustrimovic, N., Rasini, E., Legnaro, M., Bombelli, R., Aleksic, I., Blandini, F., Comi, C., Mauri, M., Minafra, B., Riboldazzi, G., Sanchez-Guajardo, V., Marino, F., & Cosentino, M. (2016). Dopaminergic Receptors on CD4+ T Naive and Memory Lymphocytes Correlate with Motor Impairment in Patients with Parkinson's Disease. *Scientific reports*, 6, 33738. <https://doi.org/10.1038/srep33738>
- Cosentino M., Ferrari M., Kustrimovic N., Rasini E., Marino F. (2015). Influence of dopamine receptor gene polymorphisms on circulating T lymphocytes: A pilot study in healthy subjects. *Human immunology*, 76, 10, 747-752. <https://doi.org/10.1016/j.humimm.2015.09.032>


BEFORE STARTING

If you need to obtain **PBMC for cell culture**, make sure you are using **sterile PBS, culture medium, filtered Lysis Buffer** and **sterile plastic disposables** as well. Moreover, **work under laminar flow hood when you are processing samples**. Otherwise, use non-sterile solutions and plastic disposables, and process samples in cell isolation laboratory.

ALL REAGENTS USED IN THIS PROTOCOL MUST BE AT ROOM TEMPERATURE!


- 1 Put the needed amount of blood sample from buffy coat into a 50 ml conical tube.


- 2 Add an equal volume of **PBS 1X** and mix well.

 SOLUTION- 02 - Phosphate Buffered Saline (PBS)
by Farmacologia Medica

- 3 Place  **3 mL** of **FICOLL** in a 15 mL conical tube.




CAREFULLY layer  **12 mL** of diluted blood on the FICOLL with a glass Pasteur Pipette to a final volume of 15 ml as shown in the figure below.

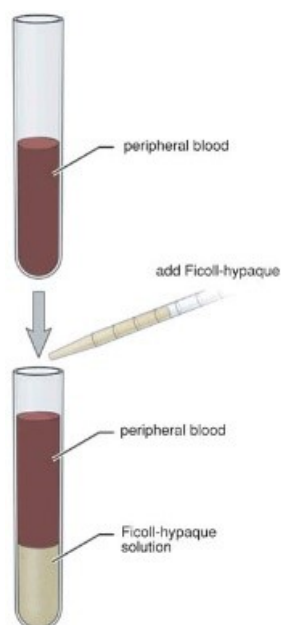

- 5 Centrifuge samples  **400 x g, 00:40:00** without break.

Allegra AVANTI 30
Centrifuge
Beckman Coulter Beckman Italy

- 6 After centrifugation, take out the tubes carefully to not disturb the mononuclear cell layer that appears as a white, cloudy band between the plasma and FICOLL as shown in the figure below.

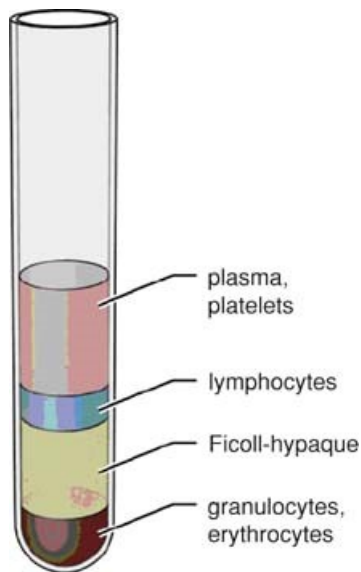

7

Carefully with a glass Pasteur pipette transfer mononuclear lymphocyte cell layer to another 15 ml conical tube.

- 8 Wash the isolated PBMC with **PBS/FBS 2%** to a final volume of **10 mL** and centrifuge at **300 x g, 00:10:00** at RT.



SOLUTION- 05 - Wash solution (PBS/FBS) for PBMC

by Elisa Storelli,

Center for Research in Medical Pharmacology, University of Insubria

Allegra AVANTI 30

Centrifuge

Beckman Coulter

Beckman Italy

- 9 Remove supernatants, resuspend pellet in **1 mL** of **Lysis Buffer** and add another **9 mL** of **Lysis Buffer**.
Immediately centrifuge the tubes at **100 x g, 00:10:00** at RT.



SOLUTION- 06 - Lysis Buffer

by **Elisa Storelli**,

Center for Research in Medical Pharmacology, University of Insubria

Allegra AVANTI 30

Centrifuge

Beckman Coulter

Beckman Italy

- 10 Remove supernatant and resuspend pellet in  **10 mL PBS/FBS 2%** and centrifuge at  **300 x g, 00:10:00** at RT.



SOLUTION- 05 - Wash solution (PBS/FBS) for PBMC

by **Elisa Storelli**,

Center for Research in Medical Pharmacology, University of Insubria

Allegra AVANTI 30

Centrifuge

Beckman Coulter

Beckman Italy

- 11 Remove supernatant and resuspend the obtained pellet in  **10 mL** of **RPMI/FBS 10%** for cell counting.



SOLUTION- 04 - Wash solution (RPMI/FBS) for PBMC

by **Farmacologia Medica**

- 12 **For manual cell count use Türk solution for checking purity.**

Follow protocol **CELL COUNT- 02**



SOLUTION- 08 - Türk solution
by Farmacologia Medica

13



OPTIONAL STEP

For automatic cell count with Cellometer machine use Trypan Blue.

Follow protocol **CELL COUNT- 03**.

Cellometer Auto T4
Automated cell counter
Nexcelom Bioscience EuroClone



SOLUTION- 09 - Trypan Blue solution
by Farmacologia Medica

14



If needed, check the purity of PBMC suspension by using morphological parameter of the flow cytometer.

For this test 0.5×10^6 PBMC in 500 μ l of PBS are enough.

BD FACS Celesta
Flow Cytometer
Becton Dickinson Milan Italy BD

15

Expected results



VIABILITY - The expected viability by Trypan Blue should be ≥ 90 %.

PURITY - The PBMC suspension obtained should contain at least 80% of lymphocytes, 10-15% of monocytes and few contaminant PMN cells ($\leq 5\%$) as confirmed by flow cytometry.

YIELD - The expected amount of PBMCs should be $\pm 100 \times 10^6$ starting from 25 ml of buffy coat.