




JAN 16, 2024

Lentivirus Production

Yiqin
Shen¹

¹University of California, San Diego

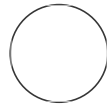

Yiqin Shen
University of California, San Diego

ABSTRACT

A protocol to produce lentivirus through cell transfection

OPEN  ACCESS



Protocol Citation: Yiqin Shen 2024. Lentivirus Production. **protocols.io** <https://protocols.io/view/lentivirus-production-c7nqzmdw>

License: This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](#), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited




Protocol status: Working
We use this protocol and it's working

Created: Jan 16, 2024

Last Modified: Jan 16, 2024

PROTOCOL integer ID:
93616



Day 1


- 1 Coat 15ml cell culture plates with D-poly-lysine (50ug/ml) for  00:30:00 to  01:00:00 .
Wash plates with cell culture water for 1 to 3 times. Let the plates to air dry.
- 2 Plate pre-grown low passage HEK cells at 6×10^6 cells per plate.
Add 20ml of media.
Incubate overnight at  37 °C 5% CO₂.

Note

Cell culture media: DMEM/F-12, GlutaMAX™ supplement + 10% FBS + 1% P/S

Day 2


- 3 Change media about 3 hours prior to transfection.
- 4 Transfect cells with viral components using CaCl₂ transfection Kit in the afternoon.
 - 4.1 Each plate = 20/25ug desired plasmid, 10.2ul PMD26, 15.6ul PsPAX, 93 ul CaCl₂; add water to 750 ul; 750 ul HBSS
 - 4.2 Add HBSS to tube 2; add plasmid, CaCl₂, H₂O to tube 1, and mix by vortex.
 - 4.3 Add tube 1 to tube 2 (on medium speed vortex) drop by drop (to form droplets).
 - 4.4 Wait  00:30:00 ; vortex again before use. 

Add  1.5 mL to each plate, distribute evenly across the plate.

Rock the plates gently to allow more distribution.

- 5 Incubate overnight at  37 °C 5% CO₂.

Day 3

- 6 Change media for all of the plates. Add  20 mL fresh media.

Day 4

- 7 Collect media from the plates (~24 hours after media change). Store at 4 degrees.

- 8 Add  20 mL of fresh media.




Day 5 & 6

5h 25m


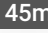

- 9 Repeat collection as day 4.

- 10 Spin down collected media at  500 x g, 4°C, 00:10:00 .

10m

- 11 Add 36ml of supernatant from step 10 to 12 ml of Lenti-X concentrator. Incubate the mixture at  4 °C in a rocker, for  04:00:00 to  Overnight .

4h 30m

- 12** Centrifuge the mixture from step 11 at  1500 x g, 4°C, 00:45:00 . After centrifugation, an off-white  pellet will be visible.
- 13** Remove supernatant from the tubes.
Resuspend the pellet in 100ul HBSS per pellet.
Aliquote 200ul per freezing tube.
Place the tubes on ice or freezing blocks immediately when you are finished, and transfer to  -80 °C for long term storage.