



OCT 18, 2023


🌐 Heart rate

Núria
Peñuelas¹

¹Vall d'Hebron Research Institute

Vilalab Public

Vilalab Public

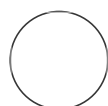

Núria Peñuelas
Vall d'Hebron Research Institute

OPEN  ACCESS



ABSTRACT

Heart rate measurement in mice.

DOI:
dx.doi.org/10.17504/protocols.io.kxygx355og8j/v1

Protocol Citation: Núria Peñuelas 2023. Heart rate.
protocols.io
<https://dx.doi.org/10.17504/protocols.io.kxygx355og8j/v1>

License: This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](#), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

Protocol status: Working
We use this protocol and it's working

Created: Oct 18, 2023

Last Modified: Oct 18, 2023

- 1 Place the animal in mild anesthesia with isoflurane at 1.5% and 1l/min O₂.
- 2 Prepare the PhysioSuite apparatus(Kent Scientific). Put the clip sensor in the right hind paw. Let it stabilize for 30 seconds before starting the recording.
- 3 When a stable heart rate is observed, perform 10 measurements.
- 4 Calculate the mean of the 10 measurements and represent the heart rate as average beats/minute.

represented as average beats/minute.