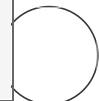




# BARseq - BARseq-styled in situ sequencing for barcoded rabies virus

Xiaoyin Chen<sup>1</sup>, Mara CP Rue<sup>1</sup>

<sup>1</sup>Allen Institute for Brain Science

 Xiaoyin Chen  
Allen Institute for Brain Science

AUG 16, 2023

## ABSTRACT

This protocol describes the application of BARseq-style in situ sequencing adapted for barcoded rabies virus. Similar procedures for both trans-synaptic tracing and retrograde tracing experiments.

## GUIDELINES

Standard precautions with RNA samples should be taken to reduce RNA degradation during tissue processing and library preparation. Pipetting and suctioning should be gentle throughout the whole procedure, and sample should not be left dried.

## MATERIALS

**DOI:**  
[dx.doi.org/10.17504/protocols.io.n2bvj82q5gk5/v1](https://dx.doi.org/10.17504/protocols.io.n2bvj82q5gk5/v1)

**Protocol Citation:** Xiaoyin Chen, Mara CP Rue 2023. BARseq - BARseq-styled in situ sequencing for barcoded rabies virus. **protocols.io** <https://dx.doi.org/10.17504/protocols.io.n2bvj82q5gk5/v1>

**MANUSCRIPT CITATION:**  
<https://doi.org/10.7554/eLife.87866.1>

**License:** This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](#), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

**Protocol status:** Working  
We use this protocol and it's working

### DNA oligos

	A	B	C
	YS220	GATCGTCGGACTGTAGAACTCTGAACCTGTCG	sequencing primer
	YS221	/5Alex594N/GATCGTCGGACTGTAGAACTCTGAACCTGTCG	Hybridization probe to visualize all genes
	XC1417	GTTCAGAGTTCTACAGTCCGACGATC	RCA primer for SNAP25 RNA padlock
	XC2757	/5AmMC12/NNNNNNNNNNNNNNNNNNNNNNNN	N20 for random priming
	XC2758	/5Alex488N/AGTCAGCGTCGAGCACGCGGCACTTATTGCA	Hybridization probe to visualize Slc17a7
	XC2759	/5Alex532N/TGAGTAGAGTTGACTAAGAGCCGTTAGATGCC	Hybridization probe to visualize Gad1
	XC2760	/5Alex647N/TCGCTGTACTAATAGTTGTCGACAGATCGTCA	Hybridization probe to visualize B19G
	XCAI5	/5phos/agctccggcattttgtattcaTCCTCTATGATTA CTGACTGCGTCTATTTAGTGGAGCCATTGCTATC TTCTTggatatacacaatccgtagattgct	pRV-4mCherry-NheI-N20 BC gapfilling padlock
	XCAI6	/5AmMC6/g+ac+at+at+tc+ga+gt+gactcataagaagt	primer 1 for pRV-4mCherry-NheI-N20 and pRV-4mCherry-CSCS2
	XCAI7	/5AmMC6/g+aa+gt+tg+aa+ta+ac+aaaatgccggagc	primer 2 for pRV-4mCherry-NheI-N20 and pRV-4mCherry-CSCS2

**Created:** Jan 12, 2023

**Last Modified:** Aug 16, 2023

**PROTOCOL integer ID:** 75215

**Keywords:** BARseq, in situ sequencing, spatial transcriptomics, barcoding-based neuroanatomy

A	B	C
XCAI63	/5AmMC6/gggagtgactgacacctccctccctg	RV B19G RT primer, Hyb XC2760 for detection
XCAI64	/5AmMC6/tatcaacatcaaggcagtcagggccc	RV B19G RT primer, Hyb XC2760 for detection
XCAI65	/5AmMC6/tcctgagacctgattgtgcacatcgg	RV B19G RT primer, Hyb XC2760 for detection
XCAI66	/5AmMC6/aactccatatgttgctggaggaggga	RV B19G RT primer, Hyb XC2760 for detection
XCAI67	/5AmMC6/tgcttttccaaacccagggaagtt	RV B19G RT primer, Hyb XC2760 for detection
XCAI68	/5AmMC6/ttcgtctgagcgaaagtcgtgcagg	RV B19G RT primer, Hyb XC2760 for detection
XCAI69	/5AmMC6/ttgcacgagacccatgttccatcca	RV B19G RT primer, Hyb XC2760 for detection
XCAI70	/5AmMC6/aggcctcttcatctacaaagccgca	RV B19G RT primer, Hyb XC2760 for detection
XCAI71	/5AmMC6/acatccctagtcctcgattctcgggc	RV B19G RT primer, Hyb XC2760 for detection
XCAI72	/5AmMC6/aagacaccgctactcctgagcacttc	RV B19G RT primer, Hyb XC2760 for detection
XCAI73	/5AmMC6/tggttttacagttcgaagccagcgg	RV B19G RT primer, Hyb XC2760 for detection
XCAI74	/5AmMC6/caccggccatcttcagttgtacgcg	RV B19G RT primer, Hyb XC2760 for detection
XCAI75	/5AmMC6/agtgtaggtttcagcctccgtcacao	RV B19G RT primer, Hyb XC2760 for detection
XCAI76	/5AmMC6/atgtaggagaacccctgacaggttggt	RV B19G RT primer, Hyb XC2760 for detection
XCAI77	/5phos/acacaatctcagaggacaggTCGCTGTACTA ATAGTTGTCGACAGATCGTCACTTCGTTCTcaat cgatcagaacctacgca	RV B19G padlock, Hyb XC2760 for detection, use with XC63
XCAI78	/5phos/gggaagtatgtattactgagtcTCGCTGTACTA ATAGTTGTCGACAGATCGTCACTTCGTTCTtgac ttgggtctcccgaactgg	RV B19G padlock, Hyb XC2760 for detection, use with XC64
XCAI79	/5phos/tgttgaagttcaccttcccgaTCGCTGTACTAAT AGTTGTCGACAGATCGTCACTTCGTTCTcggtga cgaggctgaggattt	RV B19G padlock, Hyb XC2760 for detection, use with XC65
XCAI80	/5phos/tcccagagatgcaatcatcccTCGCTGTACTAA TAGTTGTCGACAGATCGTCACTTCGTTCTgacct gacggcaatgtcttaa	RV B19G padlock, Hyb XC2760 for detection, use with XC66
XCAI81	/5phos/gtttcagacgtctcagtcattTCGCTGTACTAAT AGTTGTCGACAGATCGTCACTTCGTTCTcatgac aaccaagtcagtga	RV B19G padlock, Hyb XC2760 for detection, use with XC67
XCAI82	/5phos/cccgataagttggtgaacctgTCGCTGTACTAA TAGTTGTCGACAGATCGTCACTTCGTTCTaatga aaccaatggtgccct	RV B19G padlock, Hyb XC2760 for detection, use with XC68
XCAI83	/5phos/tggagttctaggacttagacttTCGCTGTACTAA TAGTTGTCGACAGATCGTCACTTCGTTCTgagca tgcaactcaagttatg	RV B19G padlock, Hyb XC2760 for detection, use with XC69
XCAI84	/5phos/ccaaaggagtgagacttgcgTCGCTGTACTA ATAGTTGTCGACAGATCGTCACTTCGTTCTatag tagaggggaagagagcat	RV B19G padlock, Hyb XC2760 for detection, use with XC70

A	B	C
XCAI85	/5phos/gattacaccatttggatgcccTCGCTGTAATAAT AGTTGTCGACAGATCGTCACTTCGTTCTacctact gctccactaaccac	RV B19G padlock, Hyb XC2760 for detection, use with XC71
XCAI86	/5phos/aggggtttccctagcgggaagtTCGCTGTAATA ATAGTTGTCGACAGATCGTCACTTCGTTCTtatg acagatcccttcaactcg	RV B19G padlock, Hyb XC2760 for detection, use with XC72
XCAI87	/5phos/caatccgtaccctgactaccgTCGCTGTAATAA TAGTTGTCGACAGATCGTCACTTCGTTCTcccag atatgaagagtctctaca	RV B19G padlock, Hyb XC2760 for detection, use with XC73
XCAI88	/5phos/ccagatgcatgtagagccgcgtTCGCTGTAATA ATAGTTGTCGACAGATCGTCACTTCGTTCTTaga aagcattccgccaaca	RV B19G padlock, Hyb XC2760 for detection, use with XC74
XCAI89	/5phos/gttcacttgcacagggcgttgcTCGCTGTAATAA AGTTGTCGACAGATCGTCACTTCGTTCTtcttagc cataaaagtgaacgg	RV B19G padlock, Hyb XC2760 for detection, use with XC75
XCAI90	/5phos/tggaggacgaaggatgcaccaTCGCTGTAATA ATAGTTGTCGACAGATCGTCACTTCGTTCTgctg cccaaacaatttgtag	RV B19G padlock, Hyb XC2760 for detection, use with XC76
XCAI13 1	TGGAGCCATTGCTATCTTCTTggatatacacaatccgt agattgct	sequencing primer for XCAI5

## Reagents



Phusion high-fidelity PCR kit Thermo Scientific Catalog  
#F553S



Tween-20 Sigma-aldrich Catalog #P-7949



BS(PEG)9, 100 mg (Note: BS(PEG)9 loses its effectiveness 1 month after reconstitution  
in DMSO. Prepare a fresh batch every month, especially if it has been frozen and  
thawed repeatedly. Thermo Scientific Catalog #21582



Formamide Thermo Fisher Scientific Catalog #AM9342



10x PBS Thermo Fisher Scientific Catalog  
#AM9624



RNase-free water Contributed by users













dNTP Mix (dATP, dCTP, dGTP, and dTTP, each at 10mM) Thermo Fisher  
Scientific Catalog #R0192



phi29 DNA Polymerase (10 U/μL) Thermo Scientific Catalog  
#EP0091



RNase H Enzymatics Catalog #Y9220L

	Glycerol Sigma Aldrich Catalog #G5516
	Ethanol Merck Millipore Catalog #100983
	Pierce™ MMTS (methyl methanethiosulfonate) Thermo Fisher Catalog #23011
	SSC (20X), RNase-free Thermo Fisher Catalog #AM9770
	RiboLock RNase Inhibitor (40 U/μL) Thermo Fisher Catalog #EO0381
	RevertAid H Minus Reverse Transcriptase (200 U/μL) Thermo Fisher Catalog #EP0452
	Paraformaldehyde 20% Electron Microscopy Sciences Catalog #15713
	BSA Molecular grade New England Biolabs Catalog #B9000S
	Ampligase DNA Ligase Kit Lucigen Catalog #A8101
	KCl (2 M) RNase-free Thermo Fisher Scientific Catalog #AM9640G
	Aminoallyl-dUTP Solution (50 mM) Thermo Fisher Scientific Catalog #R1101
	Tris (1 M) pH 8.0 RNase-free Thermo Fisher Scientific Catalog #AM9855G
	HiSeq SBS Kit v4 illumina Catalog #FC-401-4003
	Grace Bio-Labs HybriWell-FL™ sealing system Fluor-friendly adhesive chamber Sigma Aldrich Catalog #GBL612204

Other equipment required include incubators set at 37 °C, 45 °C, and 60 °C. All tubes should be RNase-free. RNase-free filter tips should be used. A Crest Xlight v3 spinning disk confocal on an Nikon Ti2E with Photometrics Kinetix, and Lumencor Celesta was used for imaging the sequencing steps. The filters and lasers used are indicated in Table 1.

A	B	C	D
Channels	Laser	Dichroic	Emission filter
G/YFP	514	Zt405/514/63 5rpc	FF01-565/24
T/RFP	561	FF421/491/5 67/659/776-Di01	FF01-441/511/593/684/817
A	640	Zt405/514/63 5rpc	FF01-676/29

A	B	C	D
C	640	Zt405/514/63 5rpc	FF01-775/140
GFP	488	FF421/491/5 72-Di01	69401m
DAPI	405	FF421/491/5 72-Di01	69401m
TexasRed	561	FF421/491/5 72-Di01	69401m
Cy5	640	Zt405/514/63 5rpc	ZET532/640m

Table 1. Laser and filter settings for sequencing imaging.

#### SAFETY WARNINGS



Use caution when handling liquids containing formaldehyde and formamide.

## Library preparation

- 1 Tissues with barcoded neurons should be cryo-sectioned to 20  $\mu\text{m}$  and mounted on slides. Slides can be stored at  $-80\text{ }^{\circ}\text{C}$  for up to a month.
- 2 **DAY 1**  
  
Take slide(s) out of  $-80\text{ }^{\circ}\text{C}$  and immerse immediately in 4% paraformaldehyde in 1x PBS (2 slides per 50mL falcon tube, back-to-back)
- 3 Incubate for 1 hour at room temperature on slow shaker
- 4 Wash the slides by immersing in 1x PBS (2 slides per 50ml falcon tube, back to back)

- 5 Wipe excess PBS off the surface of the chamber, then stick on the Hybriwell-FL chambers. Note that the ports on the chamber should be placed as far away from the tissue slices as possible.
- 6 Wash twice in PBST (1x PBS + 0.5% Tween-20)
- 7 Wash in 70% Ethanol for 5 mins
- 8 Wash in 85% Ethanol for 5 mins
- 9 Wash in 100% Ethanol for 5 mins
- 10 Replace with new 100% Ethanol, drop extra 100% Ethanol on top of slides and cover with ParaFilm to avoid evaporation. Incubate for at least 1.5 hrs at 4 °C (up to 3 hours)
- 11 Wash in PBST for 4-6 times, until all bubbles are cleared in the chamber and PBST flows in and out of the chamber smoothly.
- 12 Make reverse transcription mix: 50 µM N20 primer (XC2757), 2 µM XCAI6, 2 µM XCAI7, 20 U/µL RevertAid H Minus M-MuLV reverse transcriptase, 500 µM dNTP, 0.2 µg/µL BSA, 1 U/µL RiboLock RNase Inhibitor, 1x RevertAid RT buffer.

For the monosynaptic tracing experiments, the RT primers additionally included 2 µM of primers for the rabies glycoprotein (XCAI63 through XCAI76).

- 13** Incubate in reverse transcription mix overnight at 37 °C. Create a humidity chamber to avoid the slides drying out using kim-wipes and DI water.
- 14** **DAY 2:**
- Wash with PBST once
- 15** Incubate in a mixture of 1 µL BS(PEG)9 per 4 µL PBST (e.g. 200ul BS(PEG)9 and 800ul PBST) for one hour at room temperature
- 16** Wash with 1M Tris pH 8.0, then incubate in new 1M Tris pH 8.0 for 30 mins
- 17** Wash twice in PBST
- 18** Make non-gap-filling ligation mix: 1x Ampligase buffer, 20 nM padlock probe each, 0.5 U/µL Ampligase, 0.4 U/µL RNase H, 1 U/µL RiboLock RNase Inhibitor, 50 mM KCl (extra of those already provided by the ampligase buffer), 20% formamide.
- In the retrograde tracing experiments, the non-gap-filling padlock probe mix included all padlock probes for endogenous genes. In the monosynaptic tracing experiments, the non-gap-filling padlock probe mix included all padlock probes for endogenous genes except *Slc30a3*, and additionally included padlocks for the rabies glycoprotein (XCAI77 – XCAI90).
- 19** Incubate in ligation mix for at least 30 mins at 37 °C (can go longer but not shorter), then at least 45 mins at 45 °C (can go longer but not shorter).
- 20** Make the gap-filling ligation mix [same as the non-gap-filling mix with the rabies barcode padlock probe (XCAI5) as the only padlock probe, and with 50 µM dNTP, 0.2 U/µL Phusion DNA polymerase, and 5% glycerol]

1x Ampligase buffer, 20 nM padlock probe each, 0.5 U/μL Ampligase, 0.4 U/μL RNase H, 1 U/μL RiboLock RNase Inhibitor, 50 mM KCl (extra of those already provided by the ampligase buffer), 20% formamide, 50 μM dNTP, 0.2 U/μL Phusion DNA polymerase, and 5% glycerol.

- 21 Wash twice in PBST
- 22 Incubate in ligation mix for 5 mins at 37 °C, then 45 mins at 45 °C. **\*\*exact timing on this step!\*\***
- 23 Wash twice in PBST, then once in FISH Wash (2x SSC with 10% formamide)
- 24 Hybridize with 1 μM RCA primer (XC1417) in FISH wash for 10 mins at room temperature
- 25 Wash twice in FISH wash, then twice in PBST
- 26 Make RCA mix: 1 U/μL phi29 DNA polymerase, 1x phi29 polymerase buffer, 0.25 mM dNTP, 0.2 μg/μL BSA, 5% glycerol (extra of those from the enzymes), 125 μM aminoallyl dUTP
- 27 Incubate in RCA mix overnight at room temperature
- 28 **DAY 3:**  
  
Wash with PBST once



- 29 Incubate in a mixture of 1  $\mu$ L BS(PEG)9 per 4  $\mu$ L PBST (e.g. 200ul BS(PEG)9 and 800ul PBST) for one hour at room temperature
- 30 Wash with 1M Tris pH 8.0, then incubate in new 1M Tris pH 8.0 for 30 mins
- 31 Wash twice in PBST

## Sequencing

- 32 **Hybridization of Gene sequencing primer:**  
Wash with FISH wash (2x SSC with 10% formamide)
- 33 Hybridize sequencing primer (YS220) with a primer concentration of 1  $\mu$ M in FISH wash for 10 mins at room temperature
- 34 Wash with FISH wash three times, 2 mins each
- 35 Wash with PBST twice
- 36 **Sequence first cycle (genes and barcodes):**  
Do the following incubations. Unless noted with incubation temperature, each step is performed at room temperature. For steps without incubation time, treat these as quick washes. large flat

metal blocks can be used to place sample slides to quickly cycle through high and low temperatures. This version uses MiSeq Nano v2 kit:

<https://www.illumina.com/products/by-type/sequencing-kits/cluster-gen-sequencing-reagents/miseq-reagent-kit-v2.html>

- 36.1**      Incorporation Buffer 60 °C 3 mins x1
  
  
  
  
  
  
  
  
  
  
- 36.2**      2% PBST x1
  
  
  
  
  
  
  
  
  
  
- 36.3**      Idoacetamide blocker: For 9.3mg vial dilute pellet in between 2.5-3.5mL 2% PBST. Make fresh tube daily and store out of light.  
  
Idoacetamide blocker 60 °C 3 mins x1
  
  
  
  
  
  
  
  
  
  
- 36.4**      2% PBST x1
  
  
  
  
  
  
  
  
  
  
- 36.5**      Incorporation Buffer x2
  
  
  
  
  
  
  
  
  
  
- 36.6**      IRM 60 °C 3 mins x2
  
  
  
  
  
  
  
  
  
  
- 36.7**      2% PBST x1
  
  
  
  
  
  
  
  
  
  
- 36.8**      2% PBST 60 °C 3 mins x4

**36.9** Replace PBST with USM and Image \*\*if slides are dirty, clean with 70% Ethanol before adding USM\*\*

**37 Sequence subsequent cycles (genes and barcodes):**

Do the following incubations. Unless noted with incubation temperature, each step is performed at room temperature. For steps without incubation time, treat these as quick washes. large flat metal blocks can be used to place sample slides to quickly cycle through high and low temperatures.

**37.1** Incorporation buffer x2

**37.2** CRM 60 °C 3 mins x2

**37.3** Incorporation buffer x1 - Wipe ports after adding the incorporation buffer, to ensure that no CRM is left on the slide's surface

**37.4** 2% PBST x1

**37.5** Idoacetamide blocker: For 9.3mg vial dilute pellet in between 2.5-3.5mL 2% PBST. Make fresh tube daily and store out of light.

Idoacetamide blocker 60 °C 3 mins x1

**37.6** 2% PBST x1

- 37.7** Incorporation buffer x2
- 37.8** IRM 60 °C 3 mins x2
- 37.9** 2% PBST x1
- 37.10** 2% PBST 60 °C 3 mins x4
- 37.11** Replace 2% PBST with USM and image \*\*if slides are dirty, clean with 70% Ethanol before adding USM\*\*

**38 After completing all gene sequencing imaging cycles:**

**Hybridization cycle**

- 38.1 Hybridize probes:**  
Make strip buffer: 60% formamide 2xSSC 0.01% Tween20  
  
Strip buffer 60 °C 5 mins x3  
Cool down quickly on metal plates between washes, place on metal plates in 60 °C oven to heat up quickly.
- 38.2** FISH wash (2x SSC with 10% formamide) 1x

- 38.3** Hybridize probes (YS221, XC2758, XC2759, XC2760) with a primer concentration of 1  $\mu$ M in FISH wash at 60 °C for 2 minutes, then for 10 mins at room temperature. Rotate plates in holder to ensure they cool down slowly.
- 38.4** FISH wash x1
- 38.5** 0.002 mg/ML DAPI in PBST, room temperature for 5 mins
- 38.6** Replace PBST with USM and image \*\*if slides are dirty, clean with 70% Ethanol before adding USM\*\*
- 39** **After completing all gene sequencing imaging cycles and hybridization cycle:**
- Hybridize barcode sequencing primers**
- 39.1** Strip buffer (60% formamide 2xSSC 0.01% Tween20)  
60 °C 5 mins x3  
Cool down quickly on metal plates between washes, place on metal plates in 60 °C oven to heat up quickly.
- 39.2** FISH wash (2x SSC with 10% formamide) 1x
- 39.3** Hybridize probe XCA1131 with a primer concentration of 1  $\mu$ M in 2x SSC and FISH wash at 60 °C for 2 minutes, then for 10 mins at room temperature. Rotate plates in holder to ensure they cool down slowly.
- 39.4** FISH wash x1

**40** Sequencing barcodes: The sequencing procedures are the same as those for gene sequencing (Steps 36 and 37)

After hybridizing barcode primers, go directly back to step 36, sequence first barcode imaging cycle. Then repeat step 37 for desired length of barcode sequence.