



NOV 13, 2023

OPEN  ACCESS**DOI:**

[dx.doi.org/10.17504/protocols.io.dm6gp356pvzp/v1](https://dx.doi.org/10.17504/protocols.io.dm6gp356pvzp/v1)

**Protocol Citation:** Ching-Chieh Chou, Judith Frydman 2023. Immunohistochemistry on free-floating and paraffin-embedded tissue sections.

**protocols.io**

<https://dx.doi.org/10.17504/protocols.io.dm6gp356pvzp/v1>

**License:** This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](#), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited


**Protocol status:** Working  
We use this protocol and it's working

**Created:** Nov 13, 2023

## Immunohistochemistry on free-floating and paraffin-embedded tissue sections

Ching-Chieh Chou<sup>1</sup>, Judith Frydman<sup>1</sup>

<sup>1</sup>Department of Biology, Stanford University

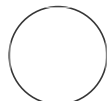

Ching-Chieh Chou  
Stanford University

### ABSTRACT

This protocol is used for free-floating frozen (30-50 microns) and paraffin-embedded (10 microns) tissue sections.

**Last Modified:** Nov 13,  
2023

**PROTOCOL integer ID:**  
90885

**Keywords:** ASAPCRN,  
Immunohistochemistry

**Funders**

**Acknowledgement:**

Aligning Science Across  
Parkinson's  
Grant ID: ASAP-000282

## Materials

### 1 Free-floating

- PBS
- 0.3% Triton-X-100 in PBS 1x, stored at 4°C
- Normal Donkey Serum (NDS), stored at -20°C
- Blocking Buffer: 10% NDS and 0.03% Triton-X-100 in PBS 1x. Diluted from 100% NDS (stored at -20°C) and 0.3% Triton-X-100. Filtered with 0.22 µM filter.
- Primary Antibody
- Secondary Antibody
- Hoechst dye
- Mounting medium

### Paraffin-embedded

- Deparaffinization station
- Slide rack
- Glass slide holder and tray
- Humidity and slide incubation chamber
- Xylene
- Ethanol
- Deionized water (diH<sub>2</sub>O)
- PBS
- Citrate buffer, pH 6.0, 10x
- 0.3% Triton-X-100 in PBS 1x, stored at 4°C
- Normal Donkey Serum (NDS), stored at -20°C
- Blocking Buffer: 10% NDS and 0.03% Triton-X-100 in PBS 1x. Diluted from 100% NDS (stored at -20°C) and 0.3% Triton-X-100. Filtered with 0.22 µM filter.
- Hydrophobic Barrier Pap Pen
- Primary Antibody
- Secondary Antibody
- Hoechst dye
- Mounting medium

## Procedure for free-floating tissue staining

- 2 Collect free-floating frozen tissue sections and process in a 24-well culture plate.
  - 2.1 **Day 1:** Wash the tissue for 3 X 5 min in PBS 1x.
  - 2.2 Permeabilize the tissue in 0.3% Triton-X-100 in PBS 1x for 20 min at room temperature on a rocking shaker.
  - 2.3 Wash the tissue for 3 X 5 min in PBS 1x.
  - 2.4 Block the tissue in Blocking Buffer for 1 hr at room temperature on the rocking shaker.
  - 2.5 Briefly rinse the tissue in PBS 1x once.
  - 2.6 Prepare Primary Antibody in Blocking Buffer.
  - 2.7 Incubate the tissue with Primary Antibody and transfer the plate to a refrigerator for overnight at 4°C.

- 2.8**      **Day 2:** Wash the tissue for 3 X 5 min in PBS 1x.
- 2.9**      Prepare Secondary Antibody in filtered 10% NDS in PBS 1x.
- 2.10**     Incubate the tissue with Secondary Antibody for 1 hr at room temperature on the rocking shaker in the dark.
- Note**

After Secondary Antibody incubation, all procedures should be performed in the dark.
- 2.11**     Wash the tissue for 2 X 5 min in PBS 1x.
- 2.12**     Incubate the tissue in Hoechst (1:2000) in PBS 1x for 10 min at room temperature on the rocking shaker.
- 2.13**     Rinse the tissue for 2 X 5 min in PBS 1x.
- 2.14**     Transfer the tissue to a slide and remove solution around the tissue by Kimwipes.
- 2.15**     Add mounting medium and put on a coverslip starting with top slowly to avoid generating any bubbles.

**2.16** Gently move the coverslip and let the mounting medium cover the edge of the slide. Let it dry overnight in the dark.

**2.17** **Day 3:** Acquire imaging by a microscope.

## Procedure for paraffin-embedded tissue staining

**3** Process paraffin-embedded tissue sections on glass slides.

**3.1** **Day 1:** Place the slides in a slide rack and rehydrate sections by sequential washes in a deparaffinization station:

1. Xylene: 2 x 5 min
2. 100% Ethanol: 2 x 10 min
3. 95% Ethanol: 10 min
4. 70% Ethanol: 10 min
5. 50% Ethanol: 10 min
6. diH<sub>2</sub>O: 2 x 5 min

**3.2** For antigen retrieval, dilute Citrate buffer in water to make 1x solution.

**3.3** Fill a beaker with Citrate buffer 1x and place the beaker on a hot plate covered with foil. Boil the buffer to 100°C.

**3.4** After the final step of rehydration, transfer the slides from the slide rack to a glass slide holder. Place the holder in the beaker and heat the slides for 20 min.

### Note

Frequently monitor the temperature of Citrate buffer by a thermometer.

- 3.5** Take off the beaker from the hot plate and let it cool for 20 min.
- 3.6** Add PBS 1x to a glass tray. Put the glass slide holder in the tray and wash for 3 X 5 min on the rocking shaker.
- 3.7** After PBS wash, carefully remove all the solution around the section with Kimwipes and draw a circle around the tissue with PAP pen a few times.
- 3.8** Add Blocking Buffer (~200  $\mu$ L) to cover the tissue without going over the PAP pen repellent border. Incubate for 2 hr at room temperature.
- 3.9** Prepare Primary Antibody in Blocking Buffer.
- 3.10** Remove all the Blocking Buffer around the tissue with Kimwipes.
- 3.11** Transfer the slides to a humidity and slide incubation chamber. Add water to bottom of incubation chamber.
- 3.12** Incubate the tissue with Primary Antibody (~200  $\mu$ L) and transfer the chamber to cold room for overnight at 4°C.

- 3.13**     **Day 2:** Knock the slides gently on Kimwipes and dry around the tissue.
- 3.14**     Stagger slides in the glase slide tray and wash with PBS 1x for 3 X 5 min on the rocking shaker.
- 3.15**     Prepare Secondary Antibody in 0.03% Triton-X-100 in PBS 1x.
- 3.16**     After PBS wash, remove all the solution around the tissue with Kimwipes.
- 3.17**     Transfer the slides to the humidity and slide incubation chamber. Add water to bottom of incubation chamber.
- 3.18**     Incubate the tissue with Secondary Antibody (~200 µL) for 1 hr at room temperature in the dark.
- Note**

After Secondary Antibody incubation, all procedures should be performed in the dark.
- 3.19**     Gently wipe off the solution around the tissue. Stagger slides in the glase slide tray and wash with PBS 1x for 2 X 5 min on the rocking shaker, covered with foil.
- 3.20**     Remove solution around the tissue. Transfer the slides to the humidity and slide incubation chamber.

- 3.21** Incubate with Hoechst dye at 1:2000 dilution in PBS 1x for 10 min.
- 3.22** Place the slides in the glase slide tray and rinse with PBS 1x for 2 X 5 min on the rocking shaker.
- 3.23** After the final wash, gently tap on the slides and remove all the solution around the tissue. Dry the slides overnight in the dark.
- 3.24** **Day 3:** Add mounting medium on the coverslip. Place the slides on the coverslip slowly to avoid generating any bubbles.
- 3.25** Flip the slides with coverslip on the top. Gently move the coverslip and let the mounting medium cover the edge of the slide.
- 3.26** Dry the slide overnight in the dark.
- 3.27** **Day 4:** Acquire imaging by a microscope.