



7

Jun 01, 2022

2-step PCR mixture and conditions (Barcoded-head primers for seqs pooling) V.7

Yin-Tse Huang¹

¹Kaohsiung Medical University



protocol .



Yin-Tse Huang

PCR mixture and condition (2X SUPERGREEN PCR MASTER MIX)

Yin-Tse Huang 2022. 2-step PCR mixture and conditions (Barcoded-head primers for seqs pooling).

protocols.io

<https://protocols.io/view/2-step-pcr-mixture-and-conditions-barcoded-head-pr-cadgsa3w>

Yin-Tse Huang

protocol

Herbold CW, Pelikan C, Kuzyk O, Hausmann B, Angel R, Berry D, Loy A. 2015. A flexible and economical barcoding approach for highly multiplexed amplicon sequencing of diverse target genes. Front. Microbiol. [Internet] 6:731. Available from: <http://dx.doi.org/10.3389/fmicb.2015.00731>

protocol ,

Jun 01, 2022

Jun 01, 2022

63624

- 1 Wear glove, clean up the working bench w. 1% bleach

For 1' PCR head-primers

- 2 Prepare 1' PCR master mixutre for **head-primers** (prepare 1.2X of solutions for pipetting error if needed)

PCR mixture for head-primers for each reaction

A	B	C	D
Component	Volume	Volume (1.2X)	Final conc.
Forward Primer (10 µM)	1.6 µl	1.9 µl	1 µM
Reverse Primer (10 µM)	1.6 µl	1.9 µl	1 µM
2X Supergreen PCR Master Mix	7.8 µl	9.4 µl	-
ddH2O	4.1 µl	4.9 µl	-
Total volume	15 µl	18 µl	-

Negative control ALWAYS NEEDED! For example, if you have 5 PCR reactions to run, prepare master mixture for 6 reactions (5 DNA template + 1 negative control).

- 3 Mix the 1' PCR master mixture gently by pippeting. Quick spin the tube.

- 4 Transfer **15 µL** 1' PCR master mixutre in 8-strip PCR tubes.

- 5

Add **0.6 µL** DNA template in 8-strip PCR tubes, resulting in a **15.6 µL** reaction mixture for 1' PCR.

Negative control contains only **15 µL** master mixture but not DNA template

6 Mix the reaction mixture gently by tapping the tubes. Quick spin the tubes.

7 Carry out PCR using the following condition:

1' PCR condition for **head-primers**

A	B	C	D
Step	Temp	Sec	Cycle
Initial denaturation	95 °C	30-180 (a)	
Denaturation	98 °C	15	20-25 cycles
Annealing	64-68 °C varied (b)	15	
Extension	72 °C	60-180 (c)	
Final extension	72 °C	210	
Preservation	Preservation	4 °C	∞

a. Varied depend on template complexity


b. Annealing varied, **62-65°C** is working based on test on 220530; Refer to 1' PCR primers for annealing temperature

c. 1kb ~ 1min extension; enough time allow full extension of sequence

7.1 1' hear-primers used in Huang lab


A	B	C	D
Name	Sequence	Tm°C	CG%
NS1B1ngs_H1	GCTATGCGCGAGCTGCcctngttgatctgccagt	71.7	60
ITS4ngs_H1	GCTATGCGCGAGCTGCtctctscgcttattgatatgc	69	55.6
LR5_H1	GCTATGCGCGAGCTGCtctgagggaaactcg	70.2	60.6
EF1-526F_H1	GCTATGCGCGAGCTGCgtctgtygtatygghcaygt	71	59.3
EF1-1567R_H1	GCTATGCGCGAGCTGCchgtccrctaccacratctt	70.6	56
EF1-2218R_H1	GCTATGCGCGAGCTGCatgacaccracrcrgrtytg	72.2	60.3
Ben2f_H1	GCTATGCGCGAGCTGCtccagactggtcagtggttaa	70.5	56.8
Bt2b_H1	GCTATGCGCGAGCTGCaccctcagtgtagtaccctggc	74.5	62.5
T22_H1	GCTATGCGCGAGCTGCtctggatgtttgggaatcc	70.3	56.8
RPB2-3bF_H1	GCTATGCGCGAGCTGCggwggwtaytytyatyaatgg	65.6	48.7
RPB2-7cR_H1	GCTATGCGCGAGCTGCcccatrgctgtttrccat	72.3	59.7
fRPB2-11aR_H1	GCTATGCGCGAGCTGCgctggatcttrtctcsacc	71.7	60.8

8 Carry out **electrophoresis** for inspection of DNA products

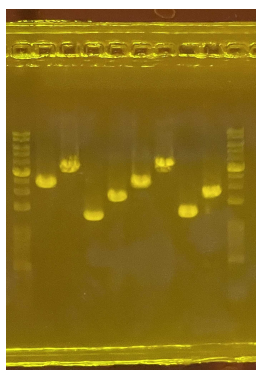

Gel before markdown

9 Markdown wells and upload the pictures to the Lab Google drive

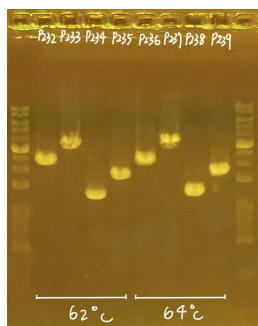

Marked gel picture go to the Lab Google drive

For 2' PCR barcoded-head primers

- 10 Prepare 2' PCR master mixutre for **barcoded-primers** (prepare 1.2X of solutions for pipetting error if needed)

PCR mixture for barcoded-primers for each reaction (**NO PRIMERS at this point!!**)

A	B	C	D
Component	Volume	Volume (1.2X)	Final conc.
2X Supergreen PCR Master Mix	10.75 µL	12.9 µL	-
ddH2O	10.75 µL	12.9 µL	-
Total volume	21.5 µL	25.8 µL	-

Negative control ALWAYS NEEDED! For example, if you have 5 PCR reactions to run, prepare master mixture for 6 reactions (5 DNA template + 1 negative control).

- 11 Mix the 2' PCR master mixture gently by pipetting. Quick spin the tube.
- 12 Transfer **21.5 µL** of the 2' PCR master mixture to 8-strip PCR tubes.
- 13 Add **2.5 µL pre-mixed barcoded-head primers** (Forward + Reverse) to each PCR tubes.

14

Add **1 µL** of **1' PCR product as template**, resulting in **25 µL** reaction mixture for 2' PCR.

Negative control contains only **24 µL** master mixture and premixed barcoded-head primers but not DNA template

- 15 Mix gently by tapping the tubes. Quick spin the tubes.

- 16 Carry out 2' PCR using the following condition:

2' PCR condition for **barcoded-head primers**


A	B	C	D
Step	Temp	Sec	Cycle
Initial denaturation	98 °C	30	
Denaturation	98 °C	15	10-15 cycles
Annealing	64-68 °C varied (a)	15	
Extension	72 °C	60 (b)	
Final extension	72 °C	210	
Preservation	Preservation	4 °C	∞

a. Annealing varied, **65 °C** is working based on test on 220531; Refer 2' PCR primers for annealing temperature
b. 1kb ~ 1min extension; enough time allow full extension of sequence

16.1 2' barcoded-head primers used in Huang lab


A	B	C	D	E	F	G	H	I	J
Name (1st)	Sequence	Tm°C	CG%	Name (2nd)	Sequence	Tm°C	CG%	Name (3rd)	Sequence
F1-1	aaatctcattttctgtagtgcgctatgctgagctgc	71.7	53.8	F2-1	ggccattctgtagtaagtagttcagctatgctgagctgc	70	53.8	F3-1	gggctagagttaaccttaattag
F1-2	aaccgcacgtctattttaagaagctatgctgagctgc	70.1	51.3	F2-2	cagatactaatctcgtcacatggctatgctgagctgc	68.5	51.3	F3-2	taatttatggctacagacgttag
F1-3	aactgtttatttcaaataatagctatgctgagctgc	61.3	35.9	F2-3	agagattgagcttatcgttttgcctatgctgagctgc	68.7	48.7	F3-3	gtaagcgtcatcacagcaaaact
F1-4	aatatgaagctccgacatatctggctatgctgagctgc	69.5	51.3	F2-4	ccattgattccagattatcatgtgctatgctgagctgc	67.7	48.7	F3-4	atccgcatactcaactcgaaagc
F1-5	acagtacaaacacgctcattaagctatgctgagctgc	70	51.3	F2-5	ccgtaggactctgtagtttaacgctatgctgagctgc	70.1	53.8	F3-5	gtgctgtcaaatcgttttgcctgc
F1-6	acagttgacctcgcgacattatgctatgctgagctgc	71.2	53.8	F2-6	gtcatgaccgctacatatctcagctatgctgagctgc	70	53.8	F3-6	gtgggataagcttgacatttag
F1-7	caccataaatgagattgctgagctatgctgagctgc	70.7	53.8	F2-7	cgataggacaagcaatgtactcagctatgctgagctgc	70.1	53.8	F3-7	gttgataggaggtcaaacaga
F1-8	accgtgtatgtatctgttaccgctatgctgagctgc	71	53.8	F2-8	cgatataatcgatccgcataccgctatgctgagctgc	71.4	56.4	F3-8	gttaaccataagtgcgacctat
R1-1	acgtcaatgctattccagctcaagctatgctgagctgc	71.4	53.8	R2-1	cagcgtttccaaagacattatggctatgctgagctgc	69.7	51.3	R3-1	taaaaaattcgtggaactccacag
R1-2	acgtcgaggtatcataaacttgcctatgctgagctgc	67.9	48.7	R2-2	cgctttatgtcttaagtagccgctatgctgagctgc	69.1	51.3	R3-2	taatccatgtgtctattcttagagc
R1-3	cgcaagtatgctcttctcatagctatgctgagctgc	70.1	53.8	R2-3	ctgctgagttataccacagtgacgctatgctgagctgc	71.5	56.4	R3-3	gtaaatctagggtgaaatgagtc
R1-4	agtcaatgctttgggtacataagctatgctgagctgc	70.1	51.3	R2-4	tacgttatattaactctagccgagctatgctgagctgc	67.4	48.7	R3-4	ctttctataataatccgggctaagc
R1-5	atatcttacacaaaagtatcgtatgctgagctgc	65.4	43.6	R2-5	tagaatgtcaacacaagtaggacgctatgctgagctgc	69.5	51.3	R3-5	ctttgttaagtgtgtccgttgtct
R1-6	cgactatttaacttcgcaacagctatgctgagctgc	69.3	51.3	R2-6	gaaatcgacaacaaattctgctctgctatgctgagctgc	69.7	51.3	R3-6	taggtatcattctcatcctatcggc
R1-7	atcgttggttctatgttcaggtatgctgagctgc	69.5	51.3	R2-7	gaatcatcaagaagggaacacagctatgctgagctgc	69.3	51.3	R3-7	tatagagacgggttcggtaaaa
R1-8	atgagatctatctagtagccgttgctatgctgagctgc	68.8	51.3	R2-8	gactatagtgaaaaatcacatagctatgctgagctgc	65.8	46.2	R3-8	tattagtagatagacactcggc
R1-9	gaggagccatagatgaatgagctatgctgagctgc	69.7	53.8	R2-9	atgacctacggagattacctgctatgctgagctgc	70.7	53.8	R3-9	tcagggttagaagactagttgta
R1-10	attgtagcattgaataggagcagctatgctgagctgc	68.7	48.7	R2-10	gagttctggatctatctgggtctgctatgctgagctgc	70	53.8	R3-10	tcgggttaaatgctaagcgtaat
R1-11	caaatctcatgctagctgacgctatgctgagctgc	70.2	53.8	R2-11	gcaacgaaacattcgttaagtagtctatgctgagctgc	68.1	48.7	R3-11	tgctgaaaaacaggaagtctcac
R1-12	caagaactagacgctgctcttaagctatgctgagctgc	70.2	53.8	R2-12	tgcccaatatgttagcaccctaagctatgctgagctgc	71.4	53.8	R3-12	accattcctaataagccaataggg

17 Carry out electrophoresis for inspection of DNA products

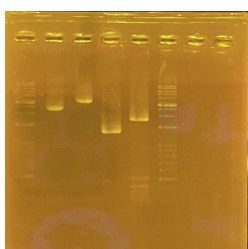

Gel before markdown

18 Markdown wells and upload the pictures to the Lab Google drive

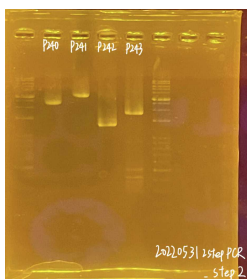

Marked gel picture go to the Lab Google drive

