



Jul 12, 2022


Primary microglial culture

Xiqun Chen¹, Qing Ye²

¹Department of Neurology Massachusetts General Hospital Harvard Medical School Charlestown MA 02129 USA, Aligning Science Across Parkinson's (ASAP) Collaborative Research Network Chevy Chase MD 20815 USA;

²Department of Neurology Massachusetts General Hospital Harvard Medical School Charlestown MA 02129 USA, Department of Neurology Longhua Hospital Shanghai University of Traditional Chinese Medicine Shanghai 200032 China

1 Works for me

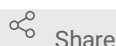

Share

dx.doi.org/10.17504/protocols.io.b4ubqwsn

Pranay Srivastava

ABSTRACT

Mixed glial cells were obtained from C57BL/6 mice embryonic day 17. Cells were cultured in high-glucose DMEM/F12 supplemented with 10% FBS in humidified air containing 5% CO₂ at 37°C. The culture medium was replaced with fresh medium 24 h after the initial preparation and every 3 days thereafter. After 1 week, microglial cells were obtained by mechanical shaking of the mixed glial cell cultures for 1 h. Cells were routinely monitored for purity by ionized calcium-binding adaptor molecule 1 (Iba1) staining and the population of Iba1+ cells was >95%

DOI

dx.doi.org/10.17504/protocols.io.b4ubqwsn

PROTOCOL CITATION

Xiqun Chen, Qing Ye 2022. Primary microglial culture. **protocols.io**
<https://dx.doi.org/10.17504/protocols.io.b4ubqwsn>



KEYWORDS

ASAPCRN

LICENSE

————— This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](https://creativecommons.org/licenses/by/4.0/), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

CREATED

Feb 08, 2022

LAST MODIFIED

Jul 12, 2022

PROTOCOL INTEGER ID

57955


- 1 Primary microglial culture - Use C57BL/6J mice at embryonic day 17
- 2 Anesthetized pregnant mice (1% sodium pentobarbital, 80mg/kg), dissect their embryos and collect the cortex.
(Separate and remove the soft membrane and blood vessels, rinse the cerebral cortex in PBS, and use the ophthalmic scissor to cut pieces of the cortex)
- 3 Collect the cortices in PBS in a 50 ml tube on ice
(The 50 ml tube contains 30 ml of PBS) ⚠ **On ice**
- 4 Transfer the cortices to 15 ml tubes containing 1.5 ml trypsin–EDTA (0.25%) and incubate it at ^{15m}
⚠ **37 °C** for ⌚ **00:15:00** Dissociate the cortices by triturating with a 10 mL serological pipette 10 – 15 times
- 5 Centrifuge the dissociated cortices (⌚ **400 x g** , ⌚ **00:05:00**). Aspirate the media and ^{5m}
resuspend the pellet in 5 ml

⌘ **DMEM, high glucose Thermo Fisher**

Scientific Catalog #11965092
- 6 Count the cells and plate them in a density of 50,000 cells/cm² into PLL coated T-75 flask.
Makeup the volume to 15 ml with the

⌘ **DMEM, high glucose Thermo Fisher**



Scientific Catalog #11965092
- 7 Change the culture medium the next day followed by the addition of a fresh culture medium ^{1w}
every ⌚ **168:00:00** till 14 days of culture
- 8 Put the flask on a shaker for ⌚ **06:00:00** at ⚠ **37 °C** (The microglia grow as a monolayer on ^{8h 5m}
the top)
Collect the detached cells, centrifuge (⌚ **400 x g** , **Room temperature** , ⌚ **00:05:00**) and

resuspended in  4 mL of

 DMEM, high glucose Thermo Fisher

Scientific Catalog #11965092

Count the cells and plate at a density of 50,000 cells/cm² in a 24 well plate


Incubate them at  37 °C in 5% CO₂. After  02:00:00 microglia usually attaches to the bottom of the plate.

The purity of the cells was monitored by staining with ionized calcium-binding adaptor molecule 1 (Iba1) antibody

(Microglia obtained by this method have a purity of 90%-95% and can be identified by immunofluorescence staining.

9 Transduction with BRAF (Optional)

1d

The purified microglial cells were seeded at a density of 80% and transduced with BRAFV^{600E}, or BRAF^{WT}, or vector lentivirus with 8 µg/ml polybrene (Sigma–Aldrich, USA) for  24:00:00 .

10 After transduction, the cells were cultured in


5d

 DMEM, high glucose Thermo Fisher

Scientific Catalog #11965092

with or

 Fetal Bovine Serum

without (FBS) ATCC Catalog #30-2020 for  120:00:00 , and the medium, as well as the cells, were collected for subsequent experiments.