

FEB 01, 2023

## OPEN ACCESS

**Protocol Citation:** Nadanai Laohakunakorn 2023. Photolithography for microfluidics. **protocols.io** <https://protocols.io/view/photolithography-for-microfluidics-b2vqqe5w>

**License:** This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](#), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

**Protocol status:** Working  
We use this protocol and it's working

**Created:** Dec 15, 2021


**Last Modified:** Feb 01, 2023

**PROTOCOL integer ID:**  
55952

## Photolithography for microfluidics

Nadanai Laohakunakorn<sup>1</sup>

<sup>1</sup>University of Edinburgh

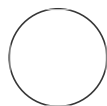

Nadanai Laohakunakorn  
University of Edinburgh

### ABSTRACT

Photolithography for microfluidic mold fabrication, successfully tested at the Scottish Microelectronics Centre. Produces rounded flow layer using positive resist SPR 220-7 with ~12μm features, and rectangular control layer using negative resist SU-8-3035 with ~35μm features).

Adapted from [dx.doi.org/10.17504/protocols.io.46wgzfe](https://doi.org/10.17504/protocols.io.46wgzfe) / Laohakunakorn et al. 2021 ([https://doi.org/10.1007/978-1-0716-1032-9\\_9](https://doi.org/10.1007/978-1-0716-1032-9_9)).

### MATERIALS

- Masks: from Compugraphics. Master-Si, NFR5009, grade 2160A2.

1X Master, Soda Lime, 5"x5"x0.09"

Minimum critical feature: 2.0μm

Defect density: 1 per in<sup>2</sup> (0.2 per cm<sup>2</sup>)

Registration: SPC

CD tolerance: 0.25μm

Defect size: 2μm

Data design: rectilinear

Write area: 1" smaller than mask size

Compugraphics product code: 2160A2

### Flow layer (SPR220, positive resist)

Digitised data = **Dark**

Data parity chrome up = **Wrong**

Title parity chrome up = **Wrong**

CD size = SPC

Data top cell = TOP

Data format = CIF

Data window = Compugraphics default

### Control layer (SU8, negative resist)

Digitised data = **Clear**

Data parity chrome up = **Wrong**

Title parity chrome up = **Wrong**

CD size = SPC  
Data top cell = TOP  
Data format = CIF  
Data window = Compugraphics default


**Double-check these settings very carefully!**

- Wafer properties: from Inseto.

Diameter: 100mm  
Material: Si  
Type: N  
Orientation: 100  
Dopant: P-doped  
Grade: Test?  
Lower resistivity: 1 ohm.cm  
Upper resistivity: 10 ohm.cm  
Thickness: 525 um  
Polish: SSP



## Flow Layer

### 1 HMDS priming

1.1 Prime a clean Si wafer in a sealed box with HMDS for  00:10:00

10m

### 2 Prebake

2.1 Transfer to hotplate and carry out pre-bake at  120 °C for  00:12:00 .

12m

### 3 Spin coat with SPR 220-7 (Megaposit) SPR\_220\_DATA\_SHEET\_RH.pdf

#### 3.1 Transfer wafer to spin coater and run following programme using SPR 220-7 resist:

2m 10s

1.  400 rpm, 00:00:30

2.  1000 rpm, 00:01:40

at 200 rpm/s. This coats the wafer with ~12  $\mu\text{m}$ .

### 4 Softbake



#### 4.1 Transfer to hotplate and carry out soft-bake at 105 °C for 00:06:00 . Remove promptly when done.

6m

### 5 Exposure

#### 5.1 Expose with the following settings (for a Karl Suss MA8 mask aligner):

1m 30s

Multiple exposure steps: [  00:01:20 exposure +  00:00:10 wait ] x 2 cycles, for a total of 160s (=638mJ/cm<sup>2</sup>)

WEC=cont

Expose type = prox

40 $\mu\text{m}$  alignment gap

5 $\mu\text{m}$  expose gap

WEC offset OFF

N2 purge NO

The Karl Suss machine has a flux of 3.99 mW/cm<sup>2</sup> as measured from last maintenance. We would like an exposure of 660 mJ/cm<sup>2</sup>.

#### 5.2 Wait: leave the wafer in a sealed wafer box for 02:00:00 .


2h

## 6 Post-exposure bake




6.1 Transfer to a hotplate and hold at  110 °C for  00:05:00 5m


6.2 Wait: leave the wafer in a sealed wafer box for  00:45:00 45m

## 7 Develop with MF26A


7.1 Develop right side up for up to  00:05:00, rinse with DI water, and dry gently with compressed N2. 5m

## 8 Reflow

8.1 Place the wafer on a hotplate at room temperature, then ramp up to  170 °C (at a ramp rate of ~0.5°C/s). Hold at  170 °C for  00:40:00, then turn off hotplate and allow to ramp down to room temperature (takes a few hours). 40m

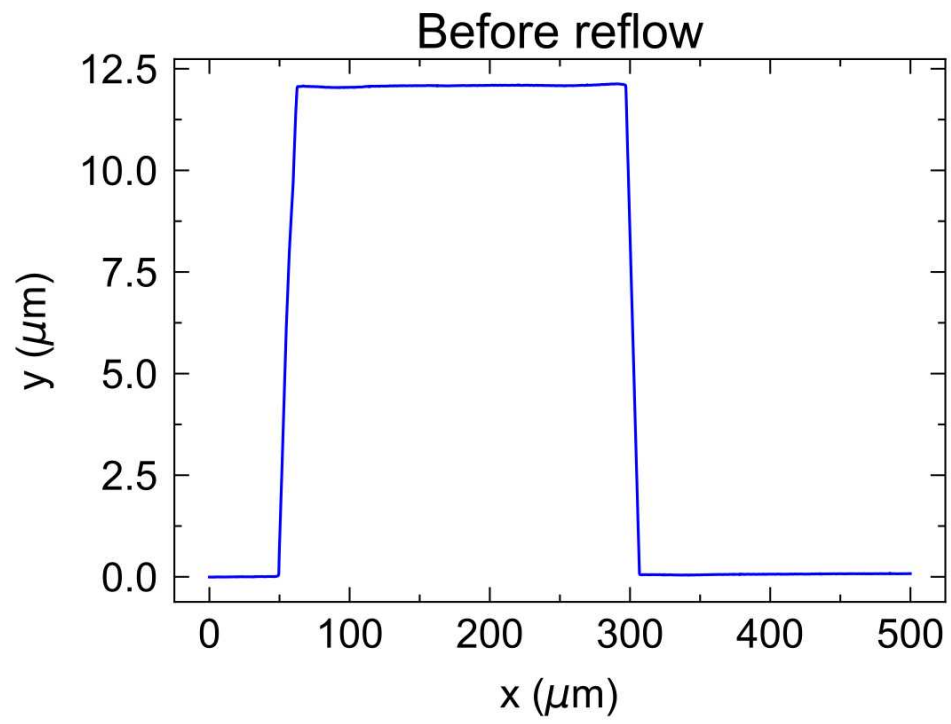

Before reflow

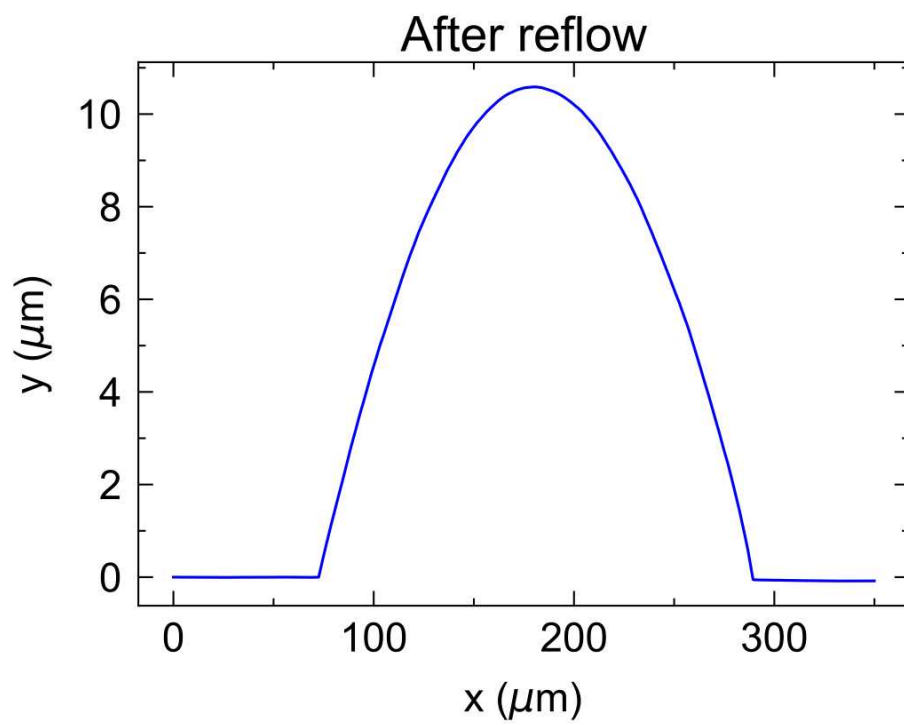






After reflow

## 9 O2 plasma treatment

- 9.1 Treat with O2 plasma for 10 minutes in barrel asher. (At SMC, Electrotech 508 with typical parameters: forward power 350W, flow 32%, pressure 0.8 torr)



## 10 Spin coat with SU-8-3035 (Kayaku) KAM-SU-8-3000-Datasheet-7.10-final.pdf

- 10.1 Spin coat with SU8-3035 using the following programme with ramp rate of 100rpm/s for all steps: 1m 51s

1.  500 rpm, 00:00:10
2.  3000 rpm, 00:01:05 = 25s ramp + 40s hold
3.  4000 rpm, 00:00:01
3.  3000 rpm, 00:00:05
3.  0 rpm, 00:00:30

Make sure the spin coater has a plastic liner installed. If necessary clean back and edge of wafer with wipe soaked in PGMEA.



## 11 Softbake

- 11.1 Place on hotplate for  00:12:30 at  95 °C 12m 30s
- Let wafer come to room temp (few mins) before exposure. If there are wrinkles, place on hotplate until wrinkles disappear.

## 12 Exposure

## 12.1 Expose with the following settings (for a Karl Suss MA8 mask aligner):

28s

Multiple exposure steps: [  00:00:18 exposure +  00:00:10 wait] x 5 cycles, for a total of 90s (=359mJ/cm<sup>2</sup>)

WEC=cont

Expose type = soft



30um alignment gap

5um expose gap



WEC offset OFF


N2 purge NO

## 13 Post-exposure bake

13.1 a. Place on hotplate at  65 °C for  00:01:00 .



5m

b. Move to second hotplate and hold at  95 °C for  00:04:00 .


14 **Wait** for a minimum of  00:10:00 .

10m


## 15 Develop

15.1 Develop upside down in PGMEA for  00:02:00 , and exchange with new PGMEA for a further  00:02:00 . Check for underdevelopment each time on the microscope (white film). If necessary continue development for up to 6 min.

4m

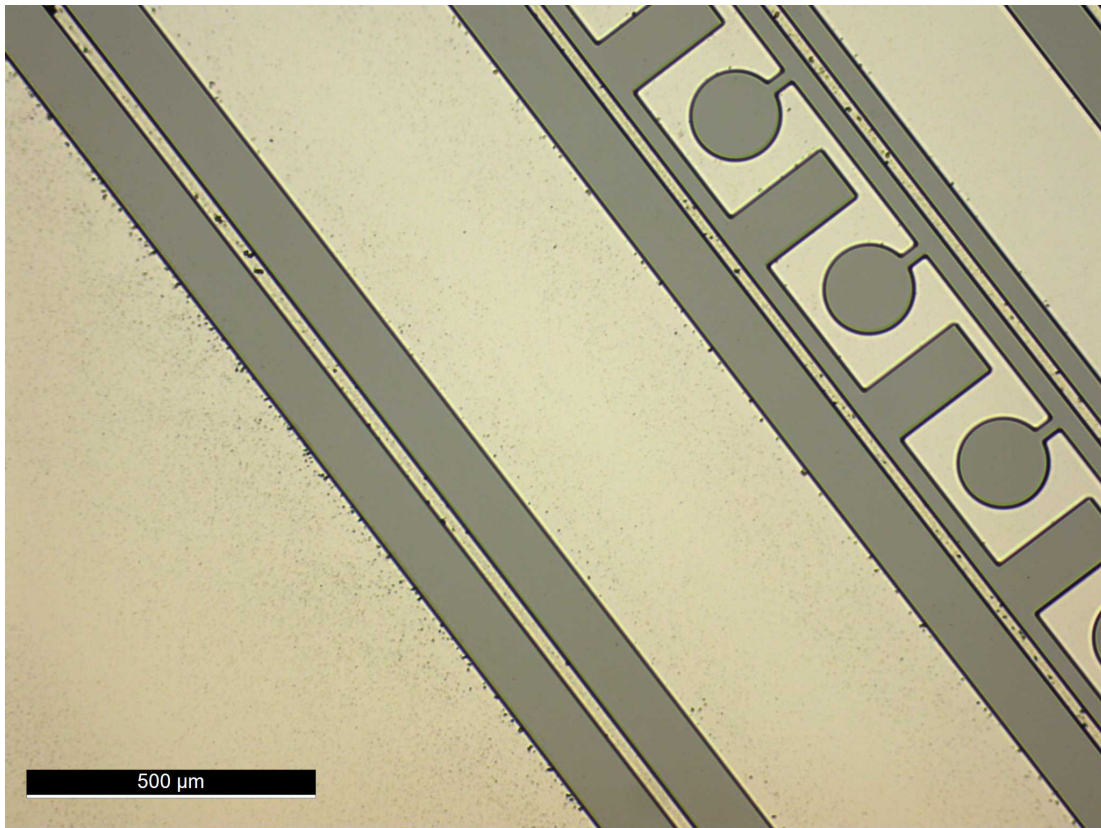

Underdevelopment after 2+1 min

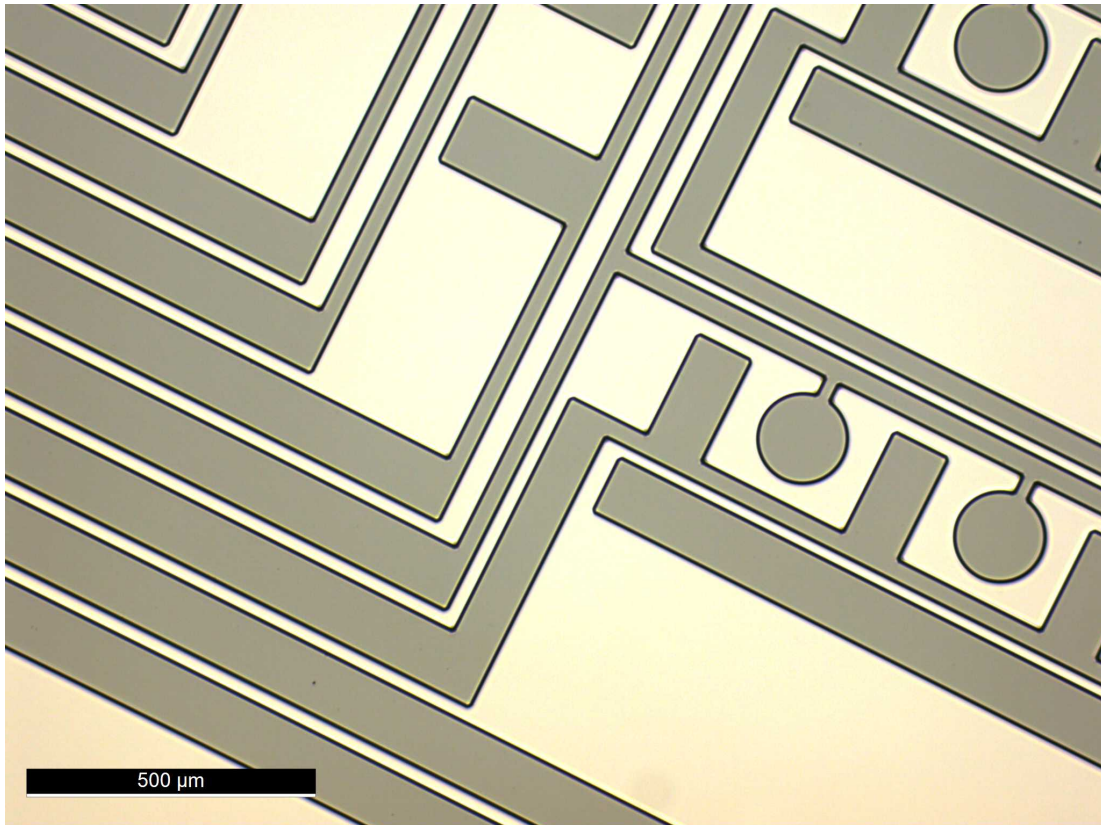




Correct development after 2+2 min




**15.2** Rinse with IPA and dry gently with compressed air.

## **16 Hardbake**

**16.1** Ramp to  135 °C , and hold at  135 °C for  02:00:00 .

2h

## **Silanization of wafers**

**17** Silanization must be carried out before first use, to prevent PDMS sticking to wafers. Leave  500  $\mu$ L of chlorotrimethylsilane (Sigma 386529) in an upturned Falcon cap, inside a sealed box with the wafers until the solvent has completely evaporated.