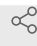





Sep 15, 2022

Fractionation of synaptosomes

Chuyu Chen¹, Ciarra Smith¹, Loukia Parisiadou¹¹Northwestern university1 *Works for me* Sharedx.doi.org/10.17504/protocols.io.x54v9d85zg3e/v1 divya.darwinarulseeli

ABSTRACT

This protocol details a step by step method to prepare pure fractions of synaptosomes for biochemical analysis.

DOI

dx.doi.org/10.17504/protocols.io.x54v9d85zg3e/v1

PROTOCOL CITATION

Chuyu Chen, Ciarra Smith, Loukia Parisiadou 2022. Fractionation of synaptosomes. **protocols.io**

<https://protocols.io/view/fractionation-of-synaptosomes-cgpntvme>



LICENSE

————— This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](https://creativecommons.org/licenses/by/4.0/), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

CREATED

Sep 15, 2022

LAST MODIFIED

Sep 15, 2022

PROTOCOL INTEGER ID

70094

- 1 Pre-chill the homogeniser and buffers on ice. Weigh the tissue
- 2 Add 4x volume of ice-cold homogenizing buffer to glass homogenizers

- 3 Apply 12 strokes of even tension with the homogeniser rod
- 4 Transfer to 1.5ml tube (original tube)
- 5 Determine which tissue has the smallest volume and use that volume for all samples
- 6 Centrifuge at 1,000g for 10mins
- 7 Transfer the supernatant into fresh 1.5ml tube and spin at 12,500g for 15mins
- 8 Discard supernatant carefully
- 9 Resuspend pellet in 1ml of H-buffer, transfer to fresh glass homogenizer and mash for 6 strokes even tension
- 10 In polypropylene/ultra-clear tube add 5ml of 1.2M sucrose
- 11 Slowly and steadily add in 5ml of 0.8M sucrose to create clear/distinct gradient
- 12 Overlay sample on top of gradient

- 13 Weigh sucrose gradient tubes (add H-buffer if weight is off by more than 0.05g)
- 14 Centrifuge with slow acceleration and zero brake in 4 degrees at 23,600g for 70 mins in a SW41 rotor
- 15 Using syringe, extract synaptosome layer between 0.8M and 1.2M sucrose and transfer into fresh 1.5ml tubes
- 16 Plate sample with appropriate dilutions into 24-well plate (dilute with H-buffer or PBS)
- 17 Spin in 4°C, max speed ~4,000 rpm for 20mins
- 18 Remove liquid and fix in 4% PFA for 20mins on rocker/shaker for 20mins
- 19 Wash x2~3 with PNBS
- 20 Either cover and place in fridge or move directly to immunofluorescence protocol