

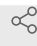



Nov 03, 2022

Passage cells

Adita Ayu Permanasari¹¹Institute of Tropical Disease, Airlangga University

1 Works for me

 Sharedx.doi.org/10.17504/protocols.io.n92ldpbd7l5b/v1 Adita Ayu Permanasari

ABSTRACT

Passage Cells

DOI

dx.doi.org/10.17504/protocols.io.n92ldpbd7l5b/v1

PROTOCOL CITATION

Adita Ayu Permanasari 2022. Passage cells. **protocols.io**
<https://dx.doi.org/10.17504/protocols.io.n92ldpbd7l5b/v1>



LICENSE

————— This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](https://creativecommons.org/licenses/by/4.0/), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

CREATED

Sep 06, 2022

LAST MODIFIED

Nov 03, 2022

PROTOCOL INTEGER ID

69611

- 1 Observe cells under microscope
- 2 Wipe outside surface of cell culture dish with ethanol cotton and bring it inside Biosafety Cabinet (BSC)

- 3 Bring medium, PBS, Trypsin-EDTA 2x to BSC
- 4 Bring autopipetter, disposable pipet (10ml, 2ml), one 15ml tube to BSC
- 5 Discard old medium
- 6 Rinse the bottom of dish with 10ml PBS
- 7 Add 1ml Trypsin-EDTA 2x to dish
- 8 Incubate the dish in CO2 Incubator 37oC 4min
- 9 Check under microscope (cell shape become rounded)
- 10 Neutralize trypsin-EDTA 2x by add 8ml medium
- 11 Suspend detached cells and transfer to 15ml tube
- 12 Centrifuge 1200 rpm, 4 min
- 13 Bring to BSC, discard medium

14 Add 20 ml DMEM and divided each 10 ml to dish (total 2 dish with @50% cells)