



FEB 01, 2024

Purification of OPTN-GST

Elias Adriaenssens¹

¹Sascha Martens lab, University of Vienna, Max Perutz Labs - Vienna

ASAP Collaborative Research Network



Elias Adriaenssens

Sascha Martens lab, University of Vienna, Max Perutz Labs - ...

ABSTRACT

This protocol describes purification of OPTN-GST.

ATTACHMENTS

[pphpb32rp.pdf](#)

OPEN ACCESS



DOI:

dx.doi.org/10.17504/protocols.io.dm6gp3nb8vzp/v1

Protocol Citation: Elias Adriaenssens 2024. Purification of OPTN-GST. **protocols.io** <https://dx.doi.org/10.17504/protocols.io.dm6gp3nb8vzp/v1>

License: This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](#), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

Protocol status: Working

Created: Jan 31, 2024

Last Modified: Feb 01, 2024

Keywords: ASAPCRN

Funders Acknowledgement:
Aligning Science Across
Parkinson's (ASAP)
Grant ID: Grant ID: ASAP-000350
Marie Skłodowska-Curie MSCA
Postdoctoral fellowship
Grant ID: 101062916

Materials:

- pETDuet-1 vector (Vector encoding OPTN-GST is available from Addgene)
- Glutathione Sepharose 4B beads (GE Healthcare)
- Amicon filter (Merck Millipore)
- Superdex 200 Increase 10/300 GL column (Cytiva)

Lysis buffer:

A	B
HEPES pH 7.4	50 mM
NaCl	300 mM
MgCl ₂	2 mM
Glycerol	5%
Imidazole	10 mM
β-mercaptoethanol	2 mM
cOmplete EDTA-free protease inhibitors (Roche)	
CIP protease inhibitor (Sigma)	
DNase (Sigma)	

Wash Buffer:

A	B
HEPES pH 7.4	50 mM
NaCl	300 mM
DTT	1 mM

High-salt wash buffer:

A	B
HEPES pH 7.4	50 mM
NaCl	700 mM
DTT	1 mM

SEC buffer:







A	B
HEPES pH 7.4	25 mM
NaCl	150 mM
DTT	1 mM

Equipment

- SORVAL RC6+ centrifuge
- Beckman Ti45 rotor

Purification of OPTN-GST

20h 45m 30s

- 1 To purify OPTN-GST, Clone human OPTN cDNA in a pETDuet-1 vector with an C-terminal GST-tag.
- 2 After the transformation of the pETDuet-1 vector encoding OPTN-GST in E. coli Rosetta pLysS cells, grow the cells in 2xTY medium at  37 °C until an OD₆₀₀ of 0.4 and then continue at  18 °C .
- 3 Once the cells reached an OD₆₀₀ of 0.8, induce protein expression with  50 micromolar (μM) IPTG for  16h  16:00:00 at  18 °C
- 4 Collect the cells by centrifugation and resuspend in lysis buffer.



Lysis buffer

A	B
HEPES pH 7.4	50 mM
NaCl	300 mM
MgCl ₂	2 mM
Glycerol	5%
Imidazole	10 mM
β-mercaptoethanol	2 mM
cOmplete EDTA-free protease inhibitors (Roche)	
CIP protease inhibitor (Sigma)	
DNase (Sigma)	

5 Sonicate cell lysates. 1m

5.1 Sonicate cell lysates twice for 00:00:30 (1/2). 30s

5.2 Sonicate cell lysates twice for 00:00:30 (2/2). 30s

6 Clear the lysates by centrifugation at 18000 rpm, 4°C, 00:45:00 in a SORVAL RC6+ centrifuge with an F21S-8x50Y rotor (Thermo Scientific). 45m

7 Collect and incubate the supernatant with pre-equilibrated Glutathione Sepharose 4B beads (GE Healthcare) for 02:00:00 at 4 °C with gentle shaking to bind OPTN-GST. 2h

8 Centrifuge the samples to pellet the beads and remove the unbound lysate.


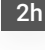



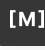
9 Wash beads twice with wash buffer, once with high salt wash buffer, and two more times with wash buffer.

Wash buffer

A	B
HEPES pH 7.4	50 mM
NaCl	300 mM
DTT	1 mM

High-salt wash buffer

A	B
HEPES pH 7.4	50 mM
NaCl	700 mM
DTT	1 mM

10 Elute the proteins  Overnight with  20 millimolar (mM) reduced L-glutathione in  2h  50 millimolar (mM) HEPES  pH 7.4,  300 millimolar (mM) NaCl,  1 millimolar (mM) DTT buffer.

11 Collect the supernatant, filter through a 0.45 µm syringe filter, concentrate using a 50 kDa cut-off Amicon filter (Merck Millipore), and load onto a pre-equilibrated Superdex 200 Increase 10/300 GL column (Cytiva).


12 Elute the proteins with SEC buffer.

SEC buffer

A	B
HEPES pH 7.4	25 mM
NaCl	150 mM
DTT	1 mM

13 Analyze fractions by SDS-PAGE and Coomassie staining.

14 Pool fractions containing purified OPTN-GST.

15 After concentrating the purified protein, aliquot the protein and snap-frozen it in liquid nitrogen. Store proteins at  -80 °C .

