



FEB 02, 2024

Very low-density lipoprotein receptor (VLDLR)-C-tag purification from HEK293E cells

Andreas Bracher¹, Patricia Yuste-Checa¹, F Ulrich Hartl¹

¹Department of Cellular Biochemistry, Max Planck Institute of Biochemistry



Patricia Yuste-Checa

ABSTRACT

This protocol details how to purify recombinant Very low-density lipoprotein receptor (VLDLR)-C-tag protein from HEK293E cells.

ATTACHMENTS

[VLDLR-C-tag purification protocol_protocols.io.docx](#)

MATERIALS

Buffers

▪ Binding buffer:

A	B
Tris-HCl pH 7.2	20 mM
NaCl	100 mM
CaCl ₂	0.5 mM

▪ Elution buffer:

A	B
Tris-HCl pH 7.0	20 mM
MgCl ₂	2 M
CaCl ₂	2 mM

FreeStyle[®]; 293 Expression Medium **Thermo Fisher Catalog #12338018**

CaptureSelect[®]; tPA Affinity Matrix **Thermo Fisher Catalog #2943430005**

Amersham NAP-25 Columns **Cytiva Catalog #17-0852-01**

OPEN ACCESS



DOI:

dx.doi.org/10.17504/protocols.io.j8nlkoqw1v5r/v1

Protocol Citation: Andreas Bracher, Patricia Yuste-Checa, F Ulrich Hartl 2024. Very low-density lipoprotein receptor (VLDLR)-C-tag purification from HEK293E cells. **protocols.io** <https://dx.doi.org/10.17504/protocols.io.j8nlkoqw1v5r/v1>

License: This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](#), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

Protocol status: Working
We use this protocol and it's working

Created: Jan 25, 2024

Last Modified: Feb 02, 2024


PROTOCOL integer ID: 94601

Keywords: ASAPCRN

Funders Acknowledgement:


Aligning Science Across
Parkinson's
Grant ID: ASAP-000282

VLDLR-C-tag expression

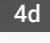1 Express VLDLR-C-tag in HEK293E cells cultured in FreeStyle 293 Expression Medium for  96:00:00 


2 Centrifuge culture and keep conditioned medium.



3 Dialyze  300 mL conditioned medium  Overnight against  10 L Binding buffer.



CaptureSelect C-tag affinity chromatography


4 Load dialyzed conditioned medium onto a CaptureSelect C-tag affinity column previously equilibrated with binding buffer (column volume, CV: 4 mL slurry for 300 mL dialyzed media) by gravity flow at  4 °C .

5 Wash the column with 5 CV of Binding buffer.



6 Elute VLDLR-C-tag protein with 6x 1 mL of Elution buffer. Collect fractions of 1 mL .

7 Analyze eluted fraction by SDS-PAGE and Coomassie blue staining.

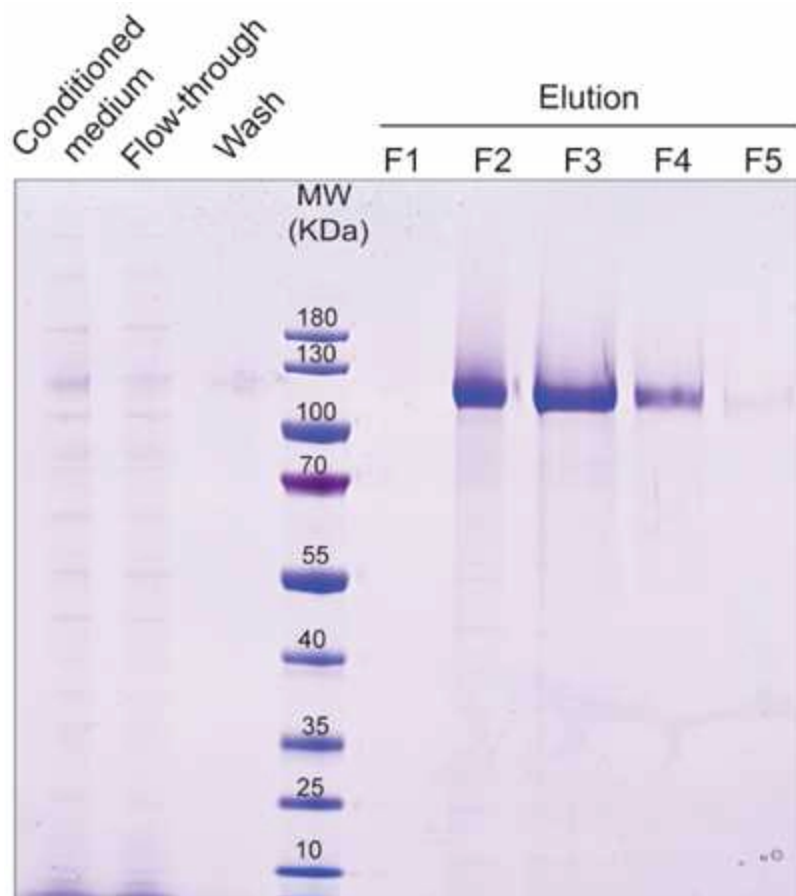

8 Pool fractions containing VLDLR-C-tag protein.

9 Exchange protein buffer to Binding buffer with a NAP-25 desalting column previously equilibrated with Binding buffer.

- 10 Collect the fractions containing protein, concentrate by ultrafiltration to $>1 \text{ mg mL}^{-1}$, aliquot and flash-freeze purified VLDLR-C-tag in liquid nitrogen for storage at -70°C .

Note

Approximate yield: From 300 ml of conditioned media around 0.6 mg of pure VLDLR-C-tag were obtained. Yield can be significantly increased if the VLDLR chaperone, RAP is co-overexpressed during protein production.