



Sep 06, 2022

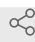
# Adapting hPSCs cultured on MEFs to feeder-free system

 In 1 collection

Hanqin Li<sup>1</sup>, Oriol Busquets<sup>2</sup>, Steven Poser<sup>2</sup>, Dirk Hockemeyer<sup>1</sup>,  
Frank Soldner<sup>2</sup>

<sup>1</sup>University of California, Berkeley; <sup>2</sup>Albert Einstein College of Medicine

1 Works for me

 Share

[dx.doi.org/10.17504/protocols.io.b4k5quy6](https://dx.doi.org/10.17504/protocols.io.b4k5quy6)



Devin E Snyder

## ABSTRACT

This protocol describes the procedure of adapting human pluripotent stem cells (hPSCs) to feeder-free culturing conditions using mTeSR-plus or StemFlex

## General Notes

1. Throughout this protocol, the term hPSC is used to collectively refer to both hiPSCs and hESCs. All described procedures have been tested and work equally well for hiPSCs and hESCs.
2. Until otherwise indicated, feeder-free hPSCs are routinely grown in a humidified cell culture incubator under “low” oxygen conditions. We have successfully maintained hPSCs using either 3% O<sub>2</sub> (3% O<sub>2</sub>, 5% CO<sub>2</sub>) or 5% O<sub>2</sub> (5% O<sub>2</sub>, 5% CO<sub>2</sub>) conditions.
3. We have routinely maintained feeder-free cells in either mTeSR-plus or StemFlex. However, these two mediums are not interchangeable. Pick one and stick to it.
4. We have routinely maintained feeder-free hPSC cultures on VTN, Matrigel and Geltrex-coated cell culture plates without observing obvious differences.

DOI

[dx.doi.org/10.17504/protocols.io.b4k5quy6](https://dx.doi.org/10.17504/protocols.io.b4k5quy6)

## PROTOCOL CITATION

Hanqin Li, Oriol Busquets, Steven Poser, Dirk Hockemeyer, Frank Soldner 2022.  
Adapting hPSCs cultured on MEFs to feeder-free system. **protocols.io**  
<https://dx.doi.org/10.17504/protocols.io.b4k5quy6>



## FUNDERS ACKNOWLEDGEMENT

Aligning Science Across Parkinson's  
Grant ID: ASAP-000486

COLLECTIONS ⓘ



## Feeder-free culturing of hPSCs

KEYWORDS

ASAPCRN

LICENSE

————— This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](#), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

CREATED

Feb 02, 2022

LAST MODIFIED

Sep 06, 2022

PROTOCOL INTEGER ID

57725

PARENT PROTOCOLS

Part of collection

[Feeder-free culturing of hPSCs](#)

MATERIALS TEXT

Item	Vendor	Catalog #
DMEM/F12	Thermo Fisher	11320082
Knockout Serum Replacement (KSR)	Thermo Fisher	10828-028
L-Glutamine	Sigma	G8540
Penicillin & Streptomycin (100x)	Thermo Fisher	15140163
MEM Non-Essential Amino Acids (100X)	Thermo Fisher	11140050
2-Mercaptoethanol	Sigma	M3148
Heat Stable Recombinant Human FGF2	Thermo Fisher	PHG0360
DPBS w/o Calcium and magnesium (DPBS)	Corning	MT21031CV
mTeSR-plus	STEMCELL Technologies	100-0276
StemFlex	Thermo Fisher	A3349401
FBS	Gibco	10437028
Vitronectin (VTN-N) Recombinant Human Protein, Truncated	Thermo Fisher	A14700
DMSO	Fisher Scientific	BP231-100
Y-27632	Chemdea	CD0141
Collagenase type IV	Thermo Fisher	17104019
Matrigel	Corning	CV40234
Geltrex	Fisher Scientific	A1413302

**Note:** This protocol makes reference to protocols in other collections. Please check for any materials found in those protocols, which might not be listed here

- 1 When MEFs-cultured hPSCs reach 50% confluency, change medium to hPSCs medium + Rock inhibitor, preparing for the feeder-free adaptation on the next day.

## 1.1 hPSCs medium

A	B
DMEM/F12	385 ml
Fetal Bovine Serum (FBS)	75 ml
Knockout Serum Replacement	25 ml
L-Glutamine (100X)	5 ml
Penicillin & Streptomycin (100X)	5 ml
MEM Non-Essential Amino Acids (100X)	5 ml
2-Mercaptoethanol (10,000X)	50 µl
Heat Stable Recombinant Human FGF2 (25ug/ml)*	80 µl

\*While we prefer Heat Stable Recombinant Human FGF2, we also have used regular FGF2.  
Final volume: 500ml

### L-Glutamine (100X)

L-Glutamine, powder	14.6 g
MilliQ H2O	500 ml

### 2-Mercaptoethanol (10,000X)

2-Mercaptoethanol	0.78 ml
MilliQ H2O	9.22 ml

### Heat Stable Recombinant Human FGF2 (25µg/ml)

A	B
Heat Stable Recombinant Human FGF2	500 µg
0.1% BSA	20 ml

Final volume: 20ml

### Y-27632 (1,000X)

Y-27632	5 mg
DMSO	1.56 ml

### hPSCs medium + Rock Inhibitor

A	B
hPSCs medium	500 ml
Y-27632 (1,000X)	500 µl

Final volume: 500ml

- 2 Coat three wells of a 6-well plate with either VTN/Matrigel/Geltrex for each cell line.

For a detailed protocol, refer to "Coating plates," which can be found in the protocol collection "Feeder-free culturing of hPSCs." This collection can be accessed using the collection link found in the title section of this protocol, located above

- 3 Wash one well of MEF-cultured hPSCs with DPBS

- 4 Add 1 ml Collagenase solution to this well.

#### 4.1 Collagenase solution



A	B
Collagenase type IV	10 mg
KSR medium	10 ml





Final volume: 10ml

#### KSR medium

A	B
DMEM/F12	385 ml
Knockout Serum Replacement	100 ml
L-Glutamine (200 mM)	5 ml
Penicillin & Streptomycin (100X)	5 ml
MEM Non-Essential Amino Acids (100X)	5 ml

Final volume: 500ml

- 5 Incubate  00:45:00  37 °C . Watch for edge curling of the colonies as this indicates that collagenase incubation is complete. <sup>45m</sup>
- 6 Add 2 ml DMEM/F12

- 7 Pipette repeatedly with 5 ml pipette to lift colonies, careful not to carry over too many MEFs.
- 8 Collect into 15 ml conical tube.
- 9 Add 7 ml DMEM/F12.
- 10 Centrifuge at  **200-300 x g, 00:05:00** 5m
- 11 Aspirate supernatant
- 12 Resuspend cell pellet in 1 ml pre-warmed Accutase
- 13 Incubate  **00:05:00**  **37 °C** 5m
- 14 Add 9 ml DMEM/F12, invert to mix
- 15 Centrifuge at  **200-300 x g, 00:05:00** 5m
- 16 Aspirate supernatant

- 17 Re-suspend cell pellet in 1 ml Feeder-Free Medium + Rock inhibitor, triturate 5-10 times to achieve single cell suspension using a P1000 tip

### 17.1 Feeder-free Medium (version A)

StemFlex basal medium	450 ml
StemFlex supplement	50 ml

Final volume: 500ml

### Feeder-free Medium (version B)

mTeSR-plus basal medium	400 ml
mTeSR-plus supplement	100 ml

Final volume: 500ml

-Feeder-free mediums (version A & B) are not interchangeable. Pick one and stick to it.  
- It is possible to include 5 ml Penicillin & Streptomycin (100X) into the feeder-free medium

### Feeder-free medium + Rock Inhibitor

A	B
Feeder-free medium	50 ml
Y-27632 (1,000X)	50 µl

Final volume: 50ml

- 18 Aspirate VTN/Matrigel/Geltrex solution from the coated plate, add 2 ml Feeder-free medium + Rock inhibitor to each well.
- 19 Dispense 20 µl, 60 µl, 200 µl cell suspension respectively into the three VTN/Matrigel/Geltrex-coated wells.

- 20 Check the cells under the microscope to get an idea of the resulting cell density.
- 21 Spread the cells by moving the plate in left-right, then backward-forward motion.
- 22 Place the plate in a low oxygen incubator
- 23 Change 2 ml pre-warmed Feeder-free medium for each well every other day.
- 24 When large colonies emerge or hPSCs density reaches 50-80%, passage the well showing the best hPSCs morphology using Accutase or ReLeSR. It usually takes 5-7 days.  
  
A detailed protocol on "Passaging of feeder-free hPSCs" can be found in the collection "Feeder-free culturing of hPSCs." This collection can be found using the collection link in the title section of this protocol, located above
- 25 It usually takes 2 passages for hPSCs to fully adapt to feeder-free culture. Differentiation and changes on growth speed are normal during the adaptation.