

Sep 09, 2024

ILT3

This protocol is a draft, published without a DOI.

Mar Roca Cugat¹, UM ILT Team¹

¹Maastricht University



Mar Roca Cugat

Maastricht University, Universitat Autònoma de Barcelona, In...

OPEN  ACCESS



Protocol Citation: Mar Roca Cugat, UM ILT Team 2024. ILT3. [protocols.io https://protocols.io/view/ilt3-dkv44w8w](https://protocols.io/view/ilt3-dkv44w8w)

License: This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](#), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

Protocol status: In development

**We are still developing and
optimizing this protocol**

Created: September 09, 2024

Last Modified: September 09, 2024

Protocol Integer ID: 107164

Abstract

-

Protocol materials

 Dral enzyme (5 U/μl) **New England Biolabs** Step 14

 MilliQ water Step 14

 NEB 10X CutSmart Buffer **New England Biolabs Catalog #B7204S** Step 14




Safety warnings

- ❗ BE SURE to discard all materials (mainly pipette tips) in the biohazard waste bin you are working with GMOs. If you spill contact the supervisors directly: stop working (note your experiment is still safe time is not a huge issue here!).



Bacterial cell destruction




10s

- 1 Resuspend the bacterial pellet in  0.3 mL of Resuspension Buffer.
 - 1.1 Mix vigorously using the vortex until no clumps remain.

The solution needs to be “clear” not “turbid” otherwise invert the tube a few times more and continue.

The lysate should appear viscous. Do not allow the lysis reaction to proceed for more than 5 min.

There may be little treads seen by opening the tube when you start step 3. That is the genomic DNA!


From now on, the bacteria are lysed so no GMO any more (lysed bacteria are no bacteria any more). However, you need to continue working under MLI conditions as the whole laboratory is MLI!
- 2 Add  0.3 mL of Lysis Buffer, mix thoroughly by vigorously inverting the sealed tube for 10 times.
- 3 Add  0.3 mL of chilled Neutralisation Buffer, mix immediately and thoroughly by vigorously shaking  00:00:10 .


If the mixture still appears viscous, more mixing is required to completely neutralize the solution.



10s

Plasmid recovery

12m 10s

- 4 Centrifuge at maximum speed in a centrifuge for  00:02:00 .

Do not disturb the pellet. You can stick your pipette directly through the potentially small white top layer. The white stuff that may stick to the outside of the pipet is no problem as long as you do not transfer this into the new tube, only all the liquid is needed.
- 5 Collect  700 μ L the supernatant in a 1.5 ml microcentrifuge tubes.


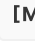




Do not disturb the pellet. You can stick your pipet directly through the potentially small white top layer. The white stuff that may stick to the outside of the pipet is no problem as long as you do not transfer this into the new tube, only all the liquid is needed.
- 6 Precipitate DNA by adding  700 μ L isopropanol to the plasmid DNA containing supernatant. Mix vigorously using the vortex.
- 7 Centrifuge at maximum speed in a centrifuge for  00:07:00 .

2m

7m












Always place the tubes with the hinge of the lid pointed outwards so you know where the pellet should be if you cannot see it in the next steps.

- 8 Pipet away all the fluid. take care NOT to touch or suck up the pellet (your plasmids).
Be care full not to touch your pellet or suck it up. (don't scratch the bottom with your pipet tip) Leave the lowest/last 10 microliters there. no need to remove them now.
- 9 Add  1 mL  70 % (v/v) ethanol to the tube with the plasmid pellet, close the tube and vortex shortly ( 00:00:10). 10s
- 10 Centrifuge at maximal speed for  00:03:00 . 3m
- 11 Remove the 70% ethanol by pipetting.
Use a 200 microliter pipet tip instead of the blue 1000 microliter tips. Pipet the last remaining microliters qualitatively. Do this really carefully to not disturb the plasmid DNA pellet.
- 12 Air dry for  00:05:00 5m
If the pellet just turns from white "opaque" it's already ok! Please set a timer to not over dry the pellet 5 minutes should be ok to evaporate the ethanol the remaining little droplets are water.
- 13 Add  50 μ L sterile water and vortex.
If the vortexing results in scattering all the water though out the tube shortly centrifugate to collect the water with plasmids again. The plasmids are now in solution and will stay in solution regardless the centrifugation.

Restriction enzyme digestion (D2)

30m

- 14 Make the mixtures for the restriction analysis, a total of  20 μ L :
-  DNA sample  5 μ L
 -  NEB 10X CutSmart Buffer **New England Biolabs Catalog #B7204S**  2 μ L
 -  DnaI enzyme (5 U/ μ L) **New England Biolabs**  1 μ L
 -  MilliQ water **Contributed by users**  12 μ L
- 14.1 Mix components by pipetting the reaction mixture up and down, or by "flicking" the reaction tube. Follow with a quick ("touch") spin-down in a microcentrifuge.
Do not vortex the reaction too vigorously -this may kill the restriction enzyme).
Do not discard the remaining plasmid product! You need to analyze this as well!
Do not forget to add the enzyme as last component.



15 Label a reaction tubes and place them in the heating block. Digestion will take place at



37 °C

for



00:30:00 .

30m

16 Add DNA gel loading dye to each sample to stop the reaction.

Agarose gel preparation


17 Prepare the casting tray by putting the sides together.

18 The silicon coated site should point towards the inside.


19 Apply the elastic band around the extruding parts in the middle of the sides to hold everything in place.

20 Flip the assembled tray in the correct position and press it lightly on the table to ensure proper alignment of the both sides.

21 Add the comb and put the tray in the box to catch any potential spills if not correctly assembled.

22 Prepare a 1% agarose gel in 1x TAE buffer. The total volume depends on the size of the casting tray; we will make  100 mL gels. Calculate and weigh the amount of agarose required for your volume. Deposit the agarose in a glass Erlenmeyer flask and add the indicated volume of 1x TAE buffer.

23 Pour the mixture in the casting tray with combs. Remove any air bubbles with a plastic filter tip.

24 Allow the gel to solidify ( 00:10:00).

24.1 Remove the comb.

25 Insert the casting tray into the electrophoresis equipment; the side of the slots in the direction of the cathode (black).

Do not add the metal casting system into the electrophoresis system, just the gel.

26 Prepare and add 1x TAE Buffer to the system until the gel is completely submerged.



27