

603.3 & 604.5_URMC_HTC_Whole Lung and Lobe Processing

Gloria S Pryhuber¹, Heidie Huyck¹

¹University of Rochester Medical Center

1 Works for me dx.doi.org/10.17504/protocols.io.biz7kf9n

Human Cell Atlas Method Development Community LungMap2 Consortium 1 more workspace



Gloria Pryhuber
University of Rochester Medical Center

ABSTRACT

Purpose and Scope of the Procedure

1. Standardize process for processing lung donations into components for storage and distribution
2. Scope: coordination of receipt and gross dissection of tissue

DOI

dx.doi.org/10.17504/protocols.io.biz7kf9n

PROTOCOL CITATION

Gloria S Pryhuber, Heidie Huyck 2020. 603.3 & 604.5_URMC_HTC_Whole Lung and Lobe Processing .
protocols.io
<https://dx.doi.org/10.17504/protocols.io.biz7kf9n>

LICENSE

This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](https://creativecommons.org/licenses/by/4.0/), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

CREATED

Jul 25, 2020

LAST MODIFIED

Aug 16, 2020

PROTOCOL INTEGER ID

39711

GUIDELINES

Scientific Principles

1. Rapid, standardized and safe processing of human tissue for the HTC requires coordination of a team of staff, materials and attention to protocols
2. Rapid processing is required to maintain the tissues in a state as close to normal as possible
3. Ultra-High Resolution CT Scan in an inflated state will provide a high-level comparative assessment of human lung structure across developmental ages

MATERIALS

NAME	CATALOG #	VENDOR
10X PBS; Diluted to 1X with DEPC H2O	Corning 46-013-CM	

NAME	CATALOG #	VENDOR
20% Paraformaldehyde (PFA); Diluted to 4% (m/v) PFA in 1X PBS	EMS 15713-S	
10% Buffered Formalin	VWR 89370-094	
Diethylpyrocarbonate (DEPC)	Sigma D5758-100ML	
Ethanol	Koptec V-1001	
OCT (inflation); 50% (v/v) OCT in 1X PBS	Tissue-Tek 4583	
OCT (mold)	Tissue-Tek 4583	
Sucrose; 30% (m/v) Sucrose in 1X PBS	Sigma S9378-5KG	
Carboxymethylcellulose (CMC); 5% (w/w) CMC in DEPC H2O	Sigma C5678-1KG	

MATERIALS TEXT

Worksheet 603.A.2 HTC_Whole_or_Partial_Lung_Processing Worksheet

- a. Grossing Station/ Biosafety Cabinet or Fume hood
- b. Gloves, face protection and clothing consistent with blood and body fluid precautions
- c. Heavy Duty Scissors and Wire Cutters
- d. Biohazards Disposal Bag
- e. Red Sharps Container
- f. Biosafe surface for manipulation of lung
- g. Small and Medium Gauze Pads – multi-pack
- h. Standard Balance and Weigh Boats: Medium and Large
- i. Airway Cannulas – Selection of sizes
 - i. 16 and 18 gauge angiocatheters (needles removed and discarded)
 - ii. tracheal cannulas 2.0-4.0 OD
 - iii. endotracheal tubes 2.5 – 8 mm OD; tubes >= 4.0 cuffed
- j. IV extension set tubing with clamp
- k. 20 ml syringe barrels for fixative
- l. Needle free suture material
- m. Zip Ties
- n. Dissection Instruments: Same or similar to Sakura Series:
 1. #4791 Scalpel Handle with
 2. #4792 or #4793 Scalpel Blades, #61 (curved tip) or #62 (pointed tip), resp
 3. #4785 Trimming Blades, Short and/or #4789 Trimming Blades, Long
 4. #4786 Trimming Handle, Short and/or #4790 Long Straight
 5. #4794 Blade Scissors with #4796 Blades Sharp/Blunt
 6. #4807 Accu-Edge® Grossing Fork 2.5 mm
 7. #4800 Accu-Edge® Grossing Board *Inches (L x W x D)*: 17 x 11.5 x 1
 8. #4802 Accu-Edge® Grossing Wells *Inches (L x W x D)*: 12.9 x 3.25 x 2.15
 9. #8031 Slide Printer
 - a. #8033 Slide Unload Station
 - b. #8035 Ink Cartridge
 - c. #8040 Slide Magazine
10. Mopec Product #AB079 ProCUT Forcep Kit
- o. Camera – iPad in water proof case works well
- p. 10% neutral buffered formalin – at least one liter with more on hand if needed for larger lungs
- q. 1xPhosphate Buffered Sodium (PBS diluted in DEPC treated water)
- r. 4% paraformaldehyde (prepared from 1:4 dilution of 16% stock non-buffered, no menthol formalin in 1x PBS in DEPC water)
- s. Liquid Nitrogen in Dewar flask
- t. Ring stand(s) and clamps (2 or more) – Fill syringe to 20 ml mark with fixative
- u. Rulers – metal 18 inch (45cm) and 12 inch (30.5cm)

SAFETY WARNINGS

Personnel will adhere to safe work processes outlined in U.S. Public Health Universal Precautions Guidelines for use of human blood and body fluids. PPE will be used including lab coat, closed shoes and gloves. All activity will be behind shield of biosafety cabinet and/or a grossing station with mask and safety glasses. Biosafety level 2 practices will be followed and the work performed in the designated lab space that is covered by annually updated IBC approved protocol. All institutional biosafety measures are followed in any manipulation of these human

tissues.

1 1. Assemble Team and Materials as needed for processing, to accomplish:

- i. CT Scan
- ii. Separation of Lobes
- iii. Tissue Dissociation
- iv. Formalin Inflation
- v. OCT /CMC Inflation
- vi. Paraformaldehyde Inflation

2 Record details of procedures in Worksheets.

Unpacking

- 3 Shipped package is opened on arrival to remove a blood sample, if provided, as it should be kept at room temperature. The remainder of the package is to be kept in cold room until processing begins to maintain approximately **4 °C** temperature of tissue.
- 4 Determine BRINDL PID to assign to tissues and images by clicking the “arrived” button in BRINDL Screening log
- 5 Record Information on Shipping Labels in Database or worksheet for entry ASAP
 - i. UNOS #; Referring Company #
 - ii. Courier and Tracking Numbers
- 6 Remove Shipping Labels from package; place aside to be added to Sample File
- 7 Open Shipping Box and Remove Styrofoam Cover
 - i. Note Condition of Packing in Database (ice sufficient, layered packaging, no leakage, etc.)
- 8 6. Open container containing tissue and move lung tissue to a clean plastic bag on ice.
 - i. Note Condition of Lung Tissue in Database
 - a. Is tissue submerged (usually in plastic bag, inside a hard plastic container)
 - b. Are Trachea and both lungs included and intact
 - c. Is Trachea Occluded, if so how
 - d. Are Thymus, Spleen, Blood Sample, Lymph Nodes included?

CT Scan

- 9 See separate protocol for CT Scan of Air Inflated Lung procedures.

Grossing of Tissues


- 10 Maintain tissues in transplant buffer on ice. Prevent open Trachea from being submerged , +/- ETT removal
- 11 Weigh intact organ as sent; Continue to Record all measurements and procedures in Worksheet or directly in BRINDL database
- 12 Determine displacement volume without inflation
- 13 Determine plan for organ lobes ie which to be inflated or infused and which to be dissociated
- 14 Remove excess tissues around trachea and large airways at hilum
Collect lymph nodes, excess large vessels, esophagus and nerves
(Best to place these back in storage buffer and Tend to Lung First)
 - a. Weigh each tissue collected
 - b. The non-lung samples are sectioned into approximately 0.5-1 cm³ portions
 - c. Divide these tissues between
 - i. 10% formalin, place in labeled tissue cassette in formalin
Store at **4 °C** to fix for 44-48 hrs,
Section at mid-time.
Volume of fixative to tissue should be approximately 10:1.
 - ii. PFA: place in 4%PFA stored at **4 °C** to fix for 20-28 hours before moving to 30%Sucrose
 - d. Record these tissues and how processed in BRINDL database
- 15 Dissect free the bronchial branches, identify pulmonary arteries and veins
- 16 Divide airways and vessels to obtain 5 independent lung lobes
 - 16.1 Weigh each lobe once isolated
 - 16.2 Proceed with processing for each lobe as Outlined in Table 1 using appropriate Protocols
- 17 Prepare accompanying tissue such as thymus, spleen, esophagus, blood, heart, in similar fashion with specific as outlined in the appropriate Protocol / SOP.

Organ	Starting Tissue	Prep	Process	BioSpecimen Type Collected from Organ
Trachea and large bronchi	Intact	Paraffin and OCT Frozen	Bisect longitudinally; Section into <1 cm long half rings	Formalin Fixed, Paraffin embedded PFASOCT embedded and frozen
Alternate Trachea and large bronchi	Intact	Cell Isolates	Enzymatic digestion and cell isolation	Frozen isolated cell aliquots
Lung: Right Upper Lobe +/- Right Middle Lobe	Fresh Tissue Blocks	Mixed Dissociated Cells	Dissociate cells by mechanical and enzymatic digestion procedure	Cryopreserve mixed cell aliquots
Sorted Dissociated Cells	Flow Activated Cell Sorting dissociated cells	Cryopreserve sorted cell aliquots "ex. Endothelial, Epithelial, Fibroblast, Bone Marrow Derived (CD45+) Separation"		
Cryo-preserved QI	Cryo-preserved cells	Recover frozen cell aliquots, test purity by PCR, culture for viability and proliferative/ differentiation potential		
Add'l Right Middle Lobe	Intact - Biopsies	Flash Frozen	Stored at -80	Protein Homogenates, RNA and/or DNA Isolation
Alternate Right Middle Lobe	Intact	Formalin Inflation, Paraffin	Inflation fixed with 10% buffered formalin x 24 hr (Intact Lobe); Sectioned by grid into 0.5-1 cm ³ pieces; Dehydrated to 70% Ethanol	Inflation fixed, paraffin embedded ready for sectioning

Alternate Right Middle Lobe	Fresh Tissue Blocks	Fresh Frozen Embedded in OCT/CMC/none	Vertical 0.1 cm slices, Section each slice by grid into 1 cm2 blocks, alternate slices to embed blocks in OCT or approx. 5% CMC (see protocol) or freeze without embedding material	Non-inflated, non-fixed flash frozen tissue blocks, mapped to originating location
Right Lower Lobe	Intact	Formalin Inflation, Paraffin	Inflation fix with 10% buffered formalin x 24 hr (Intact Lobe) Sectioned by grid into 0.5-1 cm3 pieces Dehydrated to 70% Ethanol	Inflation fixed, paraffin embedded ready for sectioning
Left Upper Lobe	Intact	PFA Fixed, Sucrose Cryopreserved, OCT embedded	Inflation fixed with 4% paraformaldehyde x 24 hr (Intact Lobe); Sectioned by grid into 0.5-1 cm3 pieces; Cryopreserved in 10% sucrose; Block in OCT	Inflation fixed, cryo-preserved and frozen for cryo-sectioning
Alternate Left Upper Lobe	Intact	OCT/PBS	Inflate with 50:50 OCT:10%PBS	Frozen for cryo-sectioning
Alternate Left Upper Lobe	Intact	Air Inflated, Vascular Contrast	To be determined	To be determined
Left Lower Lobe	Intact	PFA Fixed, Sucrose Cryopreserved, OCT embedded	Inflation fixed with 4% paraformaldehyde x 24 hr (Intact Lobe); Sectioned by grid into 0.5-1 cm3 pieces; Cryopreserved in 10% sucrose; Block in OCT	Inflation fixed, cryo-preserved and frozen for cryo-sectioning
Additional Tissue	Intact	Similar to Lung, see specific SOPs	These include thymus, spleen, esophagus, heart, nerve, etc; Block as appropriate for tissue	Formalin Fixed, Paraffin embedded , PFA/OCT embedded -frozen Single cell dissociation
Peripheral Blood	Whole	Cryopreserved, see specific SOPs	PBMC isolation, Plasma and Serum Storage	Cryopreserved

Table 1: Examples of LungMAP Human Tissue Core BioSpecimen Collection, Processing and Handling at Collection Site


19

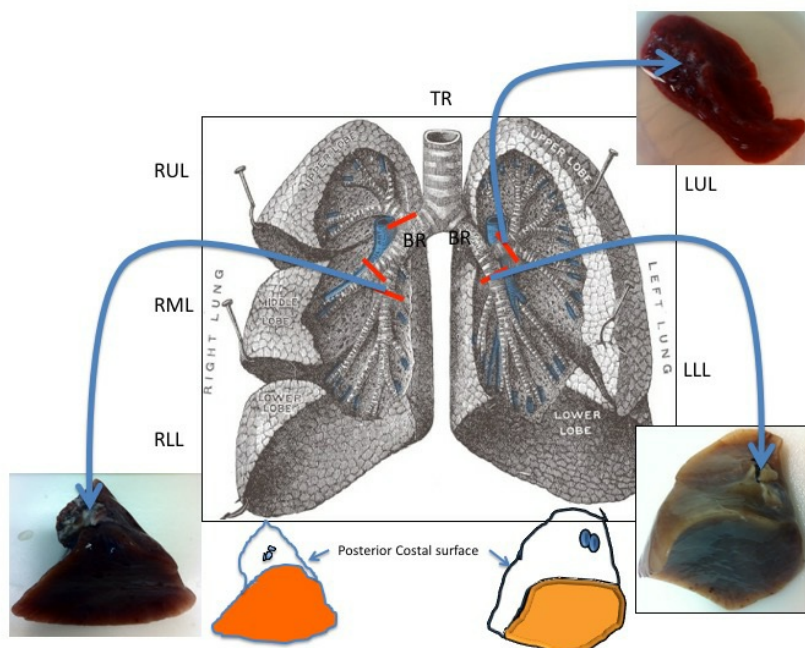

Separation into lobes at bronchi


20

Each Lobe has a characteristic shape, that can be diagramed like this

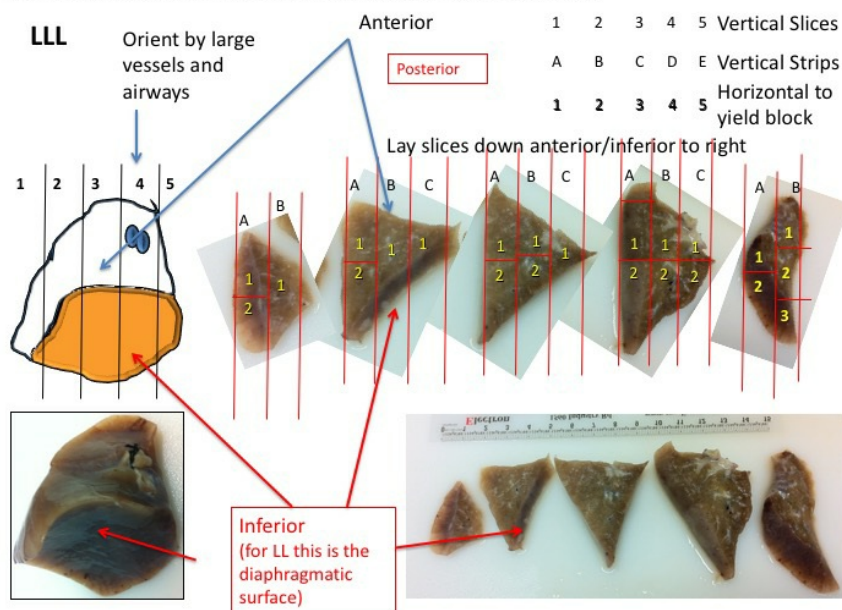


Left Lower Lobe Tissue Blocking Strategy

21

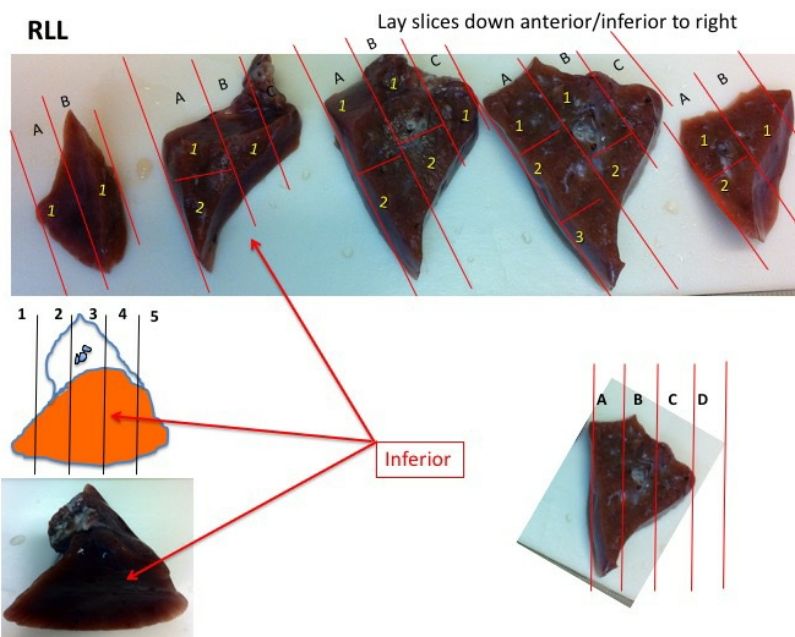


Right Lower Lobe Tissue Blocking Strategy


22

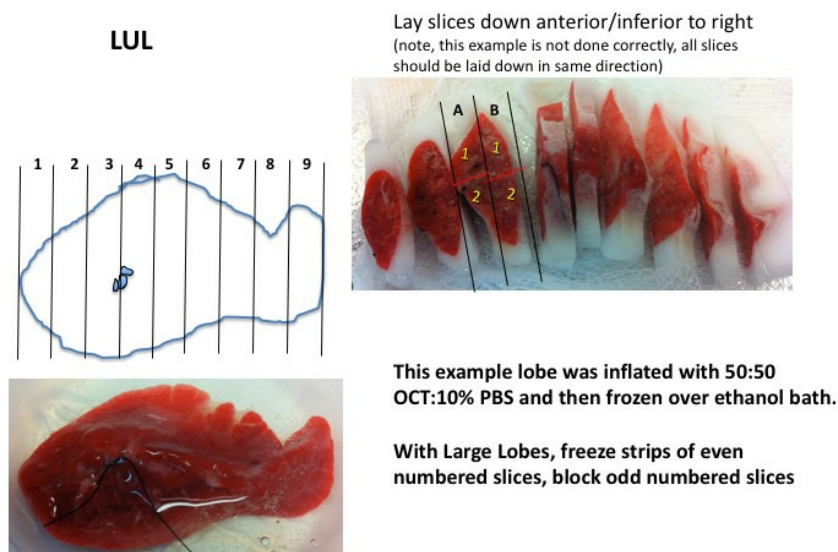


Left Upper Lobe Tissue Blocking Strategy

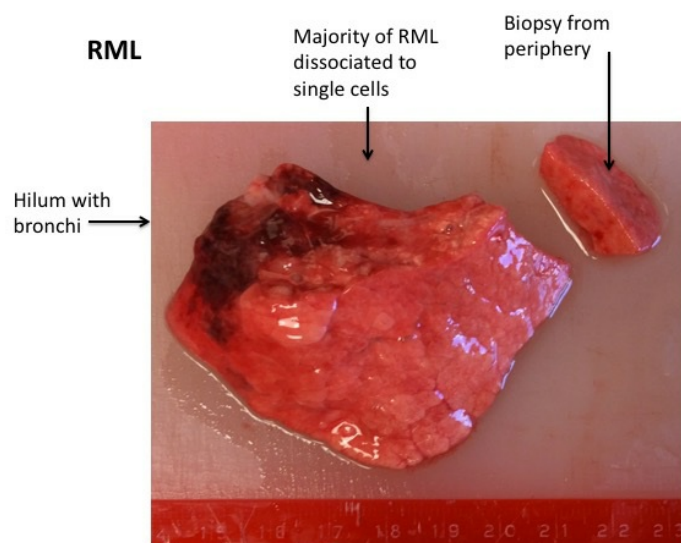


Right Middle Lobe Biopsy and Dissociation Strategy

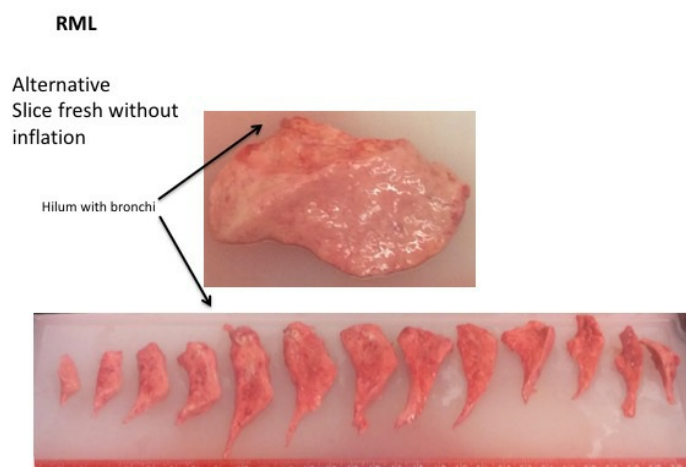

Right Middle Lobe Tissue Blocking Strategy - alternating slices flash frozen / FFPE / PFAFrozenOCT

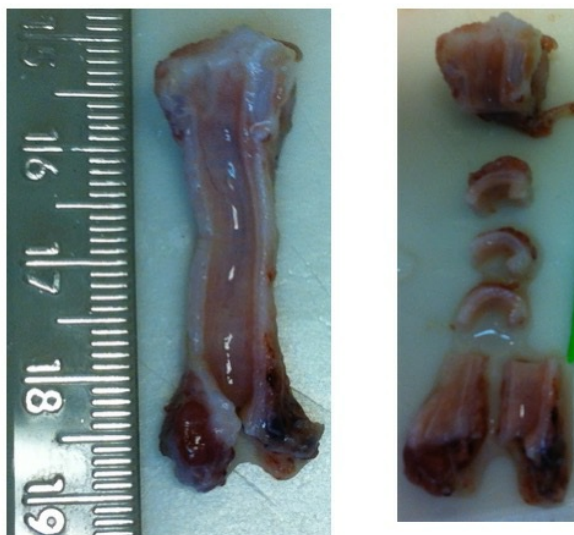

Trachea and Bronchi Tissue Blocking Strategy - image shows one side of bisected trachea only

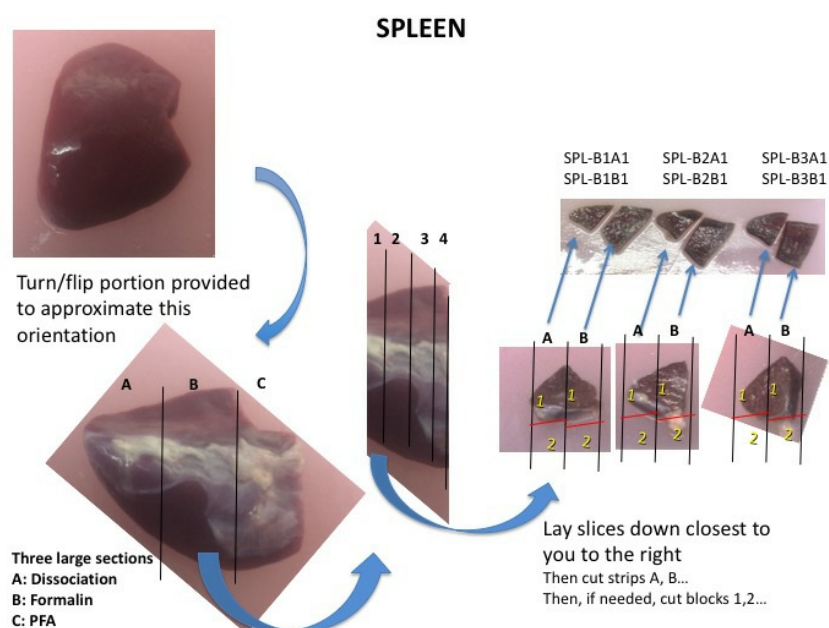

Spleen Tissue Blocking Strategy