

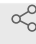



Sep 21, 2022

Flow cytometry

[gurvir.virdi](#)¹¹UCL Institute of Neurology

1 Works for me

 Sharedx.doi.org/10.17504/protocols.io.kqdg3956eg25/v1 gurvirdi

ABSTRACT

This protocol is for immunolabelling fixed midbrain dopaminergic neurons for flow cytometry analysis and downstream flow cytometry acquisition.

DOI

dx.doi.org/10.17504/protocols.io.kqdg3956eg25/v1

PROTOCOL CITATION

gurvir.virdi 2022. Flow cytometry. **protocols.io**
<https://protocols.io/view/flow-cytometry-cgxgtxjw>



LICENSE

————— This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](#), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

CREATED

Sep 21, 2022

LAST MODIFIED

Sep 21, 2022

PROTOCOL INTEGER ID

70344

Sample preparation

- 1 Cells were washed once with PBS
- 2 The cells were incubated with Accutase (Gibco) to generate a single-cell suspension for 5m

🕒 00:05:00 at 🌡 37 °C

- 3 A cell suspension of 500k/ml was prepared in media
- 4 Cells were then spun down at 🌀 200 x g, 00:05:00 , and the supernatant was removed. ^{5m}
- 5 Cell pellet was resuspended gently in 4ml of 4% paraformaldehyde and briefly vortexed at a ^{10m} low speed before being spun on a rotation spinner for 🕒 00:10:00 at room temperature.
- 6 After fixation, samples were spun down and supernatant removed
- 7 Cells were resuspended in 2ml of 0.1% BSA (Sigma) in PBS.
- 8 After resuspension, cells were filtered through a 70µm strainer (Miltenyi Biotec) to filter out any cell clumps
- 9 Cells were then centrifuged 🌀 200 x g, 00:05:00 , and the supernatant was removed. ^{5m}
- 10 Cell pellets were then resuspended in 1ml of permeabilization/blocking buffer (0.1% Triton X-100, 1% BSA, 10% normal goat serum (Sigma) in PBS), and incubated on a rotation spinner for ^{30m} 🕒 00:30:00 at room temperature.
- 11 After permeabilization/blocking, cells were centrifuged 🌀 200 x g, 00:05:00 and the ^{5m} supernatant was removed.
- 12 Cells were then resuspended in the primary antibodies (1:200) made up in 0.1% BSA in PBS, ^{1h}

and incubated on the rotation spinner for 🕒01:00:00 at room temperature.

- 13 After primary antibody incubation, cells were centrifuged 🌀200 x g, 00:05:00 , supernatant^{5m} removed and washed once in 0.1% BSA in PBS.
- 14 Cells were then resuspended in the species-specific secondary antibodies (AlexaFluor 488,^{30m} 647) at a dilution of 1:500 made up in 0.1% BSA in PBS and incubated in the dark on a rotation spinner for 🕒00:30:00 .
- 15 After incubation, cells were centrifuged 🌀200 x g, 00:05:00 , supernatant removed and^{10m} washed once in PBS, followed by incubation with DAPI made up in PBS for 🕒00:05:00 .
- 16 The DAPI + PBS was then removed, followed by one wash in PBS, before being analysed on the flow cytometer.

Cell sorting and analysis

- 17 The samples were run on the LSRii (BD) cell sorter. Scattering was initially used to discard debris as well as cell doublets and larger clumps.
- 18 The single-cell population was then gated to include DAPI positive only cells (negative control).
- 19 The gating threshold for measured channels was determined using the control lacking the antibody of interest (Fluorescence minus one (FMO) control), for both channels being recorded.
- 20 Once the parameters had been set, 10,000 cell events were recorded, and data were processed and analysed on the FlowJo software.