



JAN 10, 2024

OPEN ACCESS



DOI:
dx.doi.org/10.17504/protocols.io.4r3l221k4l1y/v1

Protocol Citation: Jonathan Tang 2024. Single action sequence reinforcement . **protocols.io** <https://dx.doi.org/10.17504/protocols.io.4r3l221k4l1y/v1>

MANUSCRIPT CITATION:
 Tang, J.C.Y., Paixao, V., Carvalho, F. *et al.* Dynamic behaviour restructuring mediates dopamine-dependent credit assignment. *Nature* (2023). <https://doi.org/10.1038/s41586-023-06941-5>

License: This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](#), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited


Single action sequence reinforcement

Jonathan Tang^{1,2,3}

¹Department of Neuroscience, Zuckerman Mind Brain Behavior Institute, Columbia University, New York, NY 10027, USA;

²Seattle Children's Research Institute, Center for Integrative Brain Research, Seattle, WA, USA;

³University of Washington School of Medicine, Department of Pediatrics, Seattle, WA, USA

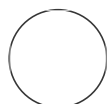

c.weber-schmidt
 Allen Institute, ASAPCRN

ABSTRACT

Single action sequence reinforcement protocol for mouse studies from Tang et al 2023.

GUIDELINES

Individual mice were subjected to a single session of protocol each day, with sessions following each other on consecutive days.

MATERIALS

The white open field box is made of 410 x 400 mm white matte acrylic walls and a 410 x 400 mm white matte acrylic base. **Note:** A white open field was used instead of a grey open field to minimize habituation effects, which would lead to reduced initial spontaneous behavior during closed loop protocol.

Protocol status: Working
We use this protocol and it's working

Created: Dec 21, 2023

Last Modified: Jan 10, 2024

PROTOCOL integer ID:
92643

Keywords: ASAPCRN,
mouse, behavior, brain

Funders

Acknowledgement:

Life Sciences Research
Fellowship

NINDS K99/R00 Award
Grant ID: 1K99NS112575

NIH
Grant ID: 5U19NS104649

Aligning Science Across
Parkinson's (ASAP)
Grant ID: ASAP-020551

Session 1 + Baseline

- 1** Mice were placed in a white open field box for closed loop reinforcement protocol. To acquire baseline behavior, individual mice were allowed to behave freely inside the box for 30 minutes on the first action A reinforcement session.
- 2** After initial behavior acquisition, closed loop reinforcement was initiated by blue laser stimulation of VTA DA neurons were made available for 60 minutes.

Sessions 2 & 3

- 3** 90 minutes of closed loop reinforcement were made available for individual mice during sessions.

Session 4

- 4** An extinction protocol was carried out comprising of 20-minute maintenance of reinforced behavior with

laser availability, followed by 60 minutes of extinction of reinforced behavior without laser availability, followed by 20-minute re-acquisition of reinforced behavior with laser availability.

Contingency Degradation Protocol

STEP CASE

Selection for Action B

1 step

To select for action B, a repeat of the protocol described above for action A was performed starting on the day following extinction protocol of action A.

- 5 Upon completion of the reinforcement and extinction protocols for action B, a contingency degradation protocol was performed comprising of 20-minute maintenance of action B with laser availability, followed by 60 minutes of contingency degradation of reinforced behavior by triggering laser randomly, followed by 40-minute re-acquisition of reinforced behavior with laser availability for action B performance.