




May 28, 2021

Isolation of Chlamydia from crocodile tissue samples

Somjit Chaiwattananarungpaissan¹, Metawee Thongdee¹, ladawan.sar¹¹The Monitoring and Surveillance Center for Zoonotic Diseases in Wildlife and Exotic Animals, Faculty of Veterinary Science, Mahidol University, Nakhon Pathom, Thailand

1 Works for me

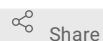

Share

dx.doi.org/10.17504/protocols.io.btcqnivw

ladawan.sar

ABSTRACT

This protocol is used for isolation of Chlamydia from crocodiles

EXTERNAL LINK

<https://doi.org/10.1371/journal.pone.0252081>

DOI

dx.doi.org/10.17504/protocols.io.btcqnivw

EXTERNAL LINK

<https://doi.org/10.1371/journal.pone.0252081>

PROTOCOL CITATION

Somjit Chaiwattananarungpaissan, Metawee Thongdee, ladawan.sar 2021. Isolation of Chlamydia from crocodile tissue samples. **protocols.io**
<https://dx.doi.org/10.17504/protocols.io.btcqnivw>



LICENSE

This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](https://creativecommons.org/licenses/by/4.0/), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

CREATED

Mar 16, 2021

LAST MODIFIED

May 28, 2021

PROTOCOL INTEGER ID

48240

- 1 Dissect tissue sample (approximate 50-100 mg) into small pieces.
- 2 Homogenize the small pieces of tissue sample with a pestle in 5 ml sucrose/phosphate/glutamate (SPG) buffer containing 500 ug/ml streptomycin, 500 ug/ml vancomycin, 50 ug/ml gentamycin, and 2.5 ug/ml fungizone and left at 4°C for 72 h.
- 3 Centrifuge the homogenized tissue sample at 250×g for 10 min and collect supernatant for cell culture inoculation.

- 4 Prepare monolayers of McCoy cells in 12-well plates by seeding 3.5×10^5 cell per well in cell culture medium (M199 medium (Gibco BRL Life Technologies Inc., NY, USA) supplemented with 10% heat-inactivated fetal calf serum (Gibco BRL Life Technologies Inc., NY, USA), 0.2% NaHCO_3 , $1 \times \text{L-glutamine}$, 100 $\mu\text{g/ml}$ streptomycin).
- 5 Inoculate 400 μl of supernatant into each well of monolayers cells.
- 6 Centrifuge 12-well plate containing infected cells at $1,000 \times g$ and 25°C for 1 h.
- 7 After centrifugation, incubate the infected cultures at 30°C in an atmosphere of 5% CO_2 for an additional 1 h.
- 8 Replace the inocula with fresh incubation medium (M199 medium (Gibco BRL Life Technologies Inc., NY, USA) containing 0.5 $\mu\text{g/ml}$ cycloheximide and supplemented with 10% heat-inactivated fetal calf serum (Gibco BRL Life Technologies Inc., NY, USA), 0.2% NaHCO_3 , 10% glucose, $1 \times \text{L-glutamine}$, 100 $\mu\text{g/ml}$ streptomycin, 100 $\mu\text{g/ml}$ vancomycin, 10 $\mu\text{g/ml}$ gentamycin, and 1 $\mu\text{g/ml}$ fungizone).
- 9 Further incubate the infected cultures at 30°C in an atmosphere of 5% CO_2 for 4–5 days.
- 10 Observe daily for the cytopathic effect (CPE) under an inverted microscope.

Chlamydia propagation

- 11
- 12 Prepare monolayers of McCoy cells in 12-well plates by seeding 3.5×10^5 cell per well in cell culture medium (M199 medium (Gibco BRL Life Technologies Inc., NY, USA) supplemented with 10% heat-inactivated fetal calf serum (Gibco BRL Life Technologies Inc., NY, USA), 0.2% NaHCO_3 , $1 \times \text{L-glutamine}$, 100 $\mu\text{g/ml}$ streptomycin).
- 13 Thaw frozen stock of Chlamydia culture and place on ice.
- 14 Dilute culture stock to 1:10 in cell culture medium.
- 15 Inoculate 200 μl of diluted culture stock into each well of monolayers cells and add 200 μl of cell culture medium.
- 16 Centrifuge 12-well plate containing infected cells at $1,000 \times g$ and 25°C for 1 h.

- 17 After centrifugation, incubate the infected cultures at 30°C in an atmosphere of 5% CO₂ for an additional 1 h.
- 18 Replace the inocula with fresh incubation medium (M199 medium (Gibco BRL Life Technologies Inc., NY, USA) containing 0.5 µg/ml cycloheximide and supplemented with 10% heat-inactivated fetal calf serum (Gibco BRL Life Technologies Inc., NY, USA), 0.2% NaHCO₃, 10% glucose, 1× L-glutamine, 100 µg/ml streptomycin, 100 µg/ml vancomycin, 10 µg/ml gentamycin, and 1 µg/ml fungizone).
- 19 Further incubate the infected cultures at 30°C in an atmosphere of 5% CO₂ for 4–5 days.
- 20 Observe daily for the CPE under an inverted microscope.

IFA

- 21
- 22 Prepare monolayer McCoy cells on glass coverslip (Ø 15 mm) in 12-well plates by seeding 3.5×10⁵ cell per well in cell culture medium (M199 medium (Gibco BRL Life Technologies Inc., NY, USA) supplemented with 10% heat-inactivated fetal calf serum (Gibco BRL Life Technologies Inc., NY, USA), 0.2% NaHCO₃, 1× L-glutamine, 100 µg/ml streptomycin).
- 23 Thaw frozen stock of Chlamydia culture and place on ice.
- 24 Dilute culture stock to 1:20 in cell culture medium.
- 25 Inoculate 400 µl of diluted culture stock into each well of monolayers cells.
- 26 Centrifuge 12-well plate containing infected cells at 1,000×g and 25°C for 1 h.
- 27 After centrifugation, incubate the infected cultures at 30°C in an atmosphere of 5% CO₂ for an additional 1 h.
- 28 Replace the inocula with fresh incubation medium (M199 medium (Gibco BRL Life Technologies Inc., NY, USA) containing 0.5 µg/ml cycloheximide and supplemented with 10% heat-inactivated fetal calf serum (Gibco BRL Life Technologies Inc., NY, USA), 0.2% NaHCO₃, 10% glucose, 1× L-glutamine, 100 µg/ml streptomycin, 100 µg/ml vancomycin, 10 µg/ml gentamycin, and 1 µg/ml fungizone).

- 29 Further incubate the infected cultures at 30°C in an atmosphere of 5% CO₂ for 3-4 days.
- 30 Observe daily for the CPE under an inverted microscope.
- 31 When CPE in cell culture occur approximately 10%, fix the infected cells with 1 ml of 4% paraformaldehyde for 24 h.
- 32 Remove 4% paraformaldehyde and wash fixed samples 3 times with 1 ml of 1× phosphate-buffered saline (PBS).
- 33 Add 1 ml of 0.2% Triton X-100 in PBS and incubate at 37°C for 30 min for permeabilization.
- 34 Wash fixed samples 3 times with 1 ml of 1×PBS.
- 35 Add 1 ml of 1% BSA and incubate at 37°C for 30 min for blocking of non-specific binding.
- 36 Transfer the coverslip to a humidity chamber and add 30 µl of 1:100 diluted primary antibody (Chlamydiaceae family-specific mouse monoclonal antibody directed against the chlamydial lipopolysaccharide (cat. no. sc-58106, Santa Cruz Biotechnology, Inc., TX, USA). Incubate the humidity chamber at 37°C for 1 h in the dark.
- 37 Transfer the coverslip to a 12-well plate and wash coverslip 3 times with 1 ml of 1×PBS.
- 38 Transfer the coverslip to the humidity chamber and add 30 µl of 1:200 diluted secondary antibody (Alexa Fluor 594-conjugated secondary goat anti-mouse antibody (cat. no. A11005, Life Technologies, CA, USA) and 2 µg/ml 4',6-diamidino-2'-phenylindole dihydrochloride (DAPI) (Molecular Probes, OR, USA). Incubate the humidity chamber at 37°C for 1 h in the dark.
- 39 Transfer the coverslip to a 12-well plate and wash coverslip 3 times with 1 ml of 1×PBS.
- 40 Mount the coverslip onto microscope glass slides using 50% glycerol.
- 41 Seal the mounted coverslip with nail polish.

42 Visualize the image under an inverted fluorescence microscope.