

Aug 26, 2024

NEBNext® Poly(A) mRNA Magnetic Isolation Module NEB #E7490S/L (Express Protocol)

 In 1 collection

DOI

dx.doi.org/10.17504/protocols.io.n92ld8627v5b/v1

Juliet Bonnevie¹

¹New England Biolabs

New England Biolabs (NEB)

NEBNext



Isabel Gautreau

New England Biolabs

OPEN  ACCESS



DOI: dx.doi.org/10.17504/protocols.io.n92ld8627v5b/v1

External link: <https://www.neb.com/en-us/-/media/nebus/files/manuals/manuale7490.pdf?rev=ed00dd803dc94459a8b437f75b82a42a&hash=CF5F02FF4E55FB9FDCA99A4D1EE58DA0>

Protocol Citation: Juliet Bonnevie 2024. NEBNext® Poly(A) mRNA Magnetic Isolation Module NEB #E7490S/L (Express Protocol). protocols.io <https://dx.doi.org/10.17504/protocols.io.n92ld8627v5b/v1>

License: This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](#), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

Protocol status: Working

We use this protocol and it's working

Created: March 08, 2024

Last Modified: August 26, 2024

Protocol Integer ID: 106283

Keywords: Oligo d(T)25, Automated high-throughput application



Abstract

The NEBNext Poly(A) mRNA Magnetic Isolation Module is designed to isolate intact poly(A)+ RNA from previously isolated total RNA. The technology is based on the coupling of Oligo d(T)25 to 1 µm paramagnetic beads which is then used as the solid support for the direct binding of poly(A)+ RNA. Thus, the procedure permits the manual processing of multiple samples and can be adapted for automated high-throughput applications. Additionally, magnetic separation technology permits elution of intact mRNA in small volumes eliminating the need for precipitating the poly(A)+ transcripts in the eluent. Intact poly(A)+ RNA which is fully representative of the mRNA population of the original sample can be obtained in less than one hour. The new protocol for poly(A) mRNA enrichment (Section 1, Express Protocol) enables comparable performance to our current protocol (Section 2, Standard Protocol) using a substantially faster and more streamlined workflow.

Guidelines

Please refer to the product page on NEB.com for FAQs about this product.

Where larger volumes, customized or bulk packaging are required, we encourage consultation with the Customized Solutions team at NEB.

Please complete the NEB Custom Contact Form at www.neb.com/CustomContactForm to learn more.

Application


Isolation of poly(A)+ RNA transcripts from Total RNA for RNA library preparation and sequencing.




Materials

The NEBNext Poly(A) mRNA Magnetic Isolation Module Includes

The volumes provided are sufficient for preparation of up to 24 reactions (NEB #E7490S) and 96 reactions (NEB #E7490L).

All reagents should be stored at  4 °C .

- NEBNext Oligo d(T)25 Beads
- NEBNext RNA Binding Buffer (2X)
- NEBNext Wash Buffer
- NEBNext Tris Buffer (10 mM Tris HCl,  7.5)
- Nuclease-free Water

Required Materials Not Included

- 96-well 0.2 ml PCR Plates and Microseal® 'B' Adhesive Sealer (Bio-Rad® MSB-1001) or 0.2 ml RNase-free tube, for example Tempassure
- PCR flex-free 8-tube strips (USA Scientific® #1402-4708)
- Magnetic Rack (NEB#S1515S, Alpaqua®, cat. #A001322 or equivalent)
- 1.5 ml Microcentrifuge tube and NEB #S1506 Magnet stand or equivalent (for washing beads only)
- Thermal cycler or heat block
- Bioanalyzer® (Agilent® Technologies, Inc.) or similar instrument and consumables



Before start

Note: When using NEBNext Poly(A) mRNA Magnetic Isolation Module E7490 with NEBNext library preparation kits, follow the protocol in the respective library prep kit manual or see Appendix A in the E7490 manual on [neb.com](https://www.neb.com). Modifications to elution buffer and volumes will need to be made at Step 26 to make the enriched mRNA directly compatible with the different library preparation workflows.

Note: For best results keep all the reagents used during the Poly(A) isolation except the NEBNext Oligo d(T)25 beads, on ice when not in use.

The protocol has been optimized using high quality Universal Human Reference Total RNA.

RNA Integrity

Assess the quality of the Input RNA by running the RNA sample on an Agilent Bioanalyzer RNA 6000 Nano/Pico Chip. For Poly(A) mRNA enrichment, high quality RNA with a RIN score > 7 is required.

RNA Sample Requirements

The RNA sample should be free of salts (e.g., Mg^{2+} , or guanidinium salts, divalent cation chelating agents (e.g., EDTA or EGTA) or organics (e.g., phenol or ethanol). RNA must be free of DNA. gDNA is a common contaminant from RNA preps. It may be carried over from the interphase of organic extractions or when the silica matrix of solid phase RNA purification methods is overloaded. If the total RNA sample may contain gDNA contamination, treat the sample with DNase I to remove all traces of DNA (DNase is not provided in this kit). After treatment with DNase I the enzyme should be removed from the sample. Any residual activity of the DNase I may degrade the oligos necessary for the enrichment. DNase I can be removed from the extraction using phenol/chloroform extraction and ethanol precipitation.






Starting Material: 1–5 µg* of DNA-free total RNA (Express Protocol)



26m



1

Note

* Inputs as low as  10 ng can be used in combination with NEBNext library preparation workflows. Please refer to the respective library prep manual for the appropriate instructions to follow. Modifications to elution buffer and volumes will need to be made here to make the enriched mRNA directly compatible with the different NEBNext library preparation workflows. Please refer to Appendix A in the manual on NEB.com and consult the appropriate manual for more information.




Dilute the total RNA with nuclease-free water to a final volume of  50 µL in a nuclease-free 0.2 ml PCR tube and keep  On ice .


2 Add  20 µL NEBNext Oligo d(T)25 beads per reaction to a 1.5 ml tube. If preparing multiple libraries, beads for up to 24 samples can be added to a single 1.5 ml tube for subsequent washes (use magnet NEB #S1506 for 1.5 ml tubes). The purpose of this step is to bring the beads from the storage buffer into the binding buffer. The NEBNext RNA Binding Buffer (2X) does not have to be diluted for this step. 

3 Place the tube on the magnetic rack at  Room temperature for  00:02:00 or until the solution is clear. 2m











4 Remove and discard all of the supernatant from the tube. Take care not to disturb the beads.



5 Remove the tube from the magnetic rack.








6 Add  50 µL NEBNext RNA Binding Buffer (2X) to the beads and wash by pipetting up and down 6 times. If preparing multiple libraries, add  50 µL RNA Binding Buffer per sample. The Binding Buffer does not have to be diluted. 

7 Place the tube on the magnet and incubate at  Room temperature until the solution is clear (~ 2 minutes).



- 8 Remove and discard the supernatant from the tube. Take care not to disturb the beads.
- 9 Remove the tube from the magnet and add  50 μ L NEBNext RNA Binding Buffer (2X) to the beads and mix by pipetting up and down until beads are homogenous. If preparing multiple libraries, add  50 μ L RNA Binding Buffer per sample. The Binding Buffer does not have to be diluted. 
- 10 Add  50 μ L beads to each RNA sample from Step 1. Mix thoroughly by pipetting up and down 6 times. This first binding step removes most of the non-target RNA.  
- 11 Heat the sample to denature the RNA and to facilitate binding of the poly(A) RNA to the beads. Place in a thermal cycler, with the heated lid set to \geq  90 °C , and run the following program: 
-  00:02:00 at  80 °C

 00:05:00 at  25 °C

Hold at  25 °C
- 12 Remove tubes from the thermal cycler when the temperature reaches hold at  25 °C .
- 13 Place the tubes on the magnetic rack at  Room temperature for  00:02:00 or until the solution is clear to separate the poly(A) RNA bound to the beads from the solution. 
- 14 Remove and discard all of the supernatant. Take care not to disturb the beads. Do not remove the tubes from the magnetic rack
- 15 While still on the magnet rinse the beads by gently adding  200 μ L of NEBNext Wash Buffer to the tubes to remove unbound RNA. 
- 16 Remove and discard all of the supernatant from each tube. Take care not to disturb the beads.

17 Remove the tubes from the magnetic rack.

18 Add the following to each tube containing mRNA bound beads to allow the poly(A) RNA to re-bind to allow the RNA to re-bind to the same beads. Mix thoroughly by gently pipetting up and down 6 times.



A	B
COMPONENT	VOLUME PER ONE LIBRARY
NEBNext Tris Buffer	50 µl
NEBNext RNA Binding Buffer (2X)	50 µl
Total Volume	100 µl

For multiple reactions, a master mix of the reaction components can be prepared On ice before addition to the sample.

19 Place the tubes in a thermal cycler, with the heated lid set to \geq 90 °C, and run the following program:

7m



00:02:00 at 80 °C
 00:05:00 at 25 °C
 Hold at 25 °C

20 Remove the tubes from the thermal cycler when the temperature reaches hold at 25 °C.

21 Place the tubes on the magnetic stand at Room temperature for 00:02:00 or until the solution is clear.

2m



22 Remove and discard all of the supernatant from each tube. Take care not to disturb the beads. Do not remove the tubes from the magnetic rack.

23 While still on the magnet, rinse the beads by gently adding 200 µL of NEBNext Wash Buffer to the tubes to remove unbound RNA.



24 Remove and discard all of the supernatant from the tubes. Take care to remove all of the NEBNext Wash Buffer and do not disturb the beads that contain the mRNA.












**Note**

Note: It is important to remove all of the supernatant to successfully use the RNA in downstream steps. Spin down the tube. Place the tube on the magnetic rack and with a 10 µl tip remove all of the NEBNext Wash Buffer. Caution: Do not disturb beads that contain the mRNA. Avoid letting the beads dry out before adding elution buffer.

- 25 Remove the tubes from the magnetic rack.


Note

Note: Modifications to elution buffer and volumes will need to be made here to make the enriched mRNA directly compatible with the different NEBNext library preparation workflows. Please refer to Appendix A in the E7490 manual on [neb.com](https://www.neb.com) and consult the appropriate manual for more information.

- 26 Elute the mRNA from the beads by adding  17 µL of the NEBNext Tris Buffer, mix by pipetting 6 times and incubate the samples at  80 °C for  00:02:00, then cool to  25 °C. Immediately, place the tubes on the magnetic rack for  00:02:00 or until the solution is clear.
- 27 Collect the purified mRNA by transferring  15 µL of the supernatant to a clean nuclease-free PCR tube. Place  On ice if using immediately. Alternatively, the sample can be placed at  -80 °C for long term storage.
- 28 Assess the yield and the size distribution of the purified mRNA. Run  1 µL on the Bioanalyzer using a RNA Pico Chip or the TapeStation using a High Sensitivity RNA ScreenTape. You may have to dilute your sample before loading.

4m



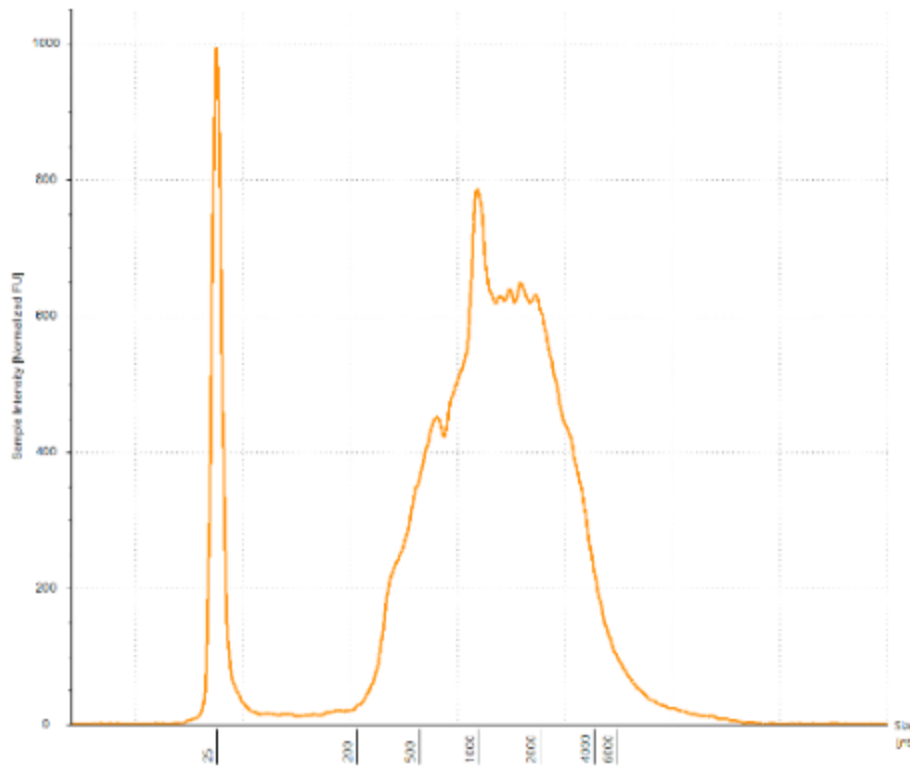

Example of mRNA distribution on a TapeStation