



VERSION 2

FEB 21, 2023

OPEN ACCESS

**Protocol Citation:** Andreas Sagen 2023. Potassium phosphate buffer (0.9 M, pH 7.7). [protocols.io](https://protocols.io/view/potassium-phosphate-buffer-0-9-m-ph-7-7-cprwvm7e) <https://protocols.io/view/potassium-phosphate-buffer-0-9-m-ph-7-7-cprwvm7e> Version created by [Andreas Sagen](#)

**License:** This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](#), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

**Protocol status:** In development  
We are still developing and optimizing this protocol

**Created:** Feb 20, 2023

**Last Modified:** Feb 21, 2023


**PROTOCOL integer ID:**  
77334

**Keywords:** potassium phosphate buffer, buffer, protein

## 🌐 Potassium phosphate buffer (0.9 M, pH 7.7) V.2

Andreas Sagen<sup>1</sup>

<sup>1</sup>University of Oslo

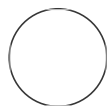

Andreas Sagen

University of Oslo, The National Institute of Occupational H...

### ABSTRACT










Potassium phosphate buffers are buffers with a buffering capacity between 5.8 and 8.0, where it is possible to adjust the pH and buffer strength with different ratios of Sodium phosphate monobasic and Sodium phosphate dibasic.

### MATERIALS

Scale  
Magnetic stirrer  
Filtering device

500 mL Phosphate buffer (0.9 M, pH 7.7)

5m

- 1 Fill a sterile  500 mL bottle with  400 mL distilled water
- 2 Measure  62.705 g Potassium phosphate dibasic and  11.568 g Potassium phosphate monobasic  
  
Materials:  
 Potassium phosphate dibasic **Sigma-aldrich Catalog #P3786**  
 Potassium phosphate monobasic **Sigma-aldrich Catalog #P5379**
- 3 Add dry ingredients and mix for  00:05:00 5m
- 4 Fill to  500 mL with distilled water
- 5 Filter sterilize through a 0.2 µm filter and store refrigerated (  4 °C )