



JUL 31, 2023

OPEN ACCESS



DOI:
dx.doi.org/10.17504/protocols.io.e6nvwjx8dlmk/v1

Protocol Citation: OLIVIA HARDING 2023. Fixing hippo neurons to assess endogenous NEMO during oxidative stress .
protocols.io
<https://dx.doi.org/10.17504/protocols.io.e6nvwjx8dlmk/v1>

License: This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](#), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

Protocol status: Working
 We use this protocol and it's working

Created: May 30, 2023

🌐 Fixing hippo neurons to assess endogenous NEMO during oxidative stress


🔗 Forked from [Investigation_of_mitophagy_in_Hippo_neurons](#)

OLIVIA
 HARDING^{1,2}

¹Department of Physiology, Perelman School of Medicine, University of Pennsylvania, Philadelphia, PA 19104;

²Aligning Science Across Parkinson's (ASAP) Collaborative Research Network, Chevy Chase, MD, 20815

Liv

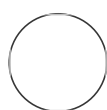

OLIVIA HARDING

ABSTRACT

We altered a previously developed method (dx.doi.org/10.17504/protocols.io.bxpfpmjn) in order to investigate the interaction of endogenous NEMO with damaged mitochondria in Hippocampal rat neurons. Mitochondrial damage was induced by the administration of 35 nM Antimycin A, a complex III inhibitor, over 1 hr. Here, we visualized neuronal mitochondria with a mito-targeted construct, Mito-SNAP, and a fluorescent SNAP ligand, JaneliaFluor 646, though other mitochondrial markers may be suitable. Endogenous NEMO was visualized with the commercially available anti-NEMO primary antibody (abcam). AntA-treated neurons exhibited fragmented mitochondria. We did not observe a change in the appearance of NEMO labeling in AntA treated cells versus non-treated cells.

MATERIALS

BSA block: 0.2% Triton-X/3% BSA in PBS.

BEFORE START INSTRUCTIONS

Prepare poly-L-lysine (PLL)-coated surfaces on which to plate hippocampal neurons by pipetting 700-1000 uL 0.5 mg/mL PLL onto glass-bottomed dishes, prioritizing coverage of the glass imaging window. Neurons may also be fixed on glass coverslips.

Plating and maintenance of hippocampal neurons

6d

- 1 Plate primary hippocampal Sprague Dawley rat neurons (embryonic day 18) on 35 mm glass-bottomed, PLL-coated dishes at 250,00 cells per dish in MEM supplemented with 10% horse serum, 33 mM D-glucose, and 1 mM sodium pyruvate. This is day 0.

Note

Information for dissection and separation of hippocampal neurons is not included because we outsource these steps to an on-campus core facility.

- 2 Incubate neurons for at least 2 but no more than 5 hours at 37 degrees C, 5% CO₂. 2h
- 3 Replace media with Neurobasal (Gibco) supplemented with 33 mM D-glucose, 2 mM GlutaMAX (Invitrogen), 100 units/ml penicillin, 100 ug/ml streptomycin, and 2% B-27 (ThermoFisher) (Maintenance Media; MM) and continue to maintain at 37 degrees C, 5% CO₂. 1d
- 4 Add AraC (5 uM) one day post-plating to prevent glial proliferation
- 5 Allow neurons to grow until Day 6.

Transfection of primary neurons

2d

- 6 On day 6, transfect neurons. 0.3 ug Mito-SNAP per 35-mm dish. Use 3.8 uL Lipofectamine 2000 2h

Transfection Reagent (ThermoFisher) per dish.

- 6.1** Equilibrate a volume of fresh MM equal to 2 mL per dish by incubating MM in a flask in the cell culture incubator. 45m
- 6.2** Replace conditioned MM with 1 mL fresh MM for each dish. Save conditioned MM in a fresh flask in the incubator.
- 6.3** Prepare transfection solutions A and B by adding 0.3 ug Mito-SNAP to 150 uL Neurobasal (A) and 3.8 uL Lipofectamine 2000 to 150 uL Neurobasal (B) for each dish. Rotate by hand 8 times and incubate at room temperature for 5 min. 5m
- 6.4** Briefly centrifuge A and B with a tabletop centrifuge. Then combine A and B and rotate by hand 8 times to mix. Incubate at room temperature for 20 min. 20m
- 6.5** Add A+B (300 uL) to each dish of primary neurons dropwise.
- 6.6** Incubate at 37 degrees 45 min. 45m
- 6.7** Replace transfection media with 50%/50% fresh and conditioned media.
- 7** Incubate transfected neurons 36-48 hours. 2d

Fluorescent labeling in neurons

1h

8 Incubate neurons in 2 μ M SNAP ligand (Janelia-Fluor 646) for 30 min.

30m

9 Wash 2X, then rest 30 min in conditioned media.

30m

Mitophagy induction in neurons

1h

10 Replace neuron media with MM containing 35 nM Antimycin A (or, for control conditions, media containing no small molecule inhibitor) and incubate for 1 hr.

1h

10.1 During this period, prepare fixation reagent (see next step)

Fix, permeabilize, and block neurons

35m

11 Prepare 4% paraformaldehyde (PFA)/4% sucrose in PBS

11.1 Thaw a 10 mL aliquot of 4% PFA in PBS using a 37 degree water bath

7m

11.2 Add 0.4 g sucrose to aliquot

- 11.3** Rotate at R.T. for 10 min to dissolve 10m
- 11.4** Warm 4% PFA/4% sucrose to 37 degrees 10m
- 12** Gently add at least 0.5 mL fixative to each dish and incubate 10 min at 37 degrees 10m
- 13** Aspirate PFA and wash samples with PBS 4X 5 min 20m
- 14** Permeabilize cells by adding at least 300 uL 0.2% Triton-X in PBS for 5 min 5m
- 15** Block cells with at least 300 uL BSA block (see Materials) for 50 min
- 15.1** Prepare primary antibody solution during this step.

Perform immunolabeling and Hoechst stain

2d

- 16** Prepare primary solution with a volume of 200 uL * number of dishes. Use anti-NEMO (abcam, raised in rabbit) at a 1:500 dilution in BSA block.

- 16.1** Centrifuge primary antibody at 4 degrees, max speed, for 4 min before use. After centrifugation, pipet from top of antibody aliquot.
- 17** Aspirate blocking solution from each sample and gently add 180 uL primary solution.
- 18** Incubate samples in primary solution overnight in a dark environment. 1d
- 19** The following day, aspirate primary solution and add 100-200 uL PBS to wash for 5 min. Carry out PBS wash for a total of 4 times. 20m
- 20** Prepare secondary solution with a volume of 110 uL * number of dishes. Use anti-rabbit 488 secondary antibody at 1:200 dilution in BSA block.
- 20.1** Centrifuge primary antibody at 4 degrees at max speed for 4 min before use. After centrifugation, pipet from top of antibody aliquot.
- 21** Aspirate PBS from final wash and add 100 uL secondary solution to each sample, making sure to submerge all areas of the glass window. Incubate 1 hr at room temperature. 1h
- 22** Aspirate secondary solution and add 100-200 uL PBS to wash for 5 min. Carry out PBS wash for a total of 2 times. 10m

- 23** Prepare Hoechst solution with a volume of 110 uL * number of dishes. Use Hoechst dye at 1:1000 dilution in PBS.
- 24** Aspirate PBS from second wash and add Hoechst solution. Incubate 5 min at room temperature. 5m
- 25** Repeat step 22. 10m
- 26** Aspirate final PBS wash and gently add enough volume of Vectashield to cover surface of glass window.

Imaging

- 27** Image samples, for example, with a Nikon Eclipse Ti Microscope with a 100X objective (Apochromat, 1.49-N.A. oil immersion) and an UltraView Vox spinning disk confocal system (PerkinElmer). Collect Z-stacks at 0.15 nm/step with Volocity acquisition software (PerkinElmer).