



DEC 27, 2022

## OPEN ACCESS

### DOI:

[dx.doi.org/10.17504/protocols.io.dm6gpjdpqzp/v1](https://dx.doi.org/10.17504/protocols.io.dm6gpjdpqzp/v1)

**Protocol Citation:** Sarah M Prostack, Edgar M Medina, Erik Kalinka, Lillian Fritz-Laylin 2022. Protocol 6: Selecting for *Spizellomyces punctatus* transformants. **protocols.io** <https://dx.doi.org/10.17504/protocols.io.dm6gpjdpqzp/v1>

**License:** This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](#), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

**Protocol status:** Working

**Created:** Nov 03, 2022

**Last Modified:** Dec 27, 2022

**PROTOCOL integer ID:** 72237


**Keywords:** *Spizellomyces punctatus*, plasmid

# Protocol 6: Selecting for *Spizellomyces punctatus* transformants

In 1 collection

Sarah M Prostack<sup>1</sup>, Edgar M Medina<sup>1</sup>, Erik Kalinka<sup>1</sup>, Lillian Fritz-Laylin<sup>1</sup>

<sup>1</sup>University of Massachusetts at Amherst

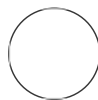

fritzlaylinlab.umass

## ABSTRACT

Once the co-culture IM plates have grown for four days, it is time to select Sp transformants. Selecting these transformants requires harvesting the cells from the co-culture plates, dislodging and removing any remaining Agro from the culture, and plating Sp onto K1 selection media. In our laboratory, we use 300 mg/L Hygromycin B (Gibco #10687010) as our selection marker. Additionally, we add 50 mg/L Carbenicillin (Carb) and 50 mg/L Tetracycline (Tet) to our K1 plates to kill *Agrobacterium*. All steps here, with the exception of centrifugation steps, must be carried out in a sterile environment, ideally in the sterile area around an open flame.

## ATTACHMENTS

[Spizellomyces transformation\\_steps.pdf](#)

## MATERIALS TEXT

### Materials

Co-culture of Agro and Sp after 4 days of growth (see *Protocol 5: Agrobacterium-mediated transformation of Spizellomyces punctatus*)


Single edge razor blades, sterile (such as Personna .009 RD #94-120-71)

K1 plates with and without selection antimicrobials


Forceps (such as AESCULAP BD232)

DS solution, sterile (see recipe)


50 mL conical tubes, sterile (such as

 Centrifuge Tubes-Bag 50mL Centrifuge Tube - Bag Sterile **Cell Treat Scientific Products Catalog #229421**

)  
100-1,000 µL micropipette such as

 Eppendorf Research Plus Single Channel pipette 100-1000 uL blue operating button for use with 100 **pipette.com Catalog #3123000063 ES-1000**

20-200 µL micropipette such as

 Eppendorf Research Plus single channel pipette 2-20 uL yellow operating button for use with 2 **pipette.com Catalog #3123000039**

Filter tips for the micropipettes, sterile such as

 TIPONE® FILTER TIPS **USA Scientific Catalog #1122-1830**

Vortexer such as

Equipment	
<b>Variable Speed Mini Vortex Mixer</b>	NAME
Vortexer	TYPE
Fisherbrand™	BRAND
14-955-163	SKU
<a href="https://www.fishersci.com/shop/products/fisher-scientific-variable-speed-mini-vortex-mixer/14955163">https://www.fishersci.com/shop/products/fisher-scientific-variable-speed-mini-vortex-mixer/14955163</a>	LINK

Centrifuge capable of holding 50 mL conicals (such

Equipment	
<b>Eppendorf™ 5810R Centrifuge</b>	NAME
Centrifuge	TYPE
Eppendorf	BRAND
02-262-8187	SKU
<a href="https://www.fishersci.com/shop/products/eppendorf-5810r-centrifuge-rotor-packages-16/022628187">https://www.fishersci.com/shop/products/eppendorf-5810r-centrifuge-rotor-packages-16/022628187</a>	LINK

)  
 Plastic storage containers, or equivalent (e.x. plastic shoe containers)  
 Open flame source

#### SAFETY WARNINGS


- ⚠ This protocol involves the use of single-edge razor blades. It is extremely important that proper safety protocols regarding sharps be taken– protective eye equipment and a sharps container must be used while conducting this protocol and disposing of razors.


#### ATTACHMENTS

[Spizellomyces transformation steps.pdf](#)



## Steps

15m

1 For each plasmid transformed, prepare a 50 mL conical with  30 mL of DS solution.

2 From the appropriate conical for the plate to be harvested, remove  1 mL of DS.

3 Distribute the DS among the four quadrants of the corresponding IM plate.

4 Incubate at  Room temperature for  00:05:00 .

5m



5 Using forceps, sterilize both sides of a single edge razor blade with the flame.

6 Tilt the rehydrated IM plate and gently scrape the surface of the agar and collect the growth into the DS pooled at the bottom.

7 Rotate the plate and scrape along the agar until most of the opaque areas on the plate are gone.

8 Use a 1000 µL pipette tip to scrap any remaining growth from the razor blade and resuspend the growth to the appropriate conical tube filled with DS.


#### Note

Once the growth is taken off the razor blade, the blade should be thoroughly sterilized using the flame and then discarded into a sharps container.

9 Aspirate all of the liquid from the scraped plate and slowly dispense the liquid back into the appropriate conical tube.

#### Note

A density gradient between the heavier cell suspension and the lighter DS solution will form. Maintaining this gradient at this point by placing retrieved co-culture in the bottom of the tube. This can help you continuously retrieve “clean” DS to continue harvesting the plate more easily.

10 From the top of the density gradient, remove  1 mL of DS.

11 Wash the surface of the scraped IM plate several times with this fresh DS.



12 Return the volume back to the conical.

13 Invert the conical 3 times.

14 Vortex the conical for 1-2 seconds to dislodge any remaining *Agro* from *Sp* cells.

15 Repeat steps 2-14 for each plasmid transformed (i.e, for each IM plate you have)

16 Centrifuge the conical(s) for  2000 rcf, Room temperature, 00:10:00 .

10m



17 Check for the presence of a pinkish pellet, this is the chytrid *Sp*.


18 Gently pour off the supernatant with the high side of the pellet facing up.






#### Note





This step can result in a great loss of *Sp* cells, thus pouring with this orientation will reduce the amount of cells lost.

19 Gently resuspend the pellet into  500  $\mu$ L of fresh DS.

20 Pipette  200  $\mu$ L of the resuspended *Sp* onto a K1 plate with selection antibiotics.




#### Note

We use  300 undetermined Hygromycin B (  Hygromycin B (50 mg/mL) **Thermo Fisher Catalog #10687010** ) and  50 undetermined Carb and  50 undetermined Tet to prevent further *Agro* growth.

21 Spread the cells using 4-5 sterile glass beads.

22 Once the plate is dry, remove the glass beads.

23 Seal and incubate plates in a humidity chamber at  Room temperature for 4 days or until colonies appear on the plate.

