




DEC 14, 2023

🌐 Open Field Protocol

Michael Lee¹

¹University of Minnesota

Team Lee

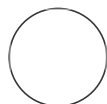

jbalster

DISCLAIMER

OPEN  ACCESS



DOI:
dx.doi.org/10.17504/protocols.io.36wgq7qjovk5/v1

Protocol Citation: Michael Lee 2023. Open Field Protocol. **protocols.io**
<https://dx.doi.org/10.17504/protocols.io.36wgq7qjovk5/v1>

License: This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](#), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

Protocol status: Working
 We use this protocol and it's working

Created: Jun 07, 2022

Last Modified: Dec 14, 2023

DISCLAIMER – FOR INFORMATIONAL PURPOSES ONLY; USE AT YOUR OWN RISK

The protocol content here is for informational purposes only and does not constitute legal, medical, clinical, or safety advice, or otherwise; content added to [protocols.io](#) is not peer reviewed and may not have undergone a formal approval of any kind. Information presented in this protocol should not substitute for independent professional judgment, advice, diagnosis, or treatment. Any action you take or refrain from taking using or relying upon the information presented here is strictly at your own risk. You agree that neither the Company nor any of the authors, contributors, administrators, or anyone else associated with [protocols.io](#), can be held responsible for your use of the information contained in or linked to this protocol or any of our Sites/Apps and Services.

ABSTRACT

This protocol details about Open Field procedure.

ATTACHMENTS

[452-957.docx](#)

GUIDELINES

Detailed records should be kept on every experiment, noting date, timing, subject test order, procedural details followed and any relevant observations pertaining to the assay (e.g. specific behavior or procedural errors).

PROTOCOL integer ID:
64066

Keywords: Open Field
Protocol, locomotor activity,
ASAPCRN

Funders
Acknowledgement:

ASAP


Set up

- 1 20 x 20 x 10 inch chamber (LxWxH).
- 2 Opaque chambers or dividers between chambers.
- 3 Light intensity = 250 lux.
- 4 Have chambers on enclosed rack to diminish outside light.

Test

30m

- 5 AnyMaze
 - 5.1 Set two zones in the chamber. Outer ~3.5 inch perimeter of box is one zone and the center of the maze is the other zone.

- 6 Record locomotor activity and track how much time the animal spends in each zone for  00:30:00 30m
- 7 Wipe down chambers with ethanol between animals.
- 8 Run males before females if running both genders on same day.
- 9 Can be used as a simple locomotion test by tracking the distance traveled and speed.
- 10 Can also be used as an anxiety test. Less anxious animals will spend more time in the center of the maze.