



NOV 29, 2023

selSeq: A method for the enrichment of non-polyadenylated RNAs including enhancer and long non-coding RNAs for sequencing

Jason D Limberis¹, John Joel Ernst¹, Metcalfe¹, Alina Nalyvayko¹

¹University of California, San Francisco



Jason D Limberis
University of California, San Francisco

OPEN ACCESS



DOI:
dx.doi.org/10.17504/protocols.io.j8nlkwpk6l5r/v1

Protocol Citation: Jason D Limberis, Joel Ernst, John Metcalfe, Alina Nalyvayko 2023. selSeq: A method for the enrichment of non-polyadenylated RNAs including enhancer and long non-coding RNAs for sequencing. **protocols.io** <https://dx.doi.org/10.17504/protocols.io.j8nlkwpk6l5r/v1>

License: This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](#), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

Protocol status: Working
We use this protocol and it's working

ABSTRACT

Non-polyadenylated RNA includes a large subset of crucial regulators of RNA expression and constitutes a substantial portion of the transcriptome, playing essential roles in gene regulation. For example, enhancer RNAs are long non-coding RNAs that perform enhancer-like functions, are bi-directionally transcribed, and usually lack polyA tails. This paper presents a novel method, *seSeq*, that selectively removes mRNA and pre-mRNA from samples to enable the selective sequencing of crucial regulatory elements, including non-polyadenylated RNA such as long non-coding RNA, enhancer RNA, and non-canonical mRNA.

Created: Mar 28, 2023

Last Modified: Nov 29, 2023

PROTOCOL integer ID: 79597

Funders

Acknowledgement:

Program for Breakthrough
Biomedical Research (PBBR)
Grant ID: NA

MATERIALS

Required

SuperScript® III First-Strand Synthesis System Thermo
Scientific Catalog #18080-051

RNase H - 1,250 units New England
Biolabs Catalog #M0297L

TURBO DNase 2 U/uL Fisher
Scientific Catalog #AM2239

Agencourt RNAClean XP Magnetic Beads Beckman Coulter Catalog #A63987

Ethanol Contributed by users

A thermocycler and a qPCR machine

A magnetic rack

Optional


Luna Universal Probe One-Step RT-qPCR Kit - 200 rxns New England
Biolabs Catalog #E3006S

Eukaryotic 18S rRNA Endogenous Control (FAM[®]/MGB probe, non-primer
limited) Thermo Fisher Catalog #4333760F

TaqMan[®]; GAPDH Control Reagents (human) Thermo
Fisher Catalog #402869

rRNA depletion oligos

BEFORE START INSTRUCTIONS




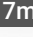

Prewarm SuperScript III 10X Buffer to  Room temperature

poly-A tailed cDNA synthesis

1 Mix the following in a 0.2ml tube

A	B
Component	Volume (µl)
Total RNA	1
Oligo dTs	1.5
10 mM dNTP mix	1.5
Nuclease-free H ₂ O	10

poly-A tailed cDNA reaction synthesis components






2 Denature sample RNA/primer mixture for  00:05:00 at  65 °C then cool to  4 °C for ≥  7m
 00:02:00

3 Spin tube briefly and add the following and mix by pipetting

55m

A	B
Component	Volume (μl)
10X SuperScript III Buffer	2
25mM MgCl ₂	4
0.1M DTT	2
Superscript III Reverse Transcriptase	2





poly-A tailed cDNA reaction synthesis components

Incubate  50 °C for  00:50:00 followed by  00:05:00 at  85 °C to deactivate the enzyme, then cool to  4 °C and proceed to the next step

Optional: rRNA depletion





4 Add in the appropriate rRNA depletion oligos for you sample

2m

Incubate  90 °C for  00:02:00 and ramp down to  Room temperature at  0.1 °C per second then proceed to the next step

poly-A tailed (and ribosomal) RNA depletion

5 Add  2 μL of RNase H

6 Incubate  37 °C for  00:20:00 followed by  00:05:00 at  65 °C to deactivate the enzyme, then cool it to  4 °C and proceed to the next step

25m

poly-A tailed (and ribosomal) DNA depletion

7 Add in the following components and mix gently by pipetting

A	B
Component	Volume (μl)
10X Turbo DNase Buffer	4
Turbo DNase	4
Nuclease-free H ₂ O	10



DNase treatment components

8 Incubate at  37 °C for  00:30:00

30m




Bead cleanup

9 Add 90 μl (1.8X) of resuspended RNAClean XP Beads to the sample
Mix by pipetting 10x

10 Incubate  00:15:00 at  On ice








15m

11 Place on the magnet, allow the beads to aggregate, and remove and discard the supernatant

12 Add  200 μL  80 % (v/v) ethanol and incubate (still on the magnet) for  00:00:30

30s

12.1 Remove the supernatant

- 12.2 Repeat  go to step #12 for a total of 2 washes
- 13 Air dry for  00:00:30, don't allow the beads to become cracked 30s
- 14 Remove the tubes from the magnetic rack
Add  50 μL H₂O (optionally add-in  1 μL RNase inhibitor) and resuspend the beads by pipetting $\geq 10\times$
- 15 Incubate  00:05:00 at  Room temperature 5m
- 16 Place on the magnet, aspirate  50 μL of the eluant into a new tube

Optional: One-step RT-qPCR quantification

17

A	B
Component	Volume (μl)
Luna Universal Probe One-Step Reaction Mix (2X)	5
Luna WarmStart RT Enzyme Mix (20X)	0.5
TaqMan GAPDH Control Reagents (human; 20x)	0.5
TaqMan 18S rRNA Control Reagents (eukaryotic; 20x)	0.5
RNA	2
Nuclease-free H ₂ O	1.5

A	B	C	D	E
Step	Temp (C)	Time (s)	Cycles	Ramp Rate (C/s)
Reverse transcription	55	600	1	2.73
Denaturation	95	60	45	2.73
Denaturation	95	10		2.73
Amplification	60	30		2.11
Capture	60	0		–

Cycle parameters for QuantStudio 3