



JUL 27, 2023

## WU sn-prep Protocol for solid tumors- joint snRNA+ATAC v2.9

Wagma

Reyka Jayasinghe<sup>1</sup>, Caravan<sup>1</sup>,


Andrew Houston<sup>1</sup>,

Nataly Naser Al

Deen<sup>1</sup>

<sup>1</sup>Department of Medicine, Washington University in St. Louis, St. Louis, MO 63110, USA

OPEN  ACCESS

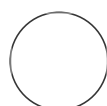

Andrew Houston

### ABSTRACT

Nuclei dissociation protocol adapted from WU sn-prep Protocol for Solid Tumors - snRNA protocol v2.8 for simultaneous profiling of genetic expression (snRNA) and chromatin accessibility (snATAC)



#### DOI:

[dx.doi.org/10.17504/protocols.io.261gednx7v47/v1](https://dx.doi.org/10.17504/protocols.io.261gednx7v47/v1)

**Protocol Citation:** Reyka Jayasinghe, Wagma Caravan, Andrew Houston, Nataly Naser Al Deen 2023. WU sn-prep Protocol for solid tumors- joint snRNA+ATAC v2.9. **protocols.io** <https://dx.doi.org/10.17504/protocols.io.261gednx7v47/v1>

**License:** This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](#), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

**Protocol status:** Working  
We use this protocol and it's working

**Created:** Jul 26, 2023

## Reagents and Tools

- 1 1x Lysis buffer (2mL):  
10mM Tris-HCl (pH 7.4) (Thermo; 15567027), 20μL  
10mM NaCl (Thermo; AM9759), 4μL  
3 mM MgCl<sub>2</sub> (Thermo; AM9530G), 6μL  
NP-40 substitute (Sigma, 74385-1L), 2μL  
1 M DTT (Sigma, 646563), 2μL  
Nuclease Free Water (Invitrogen, AM9937), 1.966mL
- 2 Lysis Dilution Buffer (10mL):  
10mM Tris-HCl (pH 7.4) (Thermo; 15567027), 100μL  
10mM NaCl (Thermo; AM9759), 20μL  
3 mM MgCl<sub>2</sub> (Thermo; AM9530G), 30μL  
1 M DTT (Sigma, 646563), 10μL  
Nuclease Free Water (Invitrogen, AM9937), 9.840mL
- 3 0.1x Lysis Buffer (10mL):  
1x Lysis Buffer (1mL) + Lysis Dilution Buffer (9mL)
- 4 Wash and Resuspension buffer (WR buffer, 10mL):  
1X PBS, 8mL  
10% Stock BSA Solution (MACS, 130-091-376), 2mL
- 5 Roche Protector RNase inhibitor-2000 U/μL (Millipore Sigma 3335399001)
- 6 Trypan blue (2X) - filtered at 0.22μm
- 7 7-AAD (7-Aminoactinomycin D) (Millipore Sigma SML1633-1ML)

- 8 Glass homogenizers (Fisher: 2mL tube- K8853030002, Small clearance pestle- K8853020002, Large clearance pestle- K8853010002)
- 9 Fine forceps and scalpels

## General Notes

- 10
- Keep everything on ice (or in the cold room)
  - Use RNase free reagents and consumables, before starting wipe down all surfaces/pipettes with RNase away and 70% ethanol
  - Avoid foam and bubbles as much as possible with gentle strokes and pipetting slowly

## Nuclei Dissociation

- 11 If using frozen tissue sample, use a scalpel (aided by a pair of fine forceps) to cut the cold samples (25-35mg) into 2mm pieces, add 1 ml of ice-cold 0.1X lysis buffer and 30  $\mu$ L RNase Inhibitor, load into the glass homogenizer. Homogenize with 4-6 push/pulls using the pestle, incubate on ice for 1 min with an additional 500-1000 $\mu$ L of 0.1X cold lysis buffer. Pipette gently for 4 times. Incubate on ice again for up to 1 min.
- 11.1 If using pulverized powder, start with 15-35 mg total, add 1 ml of ice-cold 0.1X lysis buffer and 30  $\mu$ L RNase Inhibitor, pipette the powder/lysis buffer mix gently for 6-8 times. Let sit on ice for 30". Pipette another 4-6 times with an additional 500-1000 $\mu$ L of 0.1X cold lysis buffer. Incubate on ice again for 1' – could be reduced (to like 20-45"). Add to the glass homogenizer, use the pestle for 3-4 push/pulls.
- 11.2 If using OCT sections, start with 300-450 $\mu$ m total sectioned into a 1.5mL tube. Add 1 ml of ice-cold 0.1X lysis buffer and 30  $\mu$ L RNase Inhibitor, pipette the mix gently for 10-12 times, will be sticky as the OCT thaws. Let sit on ice for 30". Pipette another 4-6 times with an additional 500-1000 $\mu$ L of 0.1X cold lysis buffer. Incubate on ice again for 1' – could be reduced (to like 20-45"). Add to the glass homogenizer, dounce with the pestle for 6-8 push/pulls.
- 12 Filter the homogenate through a 40 $\mu$ M cell strainer on ice on top of a 50ml conical tube. Wash the filter with 1ml WR buffer. Collect this into the same tube, the total filtrate is ~3 ml.

If there is still tissue on the strainer, backwash with 2 mL 0.1x Lysis buffer, follow previous steps again to dissociate completely. If going to FACS with the backwash, proceed with this sample as if it were a different tissue but sort into same collection tube.

- 13** Transfer the filtrate to a 5ml Eppendorf tube. Centrifuge at 500g for 6' at 4°C, resuspend with 100-400µL WR buffer (depending on pellet size) and 10 µL RNase inhibitor.
- 14** Transfer into a FACS sorting tube and add 3µL 7-AAD per 500 µL of resuspended sample, incubate for 10 minutes before sorting at FACS. After resuspending sample in wash buffer, if small chunks are still visible (after ~3 minutes of resuspension) use 40µM mini-strainer over FACS tube to remove chunks (proceeding to FACS with sample in current condition will clog machine and will result in additional lost sample).
- 15** Add 50µL 10% BSA solution and 10µL RNase inhibitor and into a 2mL nonbinding tube for collection. Sort 300-400K of the nuclei into the collection tube.
- 16** After FACS sorting is done centrifuge the 2ml collection tube at 500g for 6' at 4°C. There will not be a visible pellet, remove all the supernatant (likely ~3µL will be left).
- 17** Resuspend in 3µL 2x Nuclei Dilution Buffer (prepared as directed by 10x Genomics, pg. 28) and 1µL RNase inhibitor, final volume should be ~10µL.
- 18** Quantify nuclei quality and quantity using a hemacytometer or Countess II utilizing Trypan blue.
- 19** Load desired concentration of nuclei, proceed with the protocol as outlined by 10x Genomics which can be found [here](#).