



JUL 12, 2023

🌐 JAX SCBL workflow for scRNA-seq of human placenta


William F

Diane Luo¹, Jessica Grassmann¹, Paul Gabriel¹, Flynn¹,
Elise Courtois¹, Paul Robson¹

¹The Jackson Laboratory for Genomic Medicine

JAX Single Cell Biology

Tech. support email: scblservice@jax.org

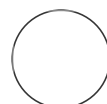

William F Flynn

OPEN  ACCESS

DOI:

dx.doi.org/10.17504/protocols.io.eq2lynw6pvx9/v1

Protocol Citation: Diane Luo, Jessica Grassmann, Paul Gabriel, William F Flynn, Elise Courtois, Paul Robson 2023. JAX SCBL workflow for scRNA-seq of human placenta.

protocols.io

<https://dx.doi.org/10.17504/protocols.io.eq2lynw6pvx9/v1>

License: This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](#), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

Protocol status: Working
We use this protocol and it's working

Created: Mar 16, 2022

Last Modified: Jul 12, 2023

PROTOCOL integer ID:
59523

ABSTRACT

Described here is the workflow used by the JAX Single Cell Biology lab to generate single cell transcriptomic libraries from human placenta samples collected at the UCSD HuBMAP Female Reproductive Tissue Mapping Center.

MATERIALS

See accompanying protocols.

Tissue collection

- 1 Tissue is collected and processed as soon as possible after Cesarean section or vaginal delivery according to the following protocol:

[Human Placenta Tissue Collection and Preservation Methods - UCSD Female Reproductive TMC](#)

Tissue is then stored in MACS storage buffer for overnight shipping.

Sample preparation

- 2 Upon receiving the sample, single cells are isolated through a digestion with cold active protease as described here:

[Human endometrium and endometriosis tissue dissociation for single-cell RNA sequencing](#)

- 2.1 The viability and cell concentration for each sample is assessed via a Countess FL II or LUNA FX7 automated cell counter.

- 2.2 The single cell suspension is then loaded onto one lane of the appropriate 10x Genomics gene expression chip as follows:

| A | B | C | D |
|--------|---------------------|------------------------|----------------------------|
| Chip | Approx Cells Loaded | Approx. Cells Targeted | Compatible Instrument |
| Chip G | 12,000 | 6,000 | Chromium Controller, Xi, X |
| Chip M | 40,000 | 20,000 | Chromium X |


Library Generation

- 3 The resulting loading, running, and library construction is performed according to the appropriate 10x Genomics protocol:

[JAX - 10x Genomics - NEXTGEM v3.1 - 3' Gene Expression](#)

Sequencing

- 4 We have typically targeted 100,000 read pairs/cell for human placenta sample, resulting in target library sizes of roughly 600M-2,000M read pairs. These are multiplexed on a NovaSeq 6000 S4



lane using v1.5 reagents and the reverse complement workflow in a 28/8/0/90 or 28/8/8/90 configuration.