

Feb 25, 2021

# NEBNext® ARTIC SARS-CoV-2 FS Library Prep Kit (Illumina®) E7658

New England Biolabs<sup>1</sup><sup>1</sup>New England Biolabs

1 Works for me [dx.doi.org/10.17504/protocols.io.bsbnak6](https://dx.doi.org/10.17504/protocols.io.bsbnak6)

New England Biolabs (NEB)

Coronavirus Method Development Community

Isabel Gautreau  
New England Biolabs

SUBMIT TO PLOS ONE

## ABSTRACT

This protocol details methods for the NEBNext® ARTIC SARS-CoV-2 FS Library Prep Kit (Illumina®), NEB #E7658S/L 24/96 reactions.

## EXTERNAL LINK

<https://www.neb.com/-/media/nebus/files/manuals/manuale7658.pdf?rev=7253d569aa3040069cb09a9aa3724f2d&hash=53C7EDF0D2FDEDAE75308379A5EB2D89>

## DOI

[dx.doi.org/10.17504/protocols.io.bsbnak6](https://dx.doi.org/10.17504/protocols.io.bsbnak6)

## EXTERNAL LINK

<https://www.neb.com/-/media/nebus/files/manuals/manuale7658.pdf?rev=7253d569aa3040069cb09a9aa3724f2d&hash=53C7EDF0D2FDEDAE75308379A5EB2D89>

## PROTOCOL CITATION

New England Biolabs 2021. NEBNext® ARTIC SARS-CoV-2 FS Library Prep Kit (Illumina®) E7658.  
**protocols.io**  
<https://dx.doi.org/10.17504/protocols.io.bsbnak6>

## KEYWORDS

NEBNext, NEB, ARTIC, SARS-CoV-2, FS, Library Prep, Illumina

## LICENSE

———— This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](https://creativecommons.org/licenses/by/4.0/), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

## CREATED

Feb 10, 2021

## LAST MODIFIED

Feb 25, 2021

## OWNERSHIP HISTORY

Feb 10, 2021 Liz Brydon Protocols.io

Feb 25, 2021 Isabel Gautreau New England Biolabs

## GUIDELINES


### Overview

The NEBNext ARTIC SARS-CoV-2 FS Library Prep Kit for Illumina contains the enzymes, buffers and oligos required to convert a broad range of total RNA into high quality, targeted, libraries for next-generation sequencing on the Illumina platform. Primers targeting the human EDF1 (NEBNext ARTIC Human Primer Mix 1) and NEDD8 (NEBNext ARTIC Human Primer Mix 2) genes are supplied as optional internal controls. The fast, user-friendly workflow also has minimal hands-on time.

Each kit component must pass rigorous quality control standards, and for each new lot the entire set of reagents is functionally validated together by construction and sequencing of indexed libraries on an Illumina sequencing platform.

For larger volume requirements, customized and bulk packaging is available by purchasing through the OEM/Bulks department at NEB. Please contact OEM@neb.com for further information.

**Figure 1. Workflow demonstrating the use of NEBNext ARTIC SARS-CoV-2 FS Library Prep Kit.**

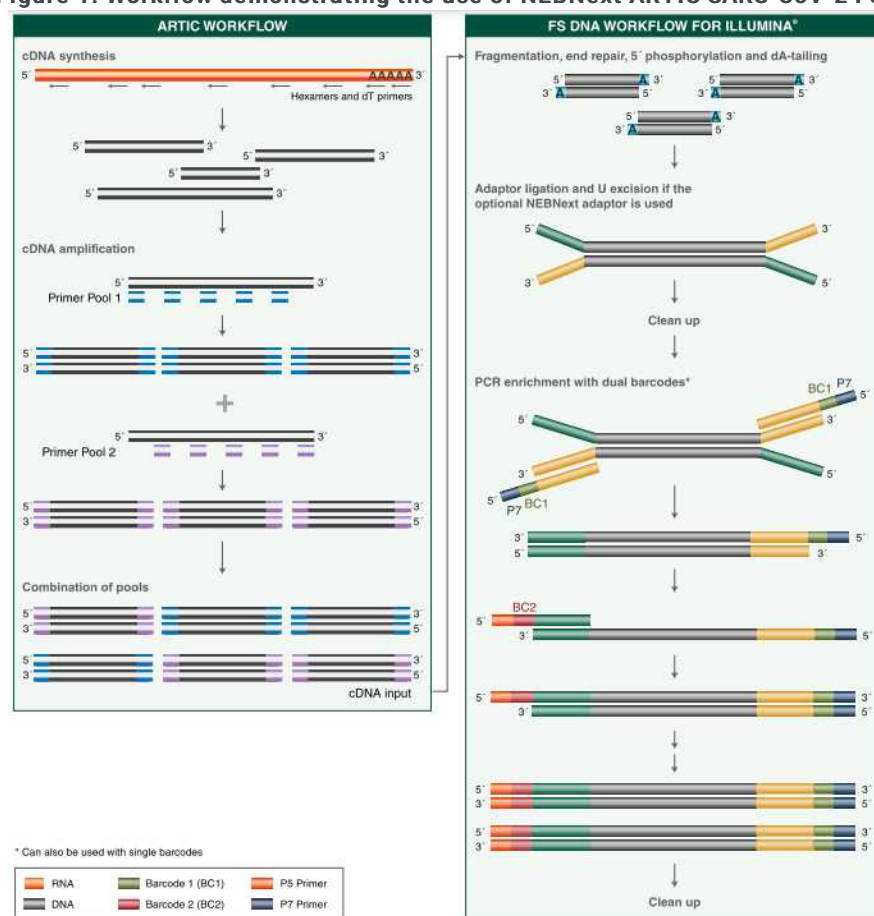

## MATERIALS TEXT

### The Library Kit Includes

The volumes provided are sufficient for preparation of up to 24 reactions (NEB #E7658S) and 96 reactions (NEB #E7658L).

### Package 1: Store at -20°C.

(lilac) LunaScript® RT SuperMix (5X)

(lilac) Q5® Hot Start High-Fidelity 2X Master Mix

(lilac) NEBNext ARTIC SARS-CoV-2 Primer Mix 1  
 (lilac) NEBNext ARTIC SARS-CoV-2 Primer Mix 2  
 (lilac) NEBNext ARTIC Human Primer Mix 1  
 (lilac) NEBNext ARTIC Human Primer Mix 2  
 (yellow) NEBNext Ultra II FS Enzyme Mix  
 (yellow) NEBNext Ultra II FS Reaction Buffer  
 (red) NEBNext Ultra II Ligation Master Mix  
 (blue) NEBNext Library PCR Master Mix  
 (white) 0.1X TE Buffer  
 (white) Nuclease-free Water

**Package 2: Store at room temperature. Do not freeze.**

NEBNext Sample Purification Beads

**Required Materials Not Included**

- NEBNext Singleplex or Multiplex Oligos for Illumina – [www.neb.com/oligos](http://www.neb.com/oligos)
- 80% Ethanol (freshly prepared)
- Nuclease-free Water
- DNA LoBind Tubes (Eppendorf® #022431021)
- Magnetic rack/stand (NEB #S1515, Alpaqua®, cat. #A001322 or equivalent)
- Thermocycler
- Vortex Mixer
- Microcentrifuge
- Agilent® Bioanalyzer® or similar fragment analyzer and associated consumables
- DNase RNase free PCR strip tubes (USA Scientific® 1402-1708)

**Kit Components**

**NEB #E7658S Table of Components**

| A      | B  | C        |
|--------|--|----------|
| NEB #  | PRODUCT                                  | VOLUME   |
| E7651A | LunaScript RT SuperMix (5X)              | 0.048 ml |
| E7652A | Q5 Hot Start High-Fidelity 2X Master Mix | 0.3 ml   |
| E7725A | NEBNext ARTIC SARS-CoV-2 Primer Mix 1    | 0.042 ml |
| E7726A | NEBNext ARTIC SARS-CoV-2 Primer Mix 2    | 0.042 ml |
| E7727A | NEBNext ARTIC Human Primer Mix 1         | 0.010 ml |
| E7728A | NEBNext ARTIC Human Primer Mix 2         | 0.010 ml |
| E7668A | NEBNext Ultra II FS Enzyme Mix           | 0.024 ml |
| E7669A | NEBNext Ultra II FS Reaction Buffer      | 0.084 ml |
| E7655A | NEBNext Ultra II Ligation Master Mix     | 0.36 ml  |
| E7656A | NEBNext Library PCR Master Mix           | 0.3 ml   |
| E7657A | 0.1X TE                                  | 1.3 ml   |
| E7667A | Nuclease free-Water                      | 1.5 ml   |
| E7659S | NEBNext Sample Purification Beads        | 2.1 ml   |

**NEB #E7658L Table of Components**

| A       | B  | C           |
|---------|--|-------------|
| NEB #   | PRODUCT                                  | VOLUME      |
| E7651AA | LunaScript RT SuperMix (5X)              | 0.192 ml    |
| E7652AA | Q5 Hot Start High-Fidelity 2X Master Mix | 1.2 ml      |
| E7725AA | NEBNext ARTIC SARS-CoV-2 Primer Mix 1    | 0.168 ml    |
| E7726AA | NEBNext ARTIC SARS-CoV-2 Primer Mix 2    | 0.168 ml    |
| E7727AA | NEBNext ARTIC Human Primer Mix 1         | 0.034 ml    |
| E7728AA | NEBNext ARTIC Human Primer Mix 2         | 0.034 ml    |
| E7668AA | NEBNext Ultra II FS Enzyme Mix           | 0.096 ml    |
| E7669AA | NEBNext Ultra II FS Reaction Buffer      | 0.336 ml    |
| E7655AA | NEBNext Ultra II Ligation Master Mix     | 2 x 0.72 ml |
| E7656AA | NEBNext Library PCR Master Mix           | 1.2 ml      |
| E7657AA | 0.1X TE                                  | 5.2 ml      |
| E7667AA | Nuclease free-Water                      | 1.5 ml      |
| E7659L  | NEBNext Sample Purification Beads        | 4 x 2.1 ml  |

#### Companion Products

| A      | B                               | C        |
|--------|---------------------------------|----------|
| NEB #  | PRODUCT                         | VOLUME   |
| T2010S | Monarch® Total RNA Miniprep Kit | 50 preps |

#### NEBNext ARTIC Human Primers

| A                                | B     | C               | D  |
|----------------------------------|-------|-----------------|--|
| PRIMER MIX                       | GENE  | POSITION        | PRIMERS                                      |
| NEBNext ARTIC Human Primer Mix 1 | EDF1  | 113 bp – 501 bp | GGCCAAATCCAAGCAGGCTA<br>GTGTTTCATTTCCGCTAGGC |
| NEBNext ARTIC Human Primer Mix 2 | NEDD8 | 110 bp – 489 bp | AAAGTGAAGACGCTGACCGG<br>GGGATCCTCACAGTCTCCCA |

Detailed information for the ARTIC Human control primers can be found at: <https://doi.org/10.5281/zenodo.4495958>

#### NEBNext ARTIC SARS-CoV-2 Primers

NEBNext ARTIC SARS-CoV-2 Primers for SARS-CoV-2 genome amplification are based on hCoV-2019/nCoV-2019 Version 3 (v3) sequences with balanced primer concentrations. Sequence information can be found at: [https://github.com/joshquick/artic-ncov2019/blob/master/primer\\_schemes/nCoV-2019/V3/nCoV-2019.tsv](https://github.com/joshquick/artic-ncov2019/blob/master/primer_schemes/nCoV-2019/V3/nCoV-2019.tsv)

#### SAFETY WARNINGS

Please refer to Safety Data Sheets (SDS) for health and environmental hazards.

This product is intended for research purposes only. This product is not intended to be used for therapeutic or diagnostic purposes in humans or animals.

#### BEFORE STARTING

**Note:** The amount of RNA required for detection depends on the abundance of the RNA of interest. In general, we recommend using > 10 copies of the (SARS-CoV-2) viral genome as input, however, results may vary depending on the quality of the input. In addition, we recommend setting up a no template control reaction and that reactions are set-up in a hood.

The presence of carry-over products can interfere with sequencing accuracy, particularly for low copy targets. Therefore, it is important to carry out the appropriate no template control (NTC) reactions to demonstrate that

positive reactions are meaningful.

## cDNA Synthesis

1 

Gently mix and spin down the LunaScript RT SuperMix reagent. Prepare the cDNA synthesis reaction as described below:

| A                              | B      |
|--------------------------------|--------|
| COMPONENT                      | VOLUME |
| RNA Sample                     | 8 µl   |
| (lilac) LunaScript RT SuperMix | 2 µl   |
| <i>Total Volume</i>            | 10 µl  |

For no template controls, mix the following components:

| A                              | B      |
|--------------------------------|--------|
| COMPONENT                      | VOLUME |
| (white) Nuclease-free Water    | 8 µl   |
| (lilac) LunaScript RT SuperMix | 2 µl   |
| <i>Total Volume</i>            | 10 µl  |

2  

Incubate reactions in a thermocycler with the following steps:

| A                 | B    | C          | D      |
|-------------------|------|------------|--------|
| CYCLE STEP        | TEMP | TIME       | CYCLES |
| Primer Annealing  | 25°C | 2 minutes  | 1      |
| cDNA Synthesis    | 55°C | 20 minutes |        |
| Heat Inactivation | 95°C | 1 minute   |        |
| Hold              | 4°C  | ∞          |        |


\*Set heated lid to 105°C

Samples can be stored at  -20 °C for up to a week.

3 Proceed to cDNA Amplification (Next Step).

## cDNA Amplification

4

 **4.5 µl cDNA input** is recommended. If using less than 4.5 µl of cDNA, add nuclease-free water to a final volume of 4.5 µl. We recommend setting up the cDNA synthesis and cDNA amplification reactions in different rooms to minimize cross-contamination of future reactions.

Use of the NEBNext ARTIC Human Primer Mix 1 and 2 are optional. If used, the appropriate ARTIC Human Primer Mix and ARTIC SARS-CoV-2 Primer Mix should be combined prior to use. More specifically, ARTIC Human Primer Mix 1 should be combined with ARTIC SARS-CoV-2 Primer Mix 1 and ARTIC Human Primer Mix 2 with ARTIC SARS-CoV-2 Primer Mix 2. Mixing directions are listed below.

## 5

Gently mix and spin down reagents. Prepare the split pool cDNA amplification reactions as described below:

### For Pool Set A:

If using the ARTIC Human Primer Mix and a *24 reaction kit*, combine **0.7 µl ARTIC Human Primer Mix 1** with **42 µl ARTIC SARS-CoV-2 Primer Mix 1**, vortex and spin down reagents. If using a *96 reaction kit*, combine **2.8 µl ARTIC Primer Mix 1** with **168 µl ARTIC SARS-CoV-2 Primer Mix 1**, vortex and spin down reagents.

Use **1.75 µl combined mix** for each Pool Set A reaction.

| A  | B       |
|--|---------|
| COMPONENT  | VOLUME  |
| cDNA (Step 2)                                    | 4.5 µl  |
| (lilac) Q5 Hot Start High-Fidelity 2X Master Mix | 6.25 µl |
| ARTIC SARS-CoV-2 Primer Mix 1 *                  | 1.75 µl |
| <i>Total Volume</i>                              | 12.5 µl |

\* If using Human Primer Mix 1, add 1.75 µl of the combined ARTIC SARS-CoV-2 Primer Mix 1 and ARTIC Human Primer Mix 1.

### For Pool Set B:

If using the ARTIC Human Primer Mix and a *24 reaction kit*, combine **0.7 µl ARTIC Human Primer Mix 2** with **42 µl ARTIC SARS-CoV-2 Primer Mix 2**, vortex and spin down reagents. If using *96 reaction kit*, combine **2.8 µl ARTIC Human Primer Mix 2** with **168 µl ARTIC SARS-CoV-2 Primer Mix 2**.

Use **1.75 µl combined mix** for each Pool Set B reaction.

| A  | B       |
|--|---------|
| COMPONENT  | VOLUME  |
| cDNA (Step 2)                                    | 4.5 µl  |
| (lilac) Q5 Hot Start High-Fidelity 2X Master Mix | 6.25 µl |
| ARTIC SARS-CoV-2 Primer Mix 2 *                  | 1.75 µl |
| <i>Total Volume</i>                              | 12.5 µl |

\* If using Human Primer Mix 2, add 1.75 µl of the combined ARTIC SARS-CoV-2 Primer Mix 2 and ARTIC Human Primer Mix 2.

## 6

Incubate reactions in a thermocycler with the following steps:

| A                    | B    | C          | D      |
|----------------------|------|------------|--------|
| CYCLE STEP           | TEMP | TIME       | CYCLES |
| Initial Denaturation | 98°C | 30 seconds | 1      |
| Denature             | 95°C | 15 seconds | 35     |
| Annealing/Extension  | 63°C | 5 minutes  |        |
| Hold                 | 4°C  | ∞          | 1      |

\*Set heated lid to 105°C

7 

Combine the paired Pool A and Pool B PCR reactions for each sample.

Samples can be stored at  -20 °C for up to a week.

8 Proceed to cDNA Amplicon cleanup (Next Step).


#### Cleanup of cDNA Amplicons

9

The volume of NEBNext Sample Purification Beads provided here are for use with the sample composition at this step (25 µl; Step 7). These bead volumes may not work properly for a cleanup at a different step in the workflow, or if this is a second cleanup at this step. For cleanups of samples contained in different buffer conditions, the volumes may need to be experimentally determined.

10 Vortex NEBNext Sample Purification Beads to resuspend.

11 

Add  20 µl (0.8X) resuspended beads to the combined PCR reaction. Mix well by pipetting up and down at least 10 times. Be careful to expel all of the liquid out of the tip during the last mix. Vortexing for 3-5 seconds on high can also be used. If centrifuging samples after mixing, be sure to stop the centrifugation before the beads start to settle out.

12 

5m

Incubate samples at  Room temperature for at least  00:05:00 .




13 Place the tube/plate on an appropriate magnetic stand to separate the beads from the supernatant. If necessary, quickly spin the sample to collect the liquid from the sides of the tube or plate wells before placing on the magnetic stand.

- 14 After 5 minutes (or when the solution is clear), carefully remove and discard the supernatant. Be careful not to disturb the beads that contain DNA targets.

**Caution: do not discard the beads.**

- 15 




30s

Add  **200 µl 80% freshly prepared ethanol** to the tube/plate while in the magnetic stand. Incubate at  **Room temperature** for  **00:00:30**, and then carefully remove and discard the supernatant. Be careful not to disturb the beads that contain DNA targets.

- 16 

30s

Repeat previous Step once for a total of two washes:



Add  **200 µl 80% freshly prepared ethanol** to the tube/plate while in the magnetic stand. Incubate at  **Room temperature** for  **00:00:30**, and then carefully remove and discard the supernatant. Be careful not to disturb the beads that contain DNA targets.

**Be sure to remove all visible liquid after the second wash. If necessary, briefly spin the tube/plate, place back on the magnet and remove traces of ethanol with a p10 pipette tip.**

- 17 Air dry the beads for up to 5 minutes while the tube/plate is on the magnetic stand with the lid open.



**Caution: Do not over-dry the beads. This may result in lower recovery of DNA. Elute the samples when the beads are still dark brown and glossy looking, but when all visible liquid has evaporated. When the beads turn lighter brown and start to crack, they are too dry.**

- 18 

Remove the tube/plate from the magnetic stand. Elute the DNA target from the beads by adding  **16 µl 0.1X TE**. If not assessing amplicons (Step 21) elute DNA in  **15 µl 0.1X TE**.

- 19 

2m

Mix well by pipetting up and down 10 times, or on a vortex mixer. Incubate for at least  **00:02:00** at  **Room temperature**. If necessary, quickly spin the sample to collect the liquid from the sides of the tube or plate wells before placing back on the magnetic stand.

- 20 

Place the tube/plate on the magnetic stand. After 5 minutes (or when the solution is clear), transfer  **14 µl** to a new



PCR tube. If not assessing amplicons (Step 21) transfer **13 µl** to a new PCR tube.

## 21

We recommend assessing cDNA amplicon (from Step 20) concentrations with a Qubit fluorometer.

**Amplicons may also be run on a Bioanalyzer or TapeStation® to confirm 400 bp size of amplicons. To run on a TapeStation, dilute amplicon 10-fold with 0.1X TE Buffer and run 2 µl on a DNA High Sensitivity ScreenTape. (See Figure 21 below for example of amplicon size profile on a Bioanalyzer).**

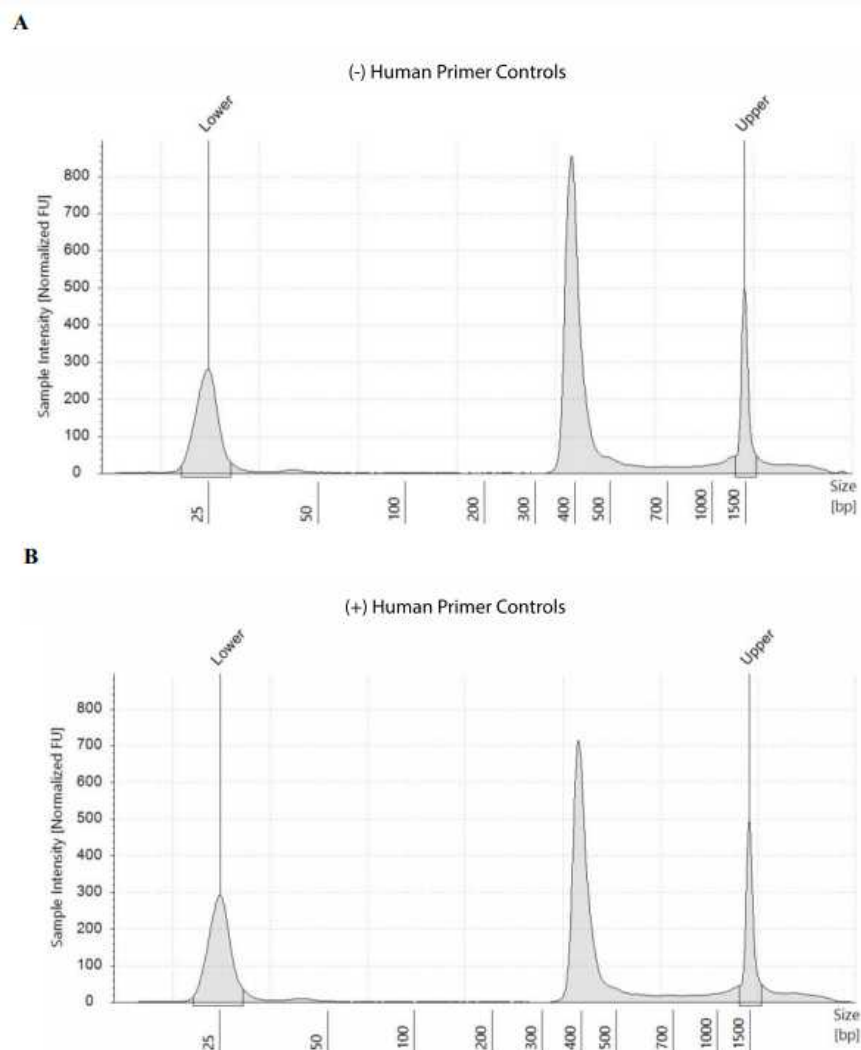


Figure 21: ARTIC SARS-CoV-2 cDNA amplicons generated from 1000 total viral copies in the absence (A) or presence (B) of the human primer controls.

**Samples can be stored at **-20 °C** for up to a week.**

## 22 Proceed to Library Prep in the Next Step.

### Fragmentation/End Prep

## 23

| A                        | B                 | C            |
|--------------------------|-------------------|--------------|
| FRAGMENTATION SIZE RANGE | INCUBATION @ 37°C | OPTIMIZATION |
| 80 bp–250 bp             | 30 min            | 30–40 min    |

**Fragmentation occurs during a 37°C incubation step. A fragmentation time of 30 minutes should be sufficient for most samples to generate mainly 100 bp fragments, however, the time may need to be optimized. For high input samples (1 µg), fragmentation time may need to be increased.**

## 24

Ensure that the Ultra II FS Reaction Buffer is completely thawed. If a precipitate is seen in the buffer, pipette up and down several times to break it up, and quickly vortex to mix. Place **On ice** until use.

There are short videos on NEB.com about handling the FS reagents:

Preparing FS Reagents: <https://www.neb.com/tools-and-resources/video-library/quick-tips--preparing-the-nebnext-ultra-ii-fs-dna-reaction-buffer-and-enzyme-mix>

Making FS master mix: <https://www.neb.com/tools-and-resources/video-library/quick-tips--preparing-nebnext-ultra-ii-fs-dna-reaction-buffer-and-enzyme-mix-master-mix>

## 25

Vortex the Ultra II FS Enzyme Mix 5-8 seconds prior to use and place **On ice**.

**It is important to vortex the enzyme mix prior to use for optimal performance**

## 26

Add the following components to a 0.2 ml thin wall PCR tube **On ice**:

| A  | B                      |
|--|------------------------|
| COMPONENT                                    | VOLUME PER ONE LIBRARY |
| ARTIC SARS-CoV-2 cDNA (Step 20)              | 13 µl                  |
| (yellow) NEBNext Ultra II FS Reaction Buffer | 3.5 µl                 |
| (yellow) NEBNext Ultra II FS Enzyme Mix      | 1 µl                   |
| <i>Total Volume</i>                          | 17.5 µl                |

27 

Vortex the reaction for  **00:00:05** and briefly spin in a microcentrifuge.

28 

In a thermocycler, run the following program:

| A    | B          |
|------|------------|
| TEMP | TIME       |
| 37°C | 30 minutes |
| 65°C | 30 minutes |
| 4°C  | ∞          |

\*Set heated lid to 75°C

If necessary, samples can be stored at  **-20 °C** ; however, a slight loss in yield (~20%) may be observed. We recommend continuing with adaptor ligation before stopping.

#### Adaptor Ligation

29 

Add the following components directly to the FS Reaction Mixture:

| A   | B        |
|---|----------|
| COMPONENT                                   | VOLUME   |
| FS Reaction Mixture (Previous Step)         | 17.5 µl  |
| (red) NEBNext Adaptor for Illumina**        | 1.25 µl  |
| (red) NEBNext Ultra II Ligation Master Mix* | 15 µl    |
| <i>Total Volume</i>                         | 33.75 µl |

\* Mix the Ultra II Ligation Master Mix by pipetting up and down several times prior to adding to the reaction.

\*\* The NEBNext adaptor is provided in NEBNext Oligo kits. NEB has several oligo options which are supplied separately from the library prep kit. Please see [www.neb.com/oligos](https://www.neb.com/oligos) for additional information

**Do not premix adaptor with the Ligation Master Mix**

30 


Set a 100 µl or 200 µl pipette to 25 µl and then pipette the entire volume up and down at least 10 times to mix thoroughly. Perform a quick spin to collect all liquid from the sides of the tube.

**Caution: The NEBNext Ultra II Ligation Master Mix is very viscous. Care should be taken to ensure adequate mixing of the ligation reaction, as incomplete mixing will result in reduced ligation efficiency. The presence of a small amount of bubbles will not interfere with performance.**

31 

15m

Incubate at  $\uparrow$  **20 °C** for  $\text{⌚}$  **00:15:00** in a thermocycler with the heated lid off.

32 

Add  $\text{📦}$  **1.5  $\mu$ l (red) USER® Enzyme** to the ligation mixture from the previous step.

**Steps 32 and 33 are only required for use with NEBNext Adaptors. USER enzyme can be found in the NEBNext Multiplex Oligos ([www.neb.com/oligos](http://www.neb.com/oligos)).**

33   

15m

Mix well and incubate at  $\uparrow$  **37 °C** for  $\text{⌚}$  **00:15:00** with the heated lid set to  $\geq$  47°C.

**Samples can be stored  $\text{⌚}$  Overnight at  $\uparrow$  -20 °C .**

#### Cleanup of Adaptor-ligated cDNA

34 The following section is for cleanup of the ligation reaction.


**The volumes of NEBNext Sample Purification Beads provided here are for use with the sample contained in the exact buffer at this step (35.25  $\mu$ l; Step 33). These bead volumes may not work properly for a cleanup at a different step in the workflow, or if this is a second cleanup at this step. For cleanups of samples contained in different buffer conditions, the volumes may need to be experimentally determined**

35 Vortex the NEBNext Sample Purification Beads to resuspend.

36 

Add  $\text{📦}$  **28  $\mu$ l (0.8X) resuspended beads** to the Adaptor Ligation reaction. Mix well by pipetting up and down at least 10 times. Be careful to expel all of the liquid out of the tip during the last mix. Vortexing for 3-5 seconds on high can also be used. If centrifuging samples after mixing, be sure to stop the centrifugation before the beads start to settle out.



37 

Incubate samples at **Room temperature** for at least **00:05:00**.

38 Place the tube/plate on an appropriate magnetic stand to separate the beads from the supernatant. If necessary, quickly spin the sample to collect the liquid from the sides of the tube or plate wells before placing on the magnetic stand.

39 After 5 minutes (or when the solution is clear), carefully remove and discard the supernatant. Be careful not to disturb the beads that contain DNA targets.

**Caution: do not discard the beads.**

40 

Add **200 µl 80% freshly prepared ethanol** to the tube/plate while in the magnetic stand. Incubate at **Room temperature** for **00:00:30**, and then carefully remove and discard the supernatant. Be careful not to disturb the beads that contain DNA targets.

41 

Repeat the previous step once for a total of two washes:

Add **200 µl 80% freshly prepared ethanol** to the tube/plate while in the magnetic stand. Incubate at **Room temperature** for **00:00:30**, and then carefully remove and discard the supernatant. Be careful not to disturb the beads that contain DNA targets.

**Be sure to remove all visible liquid after the second wash. If necessary, briefly spin the tube/plate, place back on the magnet and remove traces of ethanol with a p10 pipette tip.**

42 Air dry the beads for up to 5 minutes while the tube/plate is on the magnetic stand with the lid open.

**Caution: Do not over-dry the beads. This may result in lower recovery of DNA. Elute the samples when the beads are still dark brown and glossy looking, but when all visible liquid has evaporated. When the beads turn lighter brown and start to crack, they are too dry.**

43 

Remove the tube/plate from the magnetic stand. Elute the DNA target from the beads by adding **10 µl 0.1X TE**.

44 

2m

Mix well by pipetting up and down 10 times, or on a vortex mixer. Incubate for at least **00:02:00** at

**Room temperature**. If necessary, quickly spin the sample to collect the liquid from the sides of the tube or plate wells before placing back on the magnetic stand.

45  

Place the tube/plate on the magnetic stand. After 5 minutes (or when the solution is clear), transfer **7.5 µl** to a new PCR tube.

Samples can be stored at **-20 °C**.

46 Proceed to PCR Enrichment of Adaptor-ligated DNA in the next step.

## PCR Enrichment of Adaptor-ligated DNA

47

Check and verify that the concentration of your oligos is 10 µM.

**Follow step-case A. if you are using the following oligos (10 µM primer):**

**Use option A for any NEBNext Oligo kit where index primers are supplied in tubes. These kits have the forward and reverse primers supplied in separate tubes.**

**Follow step-case B. if you are using the following oligos (10 µM primer):**

**Use Option B for any NEBNext Oligo kit where index primers are supplied in a 96-well plate format. These kits have the forward and reverse (i7 and i5) primers combined.**

Step 47 includes a Step case.

**Primers: Tubes**

**96-well plate**

step case

### Primers: Tubes

**Index primers are supplied in tubes.** These kits have the forward and reverse primers supplied in separate tubes.

48 Add the following components to a sterile strip tube:

| A   | B       |
|---|---------|
| COMPONENT                                 | VOLUME  |
| Adaptor Ligated DNA Fragments (Step 45)   | 7.5 µl  |
| (blue) NEBNext Library PCR Master Mix     | 12.5 µl |
| (blue) Index Primer/i7 Primer*,**         | 2.5 µl  |
| (blue) Universal PCR Primer/i5 Primer*,** | 2.5 µl  |
| <i>Total Volume</i>                       | 25 µl   |

\* NEBNext Oligos must be purchased separately from the library prep kit. Refer to the corresponding NEBNext Oligo kit manual for determining valid barcode combinations.

\*\* Use only one i7/primer/index primer per sample. Use only one i5 primer (or the universal primer for single index kits) per sample.

49 

Set a 100 µl pipette to 20 µl and then pipette the entire volume up and down at least 10 times to mix thoroughly. Perform a quick spin to collect all liquid from the sides of the tube.

50 

Place the tube on a thermocycler and perform PCR amplification using the following PCR cycling conditions:

| A                    | B    | C          | D      |
|----------------------|------|------------|--------|
| CYCLE STEP           | TEMP | TIME       | CYCLES |
| Initial Denaturation | 98°C | 30 seconds | 1      |
| Denaturation         | 98°C | 10 seconds | 5*     |
| Annealing/Extension  | 65°C | 75 seconds |        |
| Final Extension      | 65°C | 5 minutes  | 1      |
| Hold                 | 4°C  | ∞          |        |

\* Set heated lid to 105°C. The number of PCR cycles recommended should be viewed as a starting point and may need to be optimized for particular sample types.

51 Proceed to Cleanup of PCR Amplification in the next step.

Cleanup of PCR Reaction 7m

52

**The NEBNext Sample Purification Beads provided here are for use with the sample contained in the exact buffer at this step. These volumes may not work properly for a cleanup at a different step in the workflow. For cleanups of samples contained in different buffer conditions, the volumes may need to be experimentally determined.**

53 Vortex NEBNext Sample Purification Beads to resuspend.

54 

Add **22.5 µl (0.9X) resuspended beads** to the PCR reaction. Mix well by pipetting up and down at least 10 times. Be careful to expel all of the liquid out of the tip during the last mix. Vortexing for 3-5 seconds on high can also be used. If centrifuging samples after mixing, be sure to stop the centrifugation before the beads start to settle out.

55 

5m




Incubate samples on bench top for at least **00:05:00** at **Room temperature**.

56 Place the tube/plate on an appropriate magnetic stand to separate the beads from the supernatant. If necessary, quickly spin the sample to collect the liquid from the sides of the tube or plate wells before placing on the magnetic stand.

- 57 After 5 minutes (or when the solution is clear), carefully remove and discard the supernatant. Be careful not to disturb the beads that contain DNA targets.




**Caution: do not discard the beads.**

- 58 

Add  **200 µl 80% freshly prepared ethanol** to the tube/plate while in the magnetic stand. Incubate at  **Room temperature** for  **00:00:30**, and then carefully remove and discard the supernatant. Be careful not to disturb the beads that contain DNA targets.

- 59 

Repeat the previous step once for a total of two washes:

Add  **200 µl 80% freshly prepared ethanol** to the tube/plate while in the magnetic stand. Incubate at  **Room temperature** for  **00:00:30**, and then carefully remove and discard the supernatant. Be careful not to disturb the beads that contain DNA targets.

Be sure to remove all visible liquid after the second wash. If necessary, briefly spin the tube/plate, place back on the magnet and remove traces of ethanol with a p10 pipette tip.

- 60 Air dry the beads for up to 5 minutes while the tube/plate is on the magnetic stand with the lid open.



**Caution: Do not over-dry the beads. This may result in lower recovery of DNA. Elute the samples when the beads are still dark brown and glossy looking, but when all visible liquid has evaporated. When the beads turn lighter brown and start to crack, they are too dry.**

- 61 



Remove the tube/plate from the magnetic stand. Elute the DNA target from the beads by adding  **17 µl 0.1X TE**.

- 62 

2m

Mix well by pipetting up and down 10 times, or on a vortex mixer. Incubate for at least  **00:02:00** at  **Room temperature**. If necessary, quickly spin the sample to collect the liquid from the sides of the tube or plate wells before placing back on the magnetic stand.


- 63 

Place the tube/plate on the magnetic stand. After 5 minutes (or when the solution is clear), transfer  **15 µl** to a new PCR tube and store at  **-20 °C**.



- 64 Assess the library size distribution with Agilent Bioanalyzer or TapeStation high sensitivity DNA reagents. The sample may need to be diluted before loading. A peak sized of 200–250 bp is expected, based on a 30-minute fragmentation time. (Figure 64)

A.

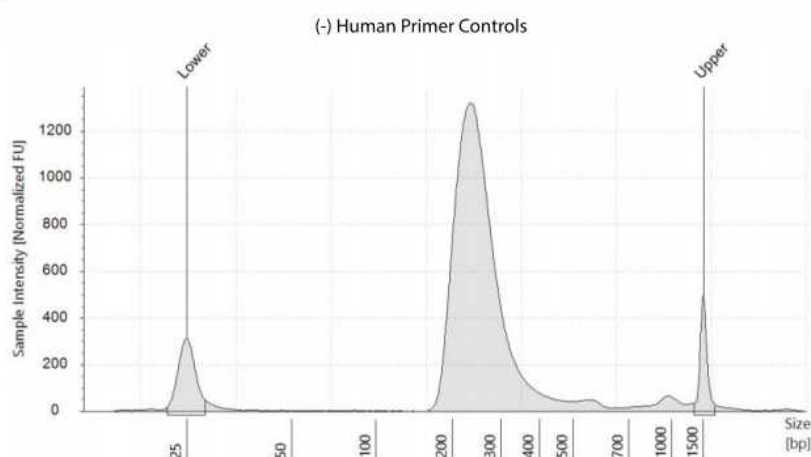

B.

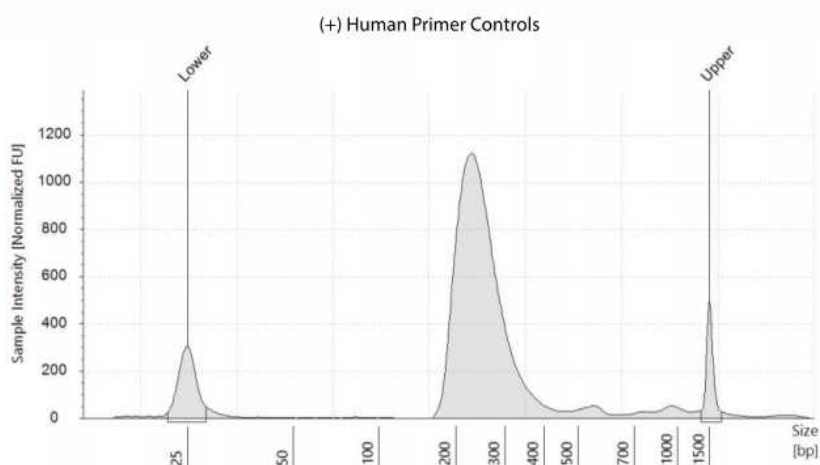


Figure 64. : Example of final library size distributions on a TapeStation. ARTIC SARS-CoV-2 libraries were generated from 1000 viral copies in the absence (A) or presence (B) of the human primer controls.