

VERSION 1

NOV 09, 2022



IN DEVELOPMENT

Monoclonal Culture (Limiting Dilution) V.1

This protocol is published without a DOI.

Binghua Darren Zhang¹

¹University of Glasgow



Binghua Darren Zhang <u>University of Glasgow</u>

COMMENTS 0

ABSTRACT

This protocol gives a step-by-step guidance on how to generate the homogenised gene-knockout (KO) monoclonal cell line from a polyclonal "KO" cells pool.

PROTOCOL CITATION

Binghua Darren Zhang 2022. Monoclonal Culture (Limiting Dilution) . **protocols.io** https://protocols.io/view/monoclonal-culture-limiting-dilution-ci2nugde

MANUSCRIPT CITATION please remember to cite the following publication along with this protocol

This protocol is generated according to the protocol established and published from Addgene (https://www.addgene.org/protocols/limiting-dilution/)

KEYWORDS

Limiting dilution, Monoclonal cell culture, Single cell clone

LICENSE

This is an open access protocol distributed under the terms of the <u>Creative Commons</u>

<u>Attribution License</u>, which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited



1

CREATED

Nov 08, 2022

LAST MODIFIED

Nov 09, 2022

PROTOCOL INTEGER ID

72494

GUIDELINES

Prepare conditioned growth medium;

Cell counting and calculation;

Culture single cell suspension in 96-well cell culture plates;

MATERIALS TEXT

Stable polyclonal cell pool;

Fresh growth medium;

DPBS Buffer;

Trypsin;

Trypan Blue;

50 mL Falcon tubes;

0.2 µm syringe filter with syringe;

0.40 µm cell strainer mesh;

Hemocytometer;

Reagent Reservoir;

Multi-channel pipette (≥ 100 μL);

96-well cell culture plate(s);

BEFORE STARTING

Generate and prepare the stable polyclonal cell pool;

Prep Conditioned Medium

 ${\bf 1} \qquad \hbox{Collect the used growth medium into a 50 mL Falcon tube}.$

Note

Do not collect the used growth medium if there is contamination or the medium turns orangish/yellowish.

5m

2 Pipette the used growth medium through a 0.2 µm syringe filter to remove the waste product and dead cells.

m protocols.io

3 Collect the supernatant (about 15 mL) used growth medium into a new 50 mL Falcon tube. And then discard the waste product. $\underline{\underline{L}}$ 15 mL For three 96-well plates

Note

Collect 15 mL used growth medium if there are three 96-well cell culture plates ready for limiting dilution. Otherwise, collect 10 mL for two plates or 5 mL for one plate only.

4 Add fresh growth medium up to 30 mL (if there are three 96-well cell culture plates).

△ 30 mL For three 96-well plates

Note

Roughly mix about 1:1 fresh and used growth medium. Add fresh medium up to 20 mL if there are two 96well cell culture plates or up to 10 mL if there is only one plate.

5 Mix the medium thoroughly and seal the tube tightly for later use as the conditioned medium.

Prep the "KO" Cells

6 Rinse cells with sufficient amount of DPBS Buffer.


7 Aspirate the DPBS with aspirator.

8 Trypsinise the cells with 1 mL of Trypsin in the culture incubator at 37°C for 2 min.

Д 1 mL



7m

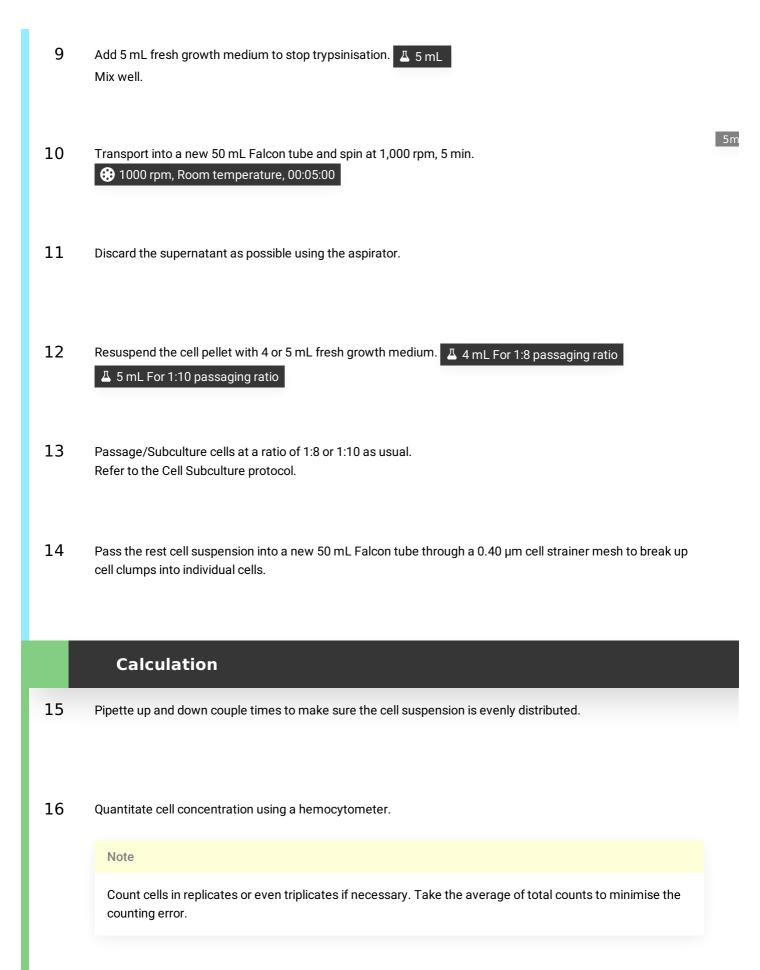

- 17 Sample Calculation (Example Calculation)
- 17.1 Expect to have a cell concentration of $\sim 3x10^6$ cells/mL.
- 17.2 $_{10X \text{ Dilution to make the concentration of } \sim 3x10^5 \text{ cells/mL}}$. Count cells in replicates.
- 17.3 Expect to have the final concentration at ~5 cells/mL. For 30mL conditioned medium, a total of ~150 cells are expected. Thus, (~150 cells)/(~3x10⁵ cells/mL) $\approx 0.5 \,\mu$ L of cell suspension is required.
- Further dilute cell suspension in 1:100 DF to obtain an expected concentration of $\sim 3x10^3$ cells/mL. (~ 150 cells)/($\sim 3x10^3$ cells/mL) ≈ 50 µL of cell suspension is required.

Monoclonal Culture Plate

- Pipette $50 \, \mu L$ which carries ~150 cells from the cell suspension into the prepared $30 \, mL$ conditioned medium (Section 1).

 Mix thoroughly.
- Transfer the medium into a tray.
- Use multi-channel pipette to transfer 100 μ L into each well, leading to the cell concentration of ~0.5 cell/well. \pm 100 μ L
- Place one control well containing ~1,000 cells in each plate. In this case, pipette ~3 μ L from the cell suspension of ~3x10⁵ cells/mL (Section 17.2) into the control well and labelled.

Place the plates into the incubator until they are ready for single cell colonies screening.