



OCT 20, 2023

Illumina Amplicon Sequencing using two step PCR

 Forked from [Illumina Amplicon Sequencing using two step PCR](#)

Lindsay Kate Newbold¹, Jonathan Joetay¹, Warren¹,

Amy

Daniel S Read¹, kerry.walsh¹, Thorpe¹

¹UK Centre for Ecology & Hydrology

UKCEH Molecular Ecology

OPEN  ACCESS



Lindsay Kate Newbold

UK Centre for Ecology & Hydrology

ABSTRACT

This document is designed to provide the researcher with all of the information required to undertake two step amplicon sequencing using the Illumina MiSeq platform at UKCEH-Wallingford.

The aim is to provide a framework upon which the researcher may choose to add their own modifications to suit novel applications.

Protocol Citation: Lindsay Kate Newbold, Joetay, Jonathan Warren, Daniel S Read, kerry.walsh, Amy Thorpe 2023. Illumina Amplicon Sequencing using two step PCR. **protocols.io** <https://protocols.io/view/illumina-amplicon-sequencing-using-two-step-pcr-cn4vvgw6>

License: This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](#), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

Protocol status: Working
We use this protocol and it's working

Created: Feb 08, 2023

Last Modified: Oct 20, 2023

PROTOCOL integer ID:
76661

Nextera index plates- Prepared prior to Step 1 PCR

- Order indexing primers direct from oligo manufacturer (For example IDT, Sigma genosys or MWG) suspended in water $0.5 \mu\text{M}$ scale, Desalt purification and $10 \mu\text{M}$ concentration. These indexing primers consist of: an Illumina Nextera adapter i5 (Forward primer) AATGATACGGCGACCACCGAGATCTACAC or i7 (reverse primer) CAAGCAGAAGACGGCATACGAGAT, Unique 8bp barcode sequence and pre-adapter i5 TCGTCGGCAGCGTC (F) or i7 GTCTCGTGGGCTCGG (R). Barcode design when demultiplexed allows for the unique assignment of an individual sequence to sample. Figure 1 gives a graphical representation of how 16F, and 24R barcodes sequences are allocated to give 384, unique barcode combinations. Full barcode sequence is given in attached document.

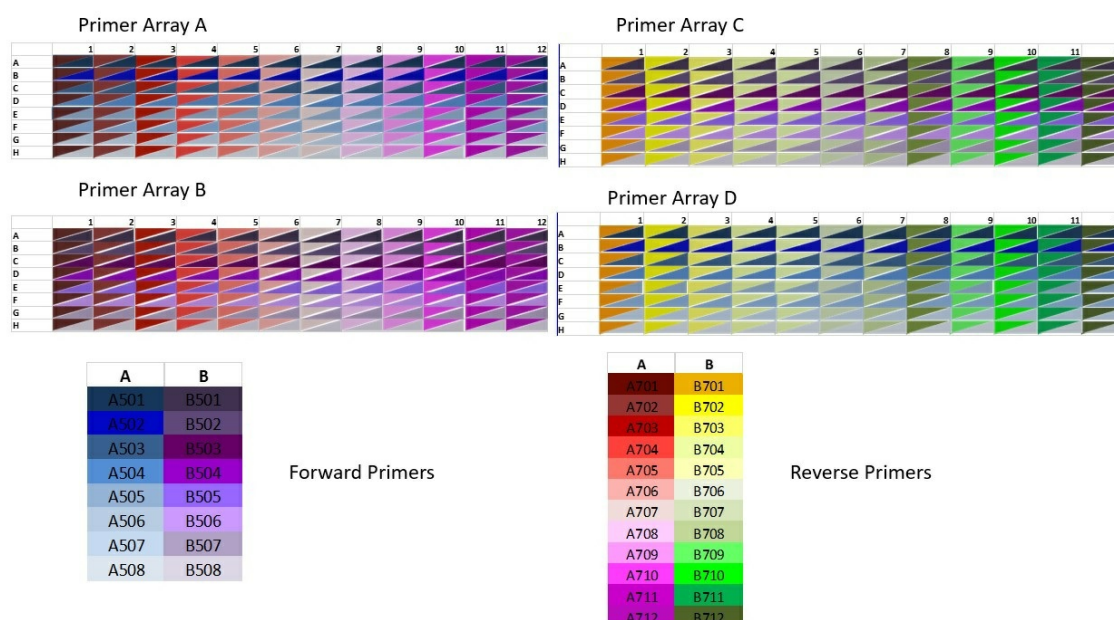


Figure 1: An illustration of Dual-index primer pair allocation

[Array Plate description.xlsx](#)

- Prepare index plate 'masters' manually allocating $250 \mu\text{L}$ of each diluted primer stock $10 \text{ micromolar } (\mu\text{M})$ per well as per plate design detailed above into Deep-Well microtiter plates. Alternatively this can be done using a liquid handling Robot.

- 3 Make PCR plate clones of the 'master' plates, each well containing 5ul of mixed primer, label, seal and store at -20 for future use.

Step I- Amplicon PCR with modified primers

- 4 Make a standardised dilution plate of template DNA from DNA extraction process (See Extraction protocol). DNA concentration is determined through either Nanodrop or Qubit BR assay of raw extract, and then each extract is diluted to a standardised concentration (this is subject to extraction success, but shouldn't be below 5ng/ul to allow for efficient PCR), using the equation:

Desired Concentration (ng/ul)

X Final desired volume = Amount of DNA to add

Raw extract Concentration (ng/ul)



To work out the amount of water its : Final desired volume - Amount of DNA to add

- 5 Step 1 Amplicon primers consist of a pre adapter, sequencing primer and specific locus primer. Target primers are designed as per Table 1:

A	B	C	D	E	F	G
Source reference	Universal primer name	Target Group	pre-adapter	Sequencing primer sequence	Specific locus primer	Combined sequence
Walters et al (2015) https://journals.asm.org/doi/10.1128/mSystems.00009-15	515F	Bacteria	TCGTCGGCAGCGTC	AGATGTGTATAAGAGACAG	GTGYCAGCMGCCGCGGTAA	TCGTCGGCAGCGTCAGATGTGTATAAGAGACAGGTGYCAGCMGCCGCGGTAA
Walters et al (2015) https://journals.asm.org/doi/10.1128/mSystems.00009-15	806R	Bacteria	GTCTCGTG GGCTCGG	AGATGTGTATAAGAGACAG	GGACTA CNVGGG TWTCTAAT	GTCTCGTG GGGCTCGGAGATGTGTATAAGAGACAGGGACTA CNVGGGTWTC TAAT
Ihrmark et al (2012) https://doi.org/10.1111/j.1574-6941.2012.01437.x	ITS7F	Fungi	TCGTCGGCAGCGTC	AGATGTGTATAAGAGACAG	GTGARTCATCGAATCTTTG	TCGTCGGCAGCGTCAGATGTGTATAAGAGACAGGTGARTCATCGAATCTTTG
Ihrmark (2012) https://doi.org/10.1111/j.1574-6941.2012.01437.x	ITS4R	Fungi	GTCTCGTG GGCTCGG	AGATGTGTATAAGAGACAG	TCCTCCGCTTAT TGATATGC	GTCTCGTG GGGCTCGGAGATGTGTATAAGAGACAGTCCTCCGCTTATTGATATGTC

A	B	C	D	E	F	G
Mangot et al (2012) https://doi.org/10.1111/1462-2920.12065	NSF563	Eukaryotes	TCGTCGGC AGCGTC	AGATGTGT ATAAGAGA CAG	CGCGGT AATTCC AGCTCC A	TCGTCGGCAG CGTCAGATGT GTATAAGAGA CAGCGCGGTA ATTCCAGCTC CA
Mangot et al (2012) https://doi.org/10.1111/1462-2920.12065	NSR951	Eukaryotes	GTCTCGTG GGCTCGG	AGATGTGT ATAAGAGA CAG	TTGGYR AATGCT TTCGC	GTCTCGTGGG CTCGGAGATG TGTATAAGAG ACAGTTGGYR AATGCTTTCG C
Kelley et al (2018) ISBN: 978-1-84911-406-6	rbcL-646F	Diatoms	TCGTCGGC AGCGTC	AGATGTGT ATAAGAGA CAG	ATGCGT TGGAGA GARGT TTC	TCGTCGGCAG CGTCAGATGT GTATAAGAGA CAGATGCGTT GGAGAGARCG TTTC
Kelley et al (2018) ISBN: 978-1-84911-406-6	rbcL-998R	Diatoms	GTCTCGTG GGCTCGG	AGATGTGT ATAAGAGA CAG	GATCAC CTTCTA ATTTAC CWACAA CTG	GTCTCGTGGG CTCGGAGATG TGTATAAGAG ACAGGATCAC CTTCTAATTTA CCWACAACG

Table 1: Step I amplicon primers employed during for the amplication of biofilm communities

- 6 Using column 'Combined sequence' given in Table 1, order primers direct from oligo manufacturer (For example IDT, Sigma genosys or MWG) suspended in water  0.5 µM scale, Desalt purification and  100 µM concentration.




7

A	B	C
Reagent	Per Sample	Per Plate
Molecular Grade Water	40.5 ul	4050 ul
5 X Buffer	1 ul	100 ul
High GC Buffer	1 ul	100 ul
10 mM dNTP's	1 ul	100 ul
Q5 Taq Polymerase	0.5 ul	50 ul
Primer F	0.5 ul	50 ul
Primer R	0.5 ul	50 ul
DNA Template (at ~10 ng/ml)	2 ul	Added separately

Table 2: 1 Step PCR Reagents in the order in which they should be added. *Note that Template is added separately after PCR Mastermix plate has been prepared.*

- 8 All reagents from above Table 2 (except template DNA) are added to a 5ml Tube (or equivalent) to form a PCR Mastermix.

Mastermix is prepared on ice, using filter tips and under sterile conditions in the PCR hood.

Once Mastermix has been made add  48 µL of Mastermix to each well of a 96 well PCR plate using a multichannel pipette (accurate for volumes between  100 µL -  300 µL).


We recommend the use of :

Axygen™ 96-Well Low Profile PCR Microplates (SKU AXY-PCR-96-LP-FLT-C)

<https://www.thistlescientific.co.uk/product/96-well-polypropylene-flat-top-pcr-microplate-low-profile-no-skirt-clear-nonsterile/>

and

Equipment	
Research Plus Multichannel pipette	NAME
Multichannel Pipette	TYPE
Eppendorf	BRAND
3125000052	SKU
https://online-shop.eppendorf.co.uk/UK-en/Manual-Liquid-Handling-44563/Manual-Pipetting--Dispensing-44564/Eppendorf-Research-plus-PF-534798.html	LINK

Add  2 µL of standard dilution from step 4 to respective position in the mastermix plate. This ensure's full traceability when samples are demultiplexed during the bioinformatics steps performed post sequencing. It is therefore imperative that experimentors make a note of which sample is added to each well of the PCR Mastermix plate.

- 9 Seal plates with PCR film, and centrifuge, to concentrate sample.

We recommend the use of:

<https://www.thermofisher.com/order/catalog/product/AB0558>

10 Initiate PCR reaction on Thermal cycler

We recommend the BIORAD C1000 touch thermal cycle <https://www.bio-rad.com/en-uk/product/c1000-touch-thermal-cycler?ID=LGTW9415>

PCR programme for 16S, 18S and ITS:

- Denature 95 °C for 02:00m
- Denature 95 °C 00:15s
- Anneal 55 °C (16S) 52 (ITS) for 00:30s
- Extension 72 °C for 00:30s
- Repeat ii. – iv. 30 cycles
- Final extension 72 °C for 10:00m

PCR programme for rbcl:

- Denature 95 °C for 02:00m
- Denature 95 °C 00:20s
- Anneal 55 °C for 00:45s
- Extension 72 °C for 00:60s
- Repeat ii. – iv. 35 cycles
- Final extension 72 °C for 5:00m


Gel Electrophoresis

11 Verify PCR products using gel electrophoresis.

11.1 Prepare the appropriate electrophoresis tank and casting tray. The number and type of gels required will depend on the number of samples.




A Mini-sub Cell GT Horizontal Electrophoresis Tank. Can run up to 28 samples + two ladder lanes. The 7 x 10 cm casting tray requires 50 ml agarose mix.

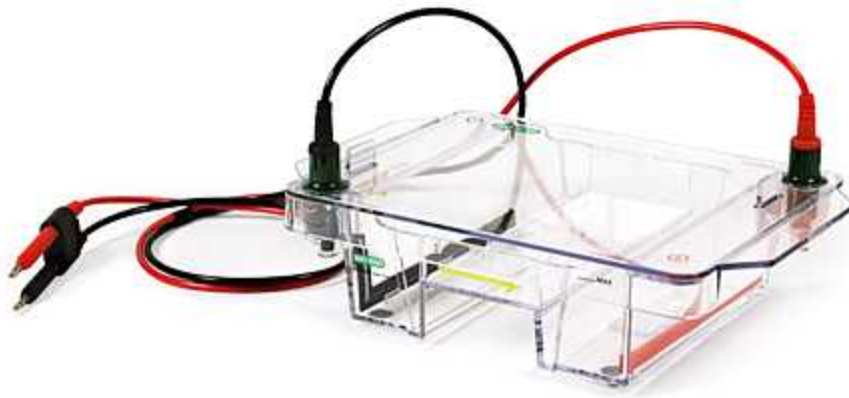



A Wide Mini-sub Cell GT Horizontal Electrophoresis Tank (Midi Tank). Can run up to 58 samples + two ladder lanes. The 15 x 10cm casting tray requires 100ml agarose mix.




Sub-Cell 96 GT Horizontal Electrophoresis Tank. Can run up to 100 samples + two ladder lanes. The 25 x 15 cm casting tray requires 250ml agarose mix.

- 11.2** To make a [M] 1 % (v/v) agarose gel : For a Midi Gel, measure 1g of molecular grade agarose powder into a glass pyrex beaker and make up to  100 mL using 1 X TBE buffer.

4m

Heat in microwave to dissolve for approx  00:02:00 . When no agarose crystals are visible, cool the gel in a sink filled with a little water, to just above room temperature,. Cooling in this way will take atleast  00:02:00 . As a rule of thumb if you're able to hold the base of the beaker without discomfort the gel is ready.


- 11.3** To a 100ml mix add 1 μL ethidium bromide, or 5 μL Gel Red
Swirl to mix and pour into gel casting tray set up as below.

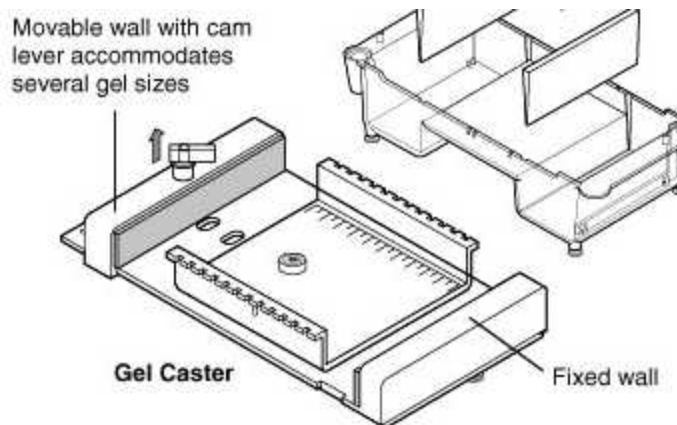



Biorad gel caster. To use push two casting walls together and tighten using lever.


- 11.4** Place combs directly into the casting tray (in above tray line up comb ends with 2nd notch and push down). If bubbles are present these can be pushed out of gel using combs before final positioning and tightening lever.

45m

Leave gel to set for atleast  00:45:00 .

- 11.5** Whilst gel is setting, prepare samples using  5 μL PCR amplicon to  1 μL Bromophenol blue gel loading buffer.

<https://www.thermofisher.com/order/catalog/product/10816015#/10816015>

- 11.6** When set load loosen casting lever and place tray and gel into Gel tank filled with 1X TBE buffer to fill line. Remove combs (doing this in the buffer prevents bubbles forming in wells). Add  5 μL Hyper ladder 1KB <https://www.bioline.com/hyperladder-1kb.html> (or similar) to first well of agarose gel by placing tip into well and gently depressing pipette control, whilst slowly removing tip from well. Repeat procedure with samples, each being placed in its an individual well.

- 11.7** When all samples are loaded put on gel tank lid, and plug in cables (red to red, black to black).

45m

Set gel to run for  00:45:00 at 90 Volts.

- 11.8 When run complete image gel using gel image system (eg <https://www.bio-rad.com/en-uk/product/gel-doc-xr-gel-documentation-system?ID=O494WJE8Z>).

A successful PCR should look similar to below.

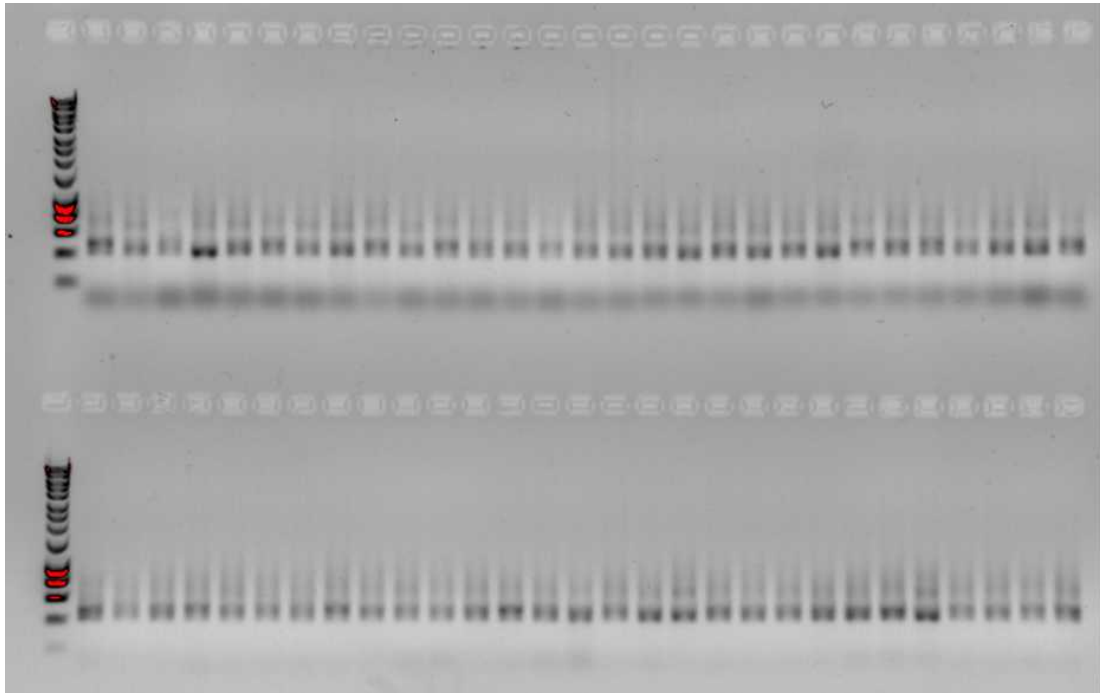


Figure 4: Result of ITS2 PCR.

PCR Clean-up

12



Clean-up Amplicon PCR step I using ZymoZR-96 kit (standard Manufacturers protocol), again using filter tips and utmost care throughout as samples are not indexed.

<https://zymoresearch.eu/collections/dcc-pcr-purification-kits/products/zr-96-dna-clean-up-kit>

Alternatively use Multiscreen-PCR₉₆ Filter plate under standard manufacturers protocol if liquid handling robot and vacuum manifold are available.

https://www.merckmillipore.com/GB/en/product/MultiScreen-PCR96-Filter-Plate,MM_NF-MSNU03050

Step II - Nextera indexing PCR

- 13 Use  1-10 µL Cleaned PCR product as template for Step II- Indexing PCR (Table 3). Water added to PCR master mix is adjusted to make each reaction total  25 µL .


A	B	C
A	B	C
Reagent	Per Sample	Per Plate
Molecular Grade Water	7.25 ul	725 ul
5 X Buffer	5 ul	500 ul
High GC Buffer	5 ul	500 ul
10 mM dNTP's	0.5 ul	50 ul
Q5 Taq Polymerase	0.25 ul	25 ul
Primer Array mix	5 ul	Added separately
PCR Template	2 ul	Added separately


Table 3: Step II PCR Reagents in the order in which they should be added. *Note that Template and primer mix are added separately after PCR Mastermix plate has been prepared.*

- 14 PCR programme :
- Denature 95 °C for 02:00m
 - Denature 95 °C for 00:15s
 - Anneal 55 °C for 00:30s
 - Extension 72 °C for 00:30s
 - Repeat to total 8 cycles



Normalisation

45m

- 15 Resultant PCR products are normalised using the Norgen NGS Normalization 96-Well Kit – 61900 under standard manufacturers protocol:
UK distributors
<https://geneflow.co.uk/product/ngs-normalization-96-well-kit-61900/>
- 16 Use ~  5 µL normalised PCR product from each sample to form a pooled library.

17 Vacuum concentrate to a volume of approx  30 µL pooled normalised plate library.
We recommend:
<https://www.eppendorf.com/gb-en/eShop-Products/Centrifugation/Concentrator/Concentrator-plus-p-PF-25748>


17.1 Whilst Vacuum concentration is occurring prepare a 50ml 1% Agarose gel, with a wide comb.

18 Add  5 µL loading buffer to concentrate and run for  00:45:00 at 90 Volts.

45m

18.1 Excise resultant band using sterile scalpel, and visualised on blue light Led trans illuminator

<https://www.smobio.com/shop/product/ve0100-b-boxtm-blue-light-led-epi-illuminator-419>

19 Extract DNA from agarose Gel using [QIAquick Gel Extraction Kit - QIAGEN](#) under manufacturers standard protocol. Samples eluted in  30 µL buffer EB.

Quantification

20 Gel Purified libraries are quantified using **Invitrogen™ Qubit™ 1X dsDNA High Sensitivity (HS) quantification kit**

<https://www.fishersci.co.uk/shop/products/qubit-1x-dsdna-hs-br-assay-kits/15850210>

Under manufacturers protocol and companion qubit fluorometer

<https://www.fishersci.co.uk/shop/products/qubit-4-fluorometer-wifi/16223001#?keyword=Qubit%20%20flurometer>

21 Calculate Molar concentration using the https://www.bioline.com/media/calculator/01_07.html

ITS Amplicon Size 488 bp

18S Amplicon Size 533 bp

16S Amplicon Size 450 bp

rbcL Amplicon Size 533 bp

If using multiple libraries mix equal concentrations of each ITS arrays: eg. ITSA+B.

Your library is now ready for sequencing!

- 22 Please note that much of the following is an ammended version of the SChloss lab SOP, therefore their work should always be referenced.

Remove a V3 600 cycle reagent cartridge from the -20°C freezer. Place in room temperature water bath for one hour. Place Hyb buffer tube in 4°C fridge. While reagent cartridge is thawing, perform steps 16-20

- 23 Prepare fresh 0.1 Molarity (M) NaOH.

- 24 To a 1.5 mL tube, add 10 μL of library and 10 μL of 0.1 Molarity (M) NaOH. To a separate tube add 2 μL PhiX, 3 μL PCR grade water, and 5 μL of 0.1 Molarity (M) NaOH. Pipette to mix. Note: NaOH concentration on the flow cell must remain under 0.001M.

- 25 Allow the tubes to incubate at room temperature for 00:05:00. Immediately add 980 μL of ice-cold Hyb to the library tube, and 990 μL Hyb to the PhiX tube.

5m

Note: the resulting 20pM PhiX can be frozen and used for subequent runs.

- 26 Use Hyb to further dilute both the library and PhiX to 6pM for a v3 kit. See example below:

- $(1.45\text{ nM library} \times 10\text{ ul}) + (0.1\text{N NaOH} \times 10\text{ ul}) + 980\text{ ul Hyb} = 14.5\text{pM Lib, } 0.002\text{N NaOH}$
- $(14.5\text{pM lib} \times 275.86\text{ ul}) + 724.14\text{ ul Hyb} = 4.0\text{pM lib, } 0.00055\text{N NaOH}$
- $[(10\text{nM PhiX} \times 2\text{ ul}) + 3\text{ ul H}_2\text{O}] + (0.2\text{N NaOH} \times 5\text{ ul}) + 990\text{ ul Hyb} = 20\text{pM PhiX, } 0.001\text{N NaOH}$

NaOH

- $(20\text{pM PhiX} \times 200\text{ ul}) + 800\text{ ul Hyb} = 4.0\text{pM PhiX, } 0.0002\text{N NaOH}$
- $(4.0\text{pM Lib} \times 900\text{ ul}) + (4.0\text{pM PhiX} \times 100\text{ ul}) = \text{solution loaded}$
- Solution loaded is 4.0pM overall with a 3.6pM Library concentration, 0.4pM PhiX concentration, and 0.000515N NaOH

- 27 For a 10% PhiX run, combine 900 μL of 6 picomolar (pM) Library and 100 μL PhiX in a final tube.

2m

Vortex. At this final stage heat to 96°C for 00:02:00. Then cool immediately on ice.

- 28 When the reagent cartridge has thawed, dry bottom with paper towel. Invert the cartridge repeatedly to check each well is thawed. This also serves to mix the reagents. Place in Fridge

- 29** Load 600 ul of the final Library/PhiX solution into well 17 on the reagent cartridge.
- 30** Unbox flow cell and PR2 bottle. Thoroughly rinse the flow cell with Milli-Q water. Carefully dry by blotting with lint free wipes (Kimwipes). Give special attention to the edges and points of intersection between the glass and plastic. Wet a new wipe with 100% alcohol and wipe the glass on both sides avoiding the rubber intake ports. Visually inspect the flow cell to ensure there are no blemishes, particles, or fibers on the glass. Follow on screen instructions to load the flow cell, reagent cartridge, and PR2 bottle. Empty and replace the waste bottle.
- 31** Transfer reagent cartridge, flow cell, PR2 bottle, and flash drive with the sample sheet to the MiSeq. Open Sample Sheet from the flash drive in the MiSeq Local run manager software and queue run. Ensure the machine recognizes the correct sample sheet and the run parameters are correct.
- 32** Wait for the MiSeq to perform its pre-run checks, and press start. NOTE: If the pre-run check fails, try wiping down the flow cell again.
- 32.1** Ideal parameters for a 90% run:
- a. Cluster density 1000-1100k/mm² for v3 kits
 - b. >85% clusters passing filter
 - c. 10% aligned (amount of PhiX)
 - d. No spikes in corrected intensity plot
 - e. All indices identified following index reads
 - f. Final >Q30 score of >70%
- 32.2** When run complete, perform post run wash as per manufacturers recommendations.