

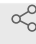


Jun 23, 2022

# QGP-1 cell line maintenance protocol

bellampalli.shreya<sup>1</sup><sup>1</sup>mayo clinic

1 Works for me

 Share[dx.doi.org/10.17504/protocols.io.bp2l61e1rvqe/v1](https://dx.doi.org/10.17504/protocols.io.bp2l61e1rvqe/v1) bellampalli.shreya

## ABSTRACT

**QGP-1 cell line maintenance protocol; thawing and passaging QGP-1 cells**

DOI

[dx.doi.org/10.17504/protocols.io.bp2l61e1rvqe/v1](https://dx.doi.org/10.17504/protocols.io.bp2l61e1rvqe/v1)

## PROTOCOL CITATION

bellampalli.shreya 2022. QGP-1 cell line maintenance protocol . **protocols.io**  
<https://dx.doi.org/10.17504/protocols.io.bp2l61e1rvqe/v1>



## LICENSE

————— This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](https://creativecommons.org/licenses/by/4.0/), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

## CREATED

Jun 23, 2022

## LAST MODIFIED

Jun 23, 2022

## PROTOCOL INTEGER ID

65196

### Thawing protocol

- 1 Warm media
- 2 Swirl cells in 37 degree water bath until vial is thawed

- 3 Add 7 mL of culture media to a T25 flask and gently pipette thawed cells into flask
- 4 Let cells adhere for 24 hours and replace half of the culture media with fresh media
- 5 Replace the media in the flask with fresh media after 48 hours
- 6 Split cells once cells reach 70-80% confluency: usually takes 3-4 days from thawing and plating in a T25 until they are ready for passage. (they won't ever become completely confluent, they will just start to grow on top of each other)
- 7 Lift cells using trypsin and plate all cells in a T75
- 8 They are usually ready to passage 3 days later.

#### passaging protocol

9

9.1 Warm media and trypsin in water bath

9.2 Clean hood with EtOH

9.3 Remove media

- 9.4 Wash with HBSS
- 9.5 Discard HBSS, add Trypsin, and place flasks in incubator for 5 minutes
- 9.6 Bring flasks back to hood, obtain lifted cells and spin at 500 rcf for 5 minutes
- 9.7 Discard supernatant
- 9.8 Resuspend in 10 mL
- 9.9 Place 250 uL of cells into a T75 flask and add 13 ml of media
- 9.10 Place flask back in incubator and split next week