




OCT 30, 2023

## Ex vivo electrophysiology

In 1 collection

Beatriz E Nielsen<sup>1</sup>

<sup>1</sup>University of Colorado Anschutz Medical Campus

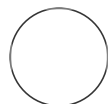

Beatriz E Nielsen

University of Colorado Anschutz Medical Campus

### ABSTRACT

This protocol describes the steps to perform whole-cell electrophysiology recordings in acute brain slices.

OPEN ACCESS



#### DOI:

[dx.doi.org/10.17504/protocols.io.j8nlkooowv5r/v1](https://dx.doi.org/10.17504/protocols.io.j8nlkooowv5r/v1)

**Protocol Citation:** Beatriz E Nielsen 2023. Ex vivo electrophysiology.

**protocols.io**

<https://dx.doi.org/10.17504/protocols.io.j8nlkooowv5r/v1>


**License:** This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](#), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited


**Protocol status:** Working

**Created:** Aug 16, 2023

**Last Modified:** Oct 30, 2023

## Set up

- 1 Rig set-up.
  - 1.1 Turn on required devices (computer, amplifier, manipulators, light sources, video camera controls, etc.)
  - 1.2 Open the N<sub>2</sub> tank connected to the rig anti-vibration table.
  - 1.3 Set 1X ACSF (add drugs needed for particular experiments) in a jug or bottle and bubble it with O<sub>2</sub>/CO<sub>2</sub>.
  - 1.4 Place the intake line into the ACSF container and allow circulation. Wait until the fluid has entered the recording chamber, then turn on the in-line heater (Warner Instruments) and set it to desired temperature  32-34 °C .
- 2 Recording electrodes or patch pipettes:
  - 2.1 Pull patch pipettes (1.5 – 2 MΩ) (World Precision Instruments) using an electrode puller (Narishige, PC-10).

2.2 Thaw an aliquot of the appropriate internal solution (stored at  -80 °C (see 'Solutions' section below) (add ATP, GTP and Phosphocreatine if solution does not already have it). Fill a syringe with filter with the internal solution and keep on ice by the rig.

3 Electrical stimulation electrodes:

3.1 Pull electrodes (World Precision Instruments) using a puller (Narishige, PC-10).

3.2 Fill electrodes with 1X ACSF using a syringe with filter.

## Whole-Cell recordings

5m

4 Transfer brain slice from incubation vial to the recording chamber and secure down the slice using a harp.

(The protocol for obtaining acute brain slices is linked below).

### Protocol



NAME

**Acute Brain Slices**


CREATED BY

kelsey.barcomb

**PREVIEW**

5 Locate and focus the desired region of the brain under the 4x objective.

- 6 Change the microscope lens to the 40x or 60x objective focus on healthy neurons selected for patching.
- 7 Fill a pipette tip with internal solution. Remove any air bubbles by gently flicking the glass pipette.
- 8 Place the pipette onto the wire electrode on the holder and tighten.
- 9 Apply a positive pressure and maintain it by quickly closing the stopcock.
- 10 Position the electrode using a micromanipulator.
- 11 Under the 40x-60x immersion objective, bring the tip of the pipette above the slice.
- 12 Approach the cell until the positive pressure create a small dimple.
- 13 Zero the pipette offset, release the positive pressure, and apply a small amount of negative pressure. The resistance should start to increase rapidly until a giga-Ohm seal is formed.
- 14 Clamp the cell at your resting potential of interest (typically -60 mV).

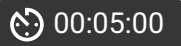- 15 Apply a few quick pulses of negative pressure to break into the cell and reach the whole-cell configuration.
- 16 Wait  00:05:00 to allow the internal solution to dialyze the cell before start your recordings. 5m  
Cells were discarded if series resistance was  $\geq 15 \text{ M}\Omega$ .  
Recordings can be made in voltage-clamp or current-clamp modes and the acquisition protocol will differ depending on the experiment.
- 16.1 Electrical stimulation: position the monopolar glass stimulating electrode filled with ACSF consistently 200  $\mu\text{m}$  away from the recorded cell. Select the appropriate number of pulses, intensity and duration of electrical stimulation for the experiment.
- 16.2 Optogenetic stimulation: deliver light pulses using using a 478 nm LED. Select the appropriate number of pulses, intensity and duration of light stimulation for the experiment.
- 17 Once finished the recording of a particular cell, set the holding potential to 0 mV and remove the electrode. Discard the electrode in sharps container.

## Solutions

- 18 Internal solutions:

■ K-Gluconate + 10 mM BAPTA

A	B	C
Drug	[mM]	g/100 mL
D-Gluconic Acid (K)	135	3.16
HEPES (K)	10	0.28
CaCl <sub>2</sub>	0.1	10 $\mu\text{L}$ (1M stock)

A	B	C
MgCl <sub>2</sub>	2	200 $\mu$ L (1M stock)
BAPTA-tetra potassium	10	0.628

+ [mM] 1 mg/mL ATP, [mM] 0.1 mg/mL GTP and [mM] 1.5 mg/mL phosphocreatine  
pH=7.35, 275 mOsm





■ K-Gluconate + 0.1 mM EGTA

A	B	C
Drug	[mM]	g/100 mL
D-Gluconic Acid (K)	135	3.16
HEPES (K)	10	0.28
CaCl <sub>2</sub>	0.1	10 $\mu$ L (1M stock)
MgCl <sub>2</sub>	2	200 $\mu$ L (1M stock)
EGTA	0.1	0.0038

+ [mM] 1 mg/mL ATP, [mM] 0.1 mg/mL GTP and [mM] 1.5 mg/mL phosphocreatine  
pH=7.35, 275 mOsm

## 19 External solutions: 1X ACSF

For 1L:

-  900 mL MilliQH2O +  100 mL stock 10X ACSF +  1.8 g NaHCO<sub>3</sub> +  2 g D-Glucose

10X ACSF stock

A	B	C
Drug	[mM]	10X Stock (g/4L)
NaCl	126	294.52
KCl	2.5	7.44
MgCl <sub>2</sub> *6H <sub>2</sub> O	1.2	9.75
NaH <sub>2</sub> PO <sub>4</sub> *H <sub>2</sub> O	1.2	6.64
CaCl <sub>2</sub> *2H <sub>2</sub> O	2.5	14.7

A	B	C
NaHCO <sub>3</sub>	21.4	
D-Glucose	11.1	