


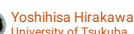

electroporation_{v.2}


Rodrigo Akiba, Liz Osoy, Nick Irwin, Patrick S. Keeling, Yoshimasa Hirakawa¹

¹University of Tsukuba, ²University of British Columbia (CA), ³University of British Columbia

1 Works for me [dx.doi.org/10.17504/protocols.io.35hgg36](https://doi.org/10.17504/protocols.io.35hgg36)

Plasmid information

Plasmid information


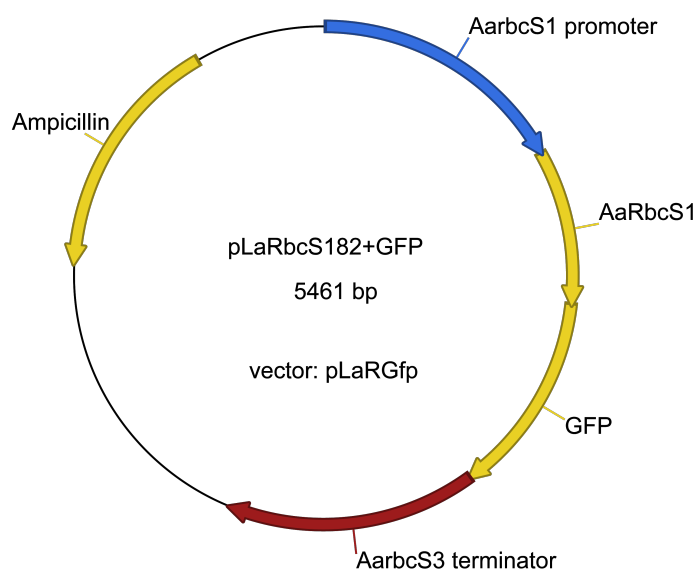


>pLaRbcS182+GFP
GAAGCTCGAGCAGCTGCACACATTAACATTCATGTGGATAATAAATATAATATATATATATATATATTATATATATATATATACAAAGCCCAAAATATTTATATATGTGTTTTAATGGAATCATGGGTCACATACATACAATTTGGTACGTGTTTACAATTCCTTACACCCCCCAAAACA

Hirakawa Y. Ishida K. (2010) Internal plastid-targeting signal found in a RubisCO small subunit protein of a chlorarachniophyte alga. The Plant Journal. 64: 402-410.

Plasmid preparation

- 1 Propagate plasmid DNA (pLaRGfp or its derivatives, such as pLaR182+GFP) in the *Escherichia coli* strain DH5 α .

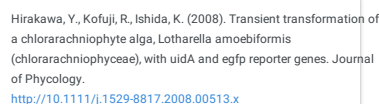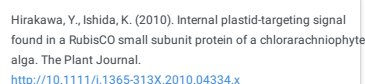

- 2 Purify plasmid DNA from 200 mL culture of *E. coli* by a Qiagen Plasmid Maxi Kit (Qiagen).
- 3 Adjust plasmid DNA concentration to 3-5 µg/µL with distilled water.

Cell culture

- 4 Culture *Amorhophclora amoebiformis* (CCMP2058) cells in 500 mL Erlenmeyer flasks containing 200 mL ESM medium at 20°C under white illumination (50-80 $\mu\text{mol photons}\cdot\text{m}^{-2}\cdot\text{s}^{-1}$) on a 14:10 hours light:dark cycle for a week.


Citation: Kodai Fukuba, Liz Cooney, Nick Irwin, Patrick J. Keeling, Yoshihisa Hirakawa (04/06/2020). Transformation of the chlorarachniophyte *Amorphochloris amoebiformis* by electroporation. <https://dx.doi.org/10.17504/protocols.io.35hgg36>

This is an open access protocol distributed under the terms of the **Creative Commons Attribution License** (<https://creativecommons.org/licenses/by/4.0/>), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

- 5 After decanting medium, resuspend cells adhered to the bottom of flasks by gentle pipetting with 2 to 3 mL ESM medium (use a glass pipette). Approximately 4×10^7 cells can be obtained from a flask.

Electroporation

- 6 Harvest a total of 5×10^6 cells by centrifugation at 2,000 g for 5 sec (use a mini centrifuge).
- 7 Resuspend cell pellet in 100 μ L of Gene Pulser electroporation buffer (Bio-Rad) with 10 μ L of plasmid DNA at the room temperature.
- 8 Transfer cell solution into electroporation cuvette with 0.2 cm gap (Bio-Rad).
- 9 Electroporate cells with a 25 ms square wave pulse at 120 V using Gene Pulser Xcell Electroporation System.
- 10 Add 0.9 mL fresh ESM medium to cuvette immediately after electroporation.
- 11 Transfer cells to glass bottom well plate/petri dish, and add an appropriate volume of ESM medium.
- 12 Incubate cells for 24 hours before observation of GFP fluorescence.

 This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](https://creativecommons.org/licenses/by/4.0/), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited