




DEC 18, 2023

## Isolation of Trophoblast Cells from Placenta

vfarm<sup>1</sup>

<sup>1</sup>Duke University

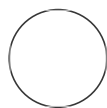

vfarm

### ABSTRACT

Protocol to derive trophoblast stem cells from placenta and culture them to form trophoblast organoids

OPEN  ACCESS



**Protocol Citation:** vfarm  
2023. Isolation of Trophoblast  
Cells from Placenta.

**protocols.io**

<https://protocols.io/view/isolation-of-trophoblast-cells-from-placenta-c6eyzbfw>

**License:** This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](#), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited


**Protocol status:** Working  
We use this protocol and it's working

**Created:** Dec 14, 2023

**Last Modified:** Dec 18, 2023



**PROTOCOL integer ID:**  
92344


## Preparations

1  0 °C the following items (and preferably have extras on hand if needed):

- 2 100mL bottles
- 2 small stir bars
- 1 large funnel
- cheese cloth cut into squares that are large enough to sit in the funnel
- dissection scissors

2 Equipment and Reagents Needed:



- stir plate that heats. alternatively could use a shaking water bath
- scalpels for dissection
- Matrigel, thawed
- Red Blood Cell Lysis Buffer (invitrogen, 00-4333-57)
- blunt P200 tips pre-chilled  4 °C
- 24-well plate warmed to  37 °C

3 Prepare the following Solutions and place all at  37 °C

The components needed to make all of the solutions below are:


- F12-K
- PennStrep
- 2.5% Trypsin
- Collagenase
- FBS

### 3.1 Wash Solution: F12K + PennStrep

 500 mL bottle of F12-K +  5 mL PennStrep

- have P/S aliquoted in 5 ml aliquots and only thaw when needed
- keep P/S aliquots in -20C
- I typically use ~400 mL of Wash Solution. No reason to need to make fresh each time, can use leftover as long as it is only a few months old


### 3.2 Trypsin Solution: 0.2% Trypsin

 75 mL of 0.2% Trypsin needed

Note: volume can be down to 25 mL, with the volume depending on the amount of starting tissue. I like to consistently use 75 mL.


 6 mL 2.5% trypsin +  69 mL PBS =  75 mL 0.2% trypsin

### 3.3 **Collagenase Solution:** 1 mg/ml collagenase in Wash Solution

 25 mL of collagenase

Dissolve  25 mg in  25 mL F12-K + P/S (Wash Solution)

### 3.4 **Stop Solution:** Final concentration of 10% FBS to inhibit Trypsin

Thaw FBS and use as needed. When using 75 mL of trypsin, I typically add  10 mL of 100% FBS to yield a final of just over 10% FBS

### 3.5 **TOM:** will need for the final plating


I typically plate 6 wells, therefore  3 mL of TOM is needed.

Will need to supplement TOM with Y-27632 so the final concentration is 10  $\mu$ M (add 1  $\mu$ l of Y-27632 for every 500  $\mu$ l of TOM).

## Isolation Protocol


39m

- 4 Cut out a lobule from the placenta. Carefully cut off the decidua to reveal chorionic villi, see image below:

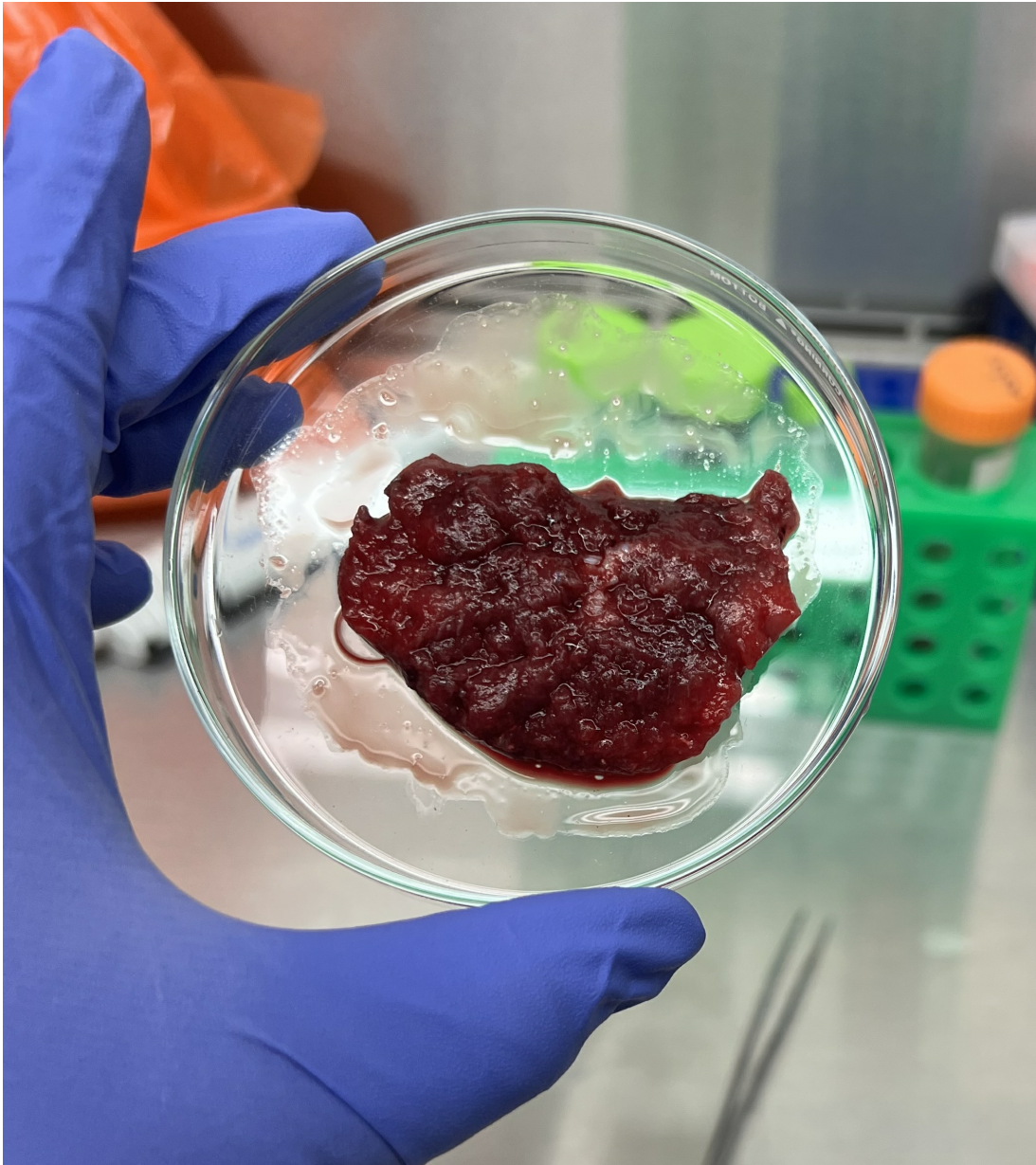


Decidua removed, villi facing up

- 5 Using a scalpel, generate very small fragments of tissue by scraping tissue off of the villus trees. The villus trees will be white veins running throughout the tissue. The smaller the tissue fragments, the more successful the isolation. Tissue should look like the tissue on the right side of the image below:

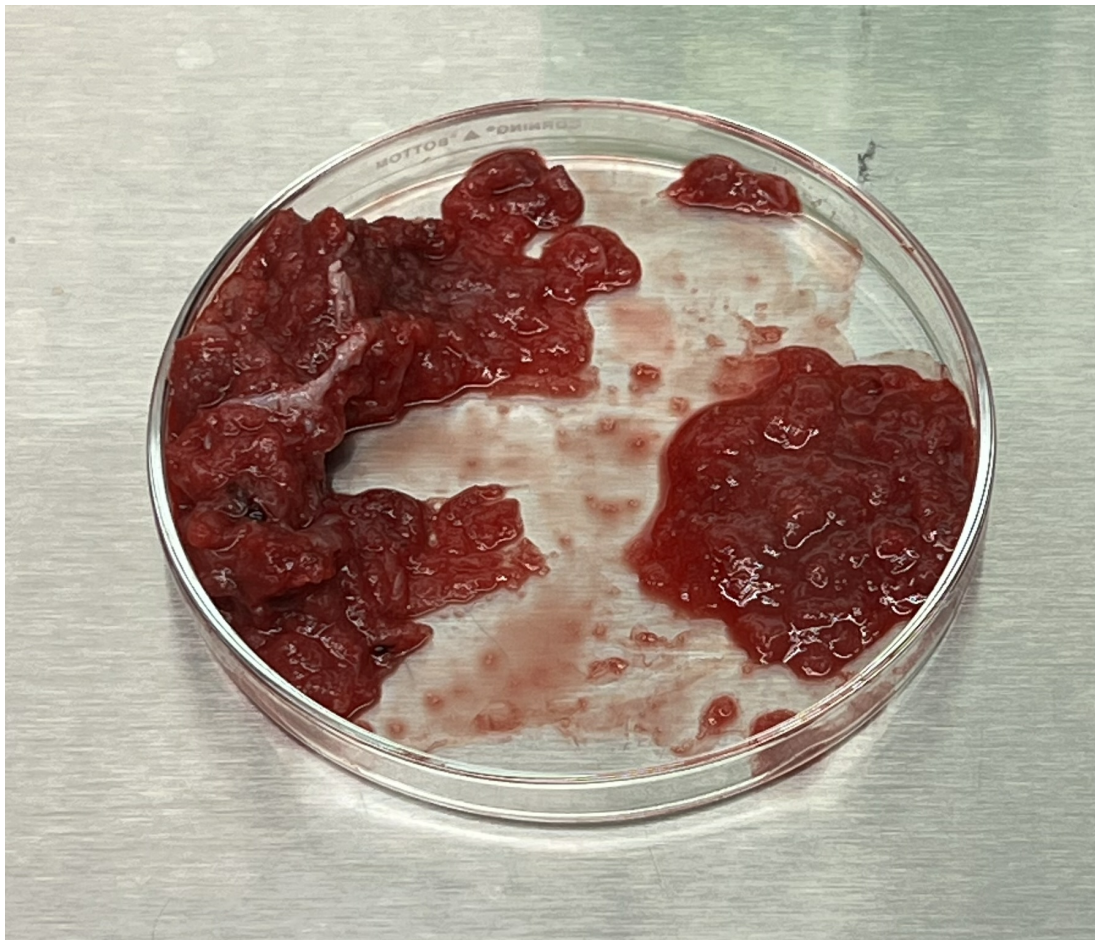

Tissue all cut up. Tissue on left will not be used, tissue on right will be used for the rest of the protocol.

**6** Place the tissue fragments into a 100 ml glass bottle and wash extensively with **Wash Buffer**.





**6.1** Add  50 mL of **Wash Buffer**, put cap on bottle and swirl bottle around ~5 times to wash tissue. Place bottle down and let tissue settle to the bottom of the bottle. Remove wash buffer using vacuum.

Note: Some tissue can get aspirated up with the tissue, depending on how well the tissue settles to the bottom of bottle. I have found that not introducing air bubbles while swirling helps increase the chance of the tissue sedimenting to the bottom of the bottle. As the wash steps go on, sedimentation gets worse and worse. Could think about centrifuging, but it hasn't been an issue yet for me.

**6.2** Repeat this at least 5 times! Wash Buffer will turn less dark pink/red throughout washes.

**6.3** When the Wash Buffer does not change from its light pink color, tissue is thoroughly washed.



- 7 Add  75 mL of Trypsin Solution and two small stir bars to the glass bottle containing the washed  8m tissue. Place bottle on  40 °C stir plate and incubate for  00:08:00 with gentle agitation.




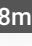
Note: I use 125 rpm on my stir plate. I use two small stir bars because I have found them bumping into each other helps the tissue disperse and not just stay in one clump.

- 7.1 Add FBS to a final concentration of about 10% to inactivate Trypsin. Add 8 mL of FBS.

Set bottle aside, will filter Collagenase/Tissue solution into this bottle later.

- 8 Place funnel on second glass bottle and place cheese cloth in funnel. Filter Trypsin/Tissue solution through cheese cloth.


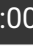
Note: Don't squeeze things through the cheese cloth; I have found an increase in the amount of debris when I do this









- 9 Gently scrape off tissue on cheese cloth and place it back in the same glass bottle with the stir bars.  25 mL of **Collagenase Solution** and incubate for  00:08:00 with gentle agitation on a  40 °C stir plate.  8m





- 10 Place funnel in the glass bottle containing the filtered cells and place new cheese cloth in funnel.

- 10.1 Remove bottle from stir plate and manually disrupt the tissue using a narrow orifice 10 mL serological pipette. Forcefully pipette up and down ~10 times to break up the tissue prior to filtering through cheese cloth.

- 10.2 Filter Collagenase/Tissue solution through cheese cloth into glass bottle containing the previous Trypsin/Tissue filtrate.

- 11 Centrifuge filtrate at  300 rcf, Room temperature,  00:04:00 . I typically have 100 ml of filtrate and split between two 50 ml conicals.

- 12 Remove supernatant using the vacuum (carefully, and don't touch pellet!). Resuspend pellets in  10 mL of F12-K and transfer to a 15 ml conical. I will resuspend both pellets in 5 mL, then wash each conical with the remaining 5 mL.
- 13 Centrifuge at  300 rcf, Room temperature, 00:04:00 4m
- 14 Carefully remove supernatant using serological pipette. Be careful not to disturb pellet!
- 15 Resuspend pellet in  1 mL of Red Blood Cell Lysis Buffer. Incubate at  Room temperature for 5m  00:05:00 .
- Note: This step is not essential but it does significantly help remove unwanted cells
- 16 Centrifuge at  300 rcf, Room temperature, 00:04:00 4m
- 17 Carefully remove supernatant using P1000, then use P200 to remove any excess media. Place conical in ice bucket.
- 18 Resuspend the pellet with  240  $\mu$ L of pre-thawed Matrigel using a pre-cooled blunt 200  $\mu$ L pipette tip (Fisher 02-707-134).
- 19 Carefully dispense  40  $\mu$ L of Matrigel-organoid suspension into a pre-warmed 24-well plate using a cold pipette tip.
- Notes: I usually set the pipette to 35  $\mu$ L (there will be lose of some Matrigel/cell mixture). Slowly and carefully lift up the pipette as dispensing the Matrigel into the well to form a dome. Do not push pipette tip fully down as this will introduce air-bubbles.

- 20** Place the 24-well plate in  37 °C incubator for  00:02:00 to allow Matrigel to pre-polymerize **2m**
- 21** Flip the plate over and incubate for an additional  00:08:00 to fully polymerize and evenly distribute the organoid fragments throughout the Matrigel **8m**
- 22** During the polarization process, prepare a stock of TOM with Y-27632 (final concentration 10 µM; add 1 µl of Y-27632 per 500 µl of TOM).
- Note: Need 500 µl of medium per well
- 23** Cover the polymerized Matrigel domes with 500 µl TOM per well and culture them in a  37 °C humidified 5% CO2 incubator.
- Note: Fill surrounding wells with 1 ml of PBS to help decrease evaporation of TOM
- 24** Keep an eye on the cells over time and replace media every 2-3 days. Keep extra Y-27632 (final of 10 µM) for the first 3 media changes.
- Notes: There will be so much debris, it will be impossible to see single cells. Debris and other cells will die overtime and clear organoids will begin to emerge after about 2-3 weeks.