

MAR 30, 2023

# OPEN ACCESS

#### יוסם

dx.doi.org/10.17504/protocol s.io.j8nlkwp56l5r/v1

**Protocol Citation:** michela.d eleidi, María José Pérez J., Hariam Raji 2023. Expansion microscopy. **protocols.io** https://dx.doi.org/10.17504/protocols.io.j8nlkwp56l5r/v1

License: This is an open access protocol distributed under the terms of the Creative Commons Attribution License, which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

**Protocol status:** Working We use this protocol and it's working

Created: Mar 28, 2023

Last Modified: Mar 30, 2023


**PROTOCOL** integer ID:

79525

Keywords: Expansion

microscopy

## Expansion microscopy

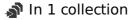

michela.deleidi<sup>1</sup>, María José Pérez J.<sup>1</sup>, Hariam Raji<sup>1</sup>

<sup>1</sup>German Center for Neurodegenerative Diseases (DZNE), Tübingen, 72076 Germany



Federico Bertoli

#### **ABSTRACT**

Expansion microscopy is a technique to visualize biological structures with higher spatial resolution than traditional microscopy methods.

**ATTACHMENTS** 

676-1424.docx

#### **Materials**

#### Stock X solution

A	В
Sodium acrylate 33% (w/v)	8.6% (w/v)
Acrylamide 50% (w/v)	2.5% (w/v)
N,N´-methylenebisacrylamide 2% (w/v)	0.15% (w/v)
5 M NaCl	11.7% (w/v)
PBS	1X

### **Digestion buffer**

A	В
Triton X-100	0.5% (w/v)
EDTA 0.5 M, pH 8	0.2% (v/v)
Tris-Cl 1 M, pH 8	5% (v/v)
NaCl	4.67% (w/v)
proteinase K	8 U/ml

- 10% (v/v) normal goat serum
- 0.1% (v/v) Triton X-100
- PBS
- secondary antibody (Alexa Fluor, Invitrogen)
- Acryloyl-X SE solution (Thermo Scientific)
- 10% (w/v) TEMED
- 10% (w/v) APS stock solution
- poly-L-ornithine-coated coverslips
- acryloyl-X SE solution
- 0.5% (w/v) 4-hydroxy-TEMPO stock solutions
- Leica TCS SP8 confocal microscope (Leica, Germany)
- GraphPad Prism version 9.0.0 (RRID:SCR\_002798)

3h

### **Expansion microscopy**

1



Note

This protocol refers to the expansion microscopy (ExM) protocol described in Asano et al., 2018 with some modifications.

Block cells with 10% (v/v) normal goat serum (NGS) in 0.1% (v/v) Triton X-100 in PBS and incubate it with primary antibodies in blocking solution (S) Overnight

2 After a 3-h incubation with the corresponding secondary antibody (Alexa Fluor, Invitrogen), wash the samples and treat with 0.1 mg/ml Acryloyl-X SE solution (Thermo Scientific) in PBS for





- 3 The freshly prepared gelling solution consisted of Stock X solution (8.6% (w/v) sodium acrylate 33% (w/v), 2.5% (w/v) acrylamide 50% (w/v), 0.15% (w/v) N,N´-methylenebisacrylamide 2% (w/v), 11.7% (w/v) NaCl 5 M, and PBS 1X), water, 10% (w/v) TEMED and 10% (w/v) APS stock solution in a 47:1:1:1 ratio.
- 4 Perform gel digestion Overnight in digestion buffer (0.5% (w/v) Triton X-100, 0.2% (v/v) EDTA 0.5 M, pH 8, 5% (v/v) Tris-Cl 1 M, pH 8, 4.67% (w/v) NaCl and 8 U/ml proteinase K).






Add the gelling solution to each well and covered by a 15-mm coverslip to ensure the formation of a smooth, flat and thin gel.



- 7 Expand the gel in water for 6001:00:00 and mount in 4000 undetermined poly-L-ornithinecoated coverslips to immobilize the gel for picture acquisition.

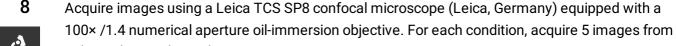

at least three independent experiments.

Analyze Images using Diffraction PSF 3D, DeconvolutionLab2, and EzColocalization plugins in Fiji-ImageJ.

- 10 Use GraphPad Prism version 9.0.0 (RRID:SCR\_002798) to calculate Spearman's rank correlation value  $(\rho)$  to identify colocalization of fluorescence signals.
- 11 The following is a variant of the protocol in case of using midbrain organoid sections:

Fix midbrain organoids and perform immunofluorescence staining as described above.

12



Treat sections with 🚨 0.1 undetermined acryloyl-X SE solution in PBS at Overnight Room temperature

- 13 Perform gelation in a 47:1:1:1 ratio of Stock X, 10% (w/v) TEMED, 10% (w/v) APS, and 0.5% (w/v) 4-hydroxy-TEMPO stock solutions.
- 14 Perform gel digestion and expansion as described above.
- 15

Acquire images using a Leica TCS SP8 confocal microscope (Leica, Germany) equipped with a 100× /1.4 numerical aperture oil-immersion objective. For each condition, acquire 5 images from at least three independent experiments.

Analyze Images using Diffraction PSF 3D, DeconvolutionLab2, and EzColocalization plugins in Fiji-ImageJ.



Use GraphPad Prism version 9.0.0 (RRID:SCR\_002798) to calculate Spearman's rank correlation value ( $\rho$ ) to identify colocalization of fluorescence signals.