

MAR 12, 2023

OPEN ACCESS

Protocol Citation: Noah Noah Snyder-Mackler 2023. Metagenomic Library Prep from fecal sample lysate . protocols.io

https://protocols.io/view/meta genomic-library-prep-fromfecal-sample-lysate-cnmnvc5e

License: This is an open access protocol distributed under the terms of the Creative Commons Attribution License, which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

Protocol status: In development

We are still developing and optimizing this protocol

Created: Jan 31, 2023

Last Modified: Mar 12, 2023

PROTOCOL integer ID: 76174

Keywords: Microbiome, Metagenome, Library, Lysate

Metagenomic Library Prep from fecal sample lysate

Noah Noah Snyder-Mackler¹

¹Arizona State University



smacklabasu

DISCLAIMER

dual index sequences here

Nextera 10bp dual index sequences

ABSTRACT

Metagenomic library prep from fecal sample lysate using Illumina DNA prep kit (1/2 reactions).

ATTACHMENTS

SMack_Lab_nxtra_index_p rimers.xlsx

GUIDELINES

Abbreviations

BLT: Bead-Linked Transposomes

TB1: Tagmentation Buffer 1

MM: Master Mix

TWB: Tagmentation Wash Buffer TSB: Tagmentation Stop Buffer

EPM: Enhanced PCR Mix

SPB: Sample Purification Beads

RSB: Resuspension Buffer

MATERIALS

- PCR plates and covers
- Nuclease-Free water
- Illumina® DNA Prep (M) Tagmentation (24 Samples IPB) Illumina, Inc. Catalog #20060060

Bead-Linked Transposomes (BLT)

Tagmentation Buffer 1 (TB1)

Tagmentation Wash Buffer (TWB)

Tagmentation Stop Buffer (TSB)

Enhanced PCR Mix (EPM)

Sample Purification Beads (SPB)

Resuspension Buffer (RSB)

- Freshly prepared 🔯 80% ethanol Fisher Scientific
- Index Adaptors (96)
- Magnetic stand

Preparation

1 Add \underline{A} 18 μ L of nuclease-free water and \underline{A} 2 μ L of sample to a PCR plate

Tagmentation Reaction

15m

- Take out TB1 and keep on On ice . Turn on thermocycler to
- Multiply by number of samples, for a 96 well plate use x100 and pipet mix together to create MM. Vortex BLT vigorously for 10 seconds to resuspend

Reagent	x1	x100
BLT	5 ul	500ul
TB1	5 ul	500ul

4 Vortex and add A 10 µL of reaction MM to each well

Note	
We found this easiest using a repeater with a $\ \ \ \ \ \ \ \ \ \ \ \ \ $	
Pipette mix slowly 5 times using multichannel (should have 20 µL of volume in each well	
Cover and seal plate. Then run "TAG" on thermocycler Choose the preheat lid option and set to Set the reaction volume to 30 µL Hold at 10°C	15m
During TAG incubation take out TSB and place at samples and vortex before addition to	
Add $\underline{\mathbb{Z}}$ 10 μL of TSB to each sample well using repeater with $\underline{\mathbb{Z}}$ 0.5 μL tip to stop enzymatic reaction	
Note	
Be careful because buffer is foamy and be sure to pipette mix after addition	

Over with thick plastic and run "PTC" protocol on thermocycler

■ Choose the preheat lid option and set to

15m

- Set the reaction volume to 🚨 40 µL
- Hold at 8 10 °C

5

6

7

8

10 Place on magnetic stand until beads pellet and discard supernatant

Note

Take out primer plates to thaw at 4 °C

Make PCR Master Mix

1h 10m 15s

- 12 Take out EPR and thaw on ice
- 1. Make master mix, multiply by number of samples (use x100 for a 96 well plate) (use a 5 ml tube for a full plate of MM)

Reagent	x1	x100
EPR	11 ul	1,056 ul
Nuclease free water	11 ul	1,056 ul

- 14 Place on magnetic stand until beads pellet + remove and discard supernatant (TWB)
- 15 Add \underline{A} 20 μL of MM using repeater and \underline{A} 5 mL attachment
- Add in $\underline{\underline{A}}$ 2.5 μ L of i5 and i7 primers

- Pipette mix, cover with thick plastic and run "BLT" 12 cycle on thermocycler (01:00:00) 1h 10m 15s

 Choose the preheat lid option and set to 100 °C
 - \$ 68 °C for \$ 00:03:00 • \$ 98 °C for \$ 00:03:00
 - 12 cycles of:
 - \$ 98 °C for \$ 00:00:45 \$ 62 °C for \$ 00:00:30 \$ 68 °C for \$ 00:02:00
 - \$\ 68 °C for \ \ \ 00:01:00
 - Hold at 8 10 °C

Note

- 1. Take out SPB and RSB while PCR is running to thaw
- 2. You may need to place the RSB on the thermocycler at 37 °C to get it to thaw


Bead Clean Up

10m

- Grab a new PCR plate and add \perp 60 μ L of nuclease-free water using a multichannel pipette and \perp 45 μ L of SPB using a repeater with a .5 ml tip to each sample well
- After PCR is complete spin down plate \rightarrow place on magnetic stand \rightarrow transfer \perp 25 μ L of supernatant to new PCR plate containing water and SPB. Discard old plate.
- Pipette mix 10 times and incubate for 00:05:00 at Room temperature

 During this incubation take out a new PCR plate and add Δ 15 μL of SPB to each well using a repeater
- 21 After the incubation place plate on magnetic stand until beads pellet

5m

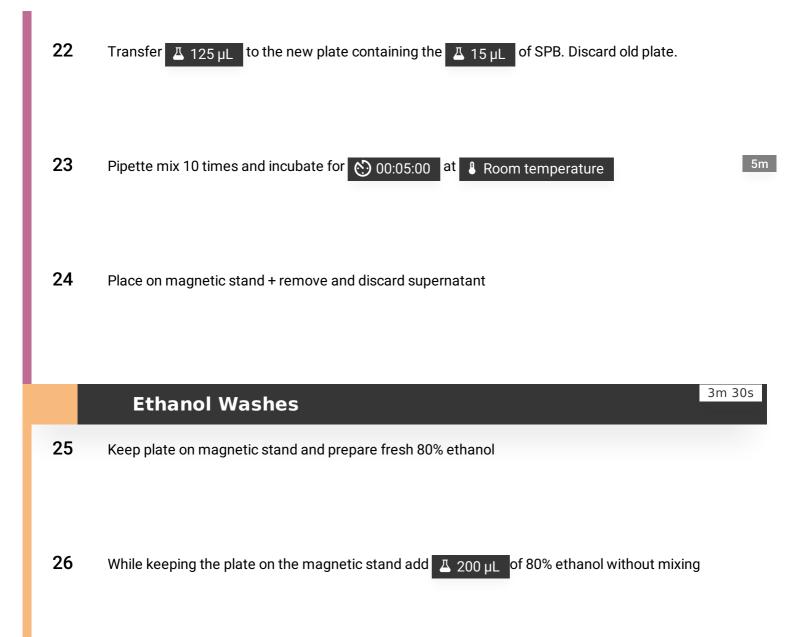

Incubate 00:00:30 + remove and discard supernatant. Repeat ethanol wash and discard

Use a 🔼 20 µL pipette to remove as much ethanol as possible

Air dry while on magnetic stand (© 00:03:00

27

28

supernatant.

Note

3m

30s

Note

be sure to not let beads crack

Elution 29 Remove beads from stand and add Δ 32 μL of RSB to beads 30 Pipette mix to resuspend and incubate at Room temperature for 00:02:00 2m 31 Place on magnetic stand and transfer supernatant to new and final sturdy PCR plate. Seal, label, and store!