



# Nuclei Extraction for tissue using Iodixanol Gradients V.1

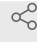
Kriegstein lab<sup>1</sup>

<sup>1</sup>UCSF


Version 1 ▾

Oct 24, 2022

1 *Works for me*

 Share

[dx.doi.org/10.17504/protocols.io.81wgby7yqvpk/v1](https://dx.doi.org/10.17504/protocols.io.81wgby7yqvpk/v1)

 Rachelle Leonard

## ABSTRACT

Nuclei isolation using Iodixanol gradients geared for multiome

## DOI

[dx.doi.org/10.17504/protocols.io.81wgby7yqvpk/v1](https://dx.doi.org/10.17504/protocols.io.81wgby7yqvpk/v1)

## PROTOCOL CITATION

Kriegstein lab 2022. Nuclei Extraction for tissue using Iodixanol Gradients.

**protocols.io**

<https://dx.doi.org/10.17504/protocols.io.81wgby7yqvpk/v1>



## LICENSE

————— This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](https://creativecommons.org/licenses/by/4.0/), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

## CREATED

Oct 06, 2022

## LAST MODIFIED

Oct 24, 2022

## PROTOCOL INTEGER ID

70981

## Creating Buffers

### 1 Stock Buffers

All stock solutions should be filtered using a .22 um PVDF/PES filter system. All solutions except 50% Iodixanol solution are stable at 4c for at least 6 months.

#### 1 M Sucrose (300 mL)

Substance:	Stock Conc.:	Amount:	Final conc. in working solution:
Sucrose	-	102.69 g	1 M
Water	-	235.5 mL	-

#### 1.0616x Homogenization Buffer Stable Solution (200 mL)

Substance:	Stock conc:	Amount:	Final conc. in working solution:
Sucrose	1 M	53.1 mL	.2653 M
KCl	2 M	2.66 mL	26.6 mM
MgCl <sub>2</sub>	1 M	1.06 mL	5.31 mM
Tricine-KOH pH 7.8	.75 M	5.67 mL	21.2 mM
Water	-	137.5 mL	-

#### Diluent Buffer (100 mL)

Substance:	Stock conc:	Amount:	Final conc. in working solution:
KCl	2 M	7.5 mL	150 mM
MgCl <sub>2</sub>	1 M	3 mL	30 mM
Tricine KOH pH 7.8	.75 M	16 mL	120 mM
Water	-	73.5 mL	-

#### 50% Iodixanol Solution (50 mL) \*\*Remake Monthly for Stability

Substance:	Stock conc:	Amount:	Final conc. in working solution:
Diluent Buffer	-	8.3 mL	-
Iodixanol	60%	41.7 mL	50%

#### 1x Homogenization Buffer Unstable Solution for 4 reactions **Prepare Fresh**

Substance:	Stock conc:	Amount:	Final conc. in working solution:
HB stable solution	1.0616X	7536 uL	1X
DTT	1 M	8 uL	1 mM
Spermidine	500 mM	8 uL	.5 mM
Spermine	150 mM	8 uL	.15 mM
NP40	10%	240 uL	.3%
cOmplete PI (diluted in HB stable)	100X	80 uL	1X
Ribolock	40U/uL	120 uL	.6 U/uL

#### 30% Iodixanol Solution per reaction **Prepare Fresh**

Substance:	Stock conc:	Amount:	Final conc. in working solution:
HB unstable	-	240 uL	-
50% Iodixanol Solution	50%	360 uL	30%

#### 40% Iodixanol Solution per reaction **Prepare Fresh**

Substance:	Stock conc:	Amount:	Final conc. in working solution:
HB unstable	-	120 uL	-
50% Iodixanol Solution	50%	480 uL	40%

#### Wash Buffer 1 mL for 4 reactions **Prepare Fresh**


Substance:	Stock conc:	Amount:	Final conc. in working solution:
Tris-HCl pH 7.4	1 M	10 uL	10 mM
NaCl	5 M	2 uL	10 mM
MgCl <sub>2</sub>	1 M	3 uL	3 mM
BSA	30%	33.3 uL	1%
Tween-20	10%	10 uL	.1%
DTT	1 M	1 uL	1 mM
Ribolock	40 U/uL	15 uL	.6 U/uL

## 2 Before starting protocol:

- 1) pre-chill swinging bucket centrifuge and a fixed angle centrifuge to 4c
- 2) Pre-chill douncers and pestles to 4c on ice
- 3) Pre-chill tubes
- 4) Fill up 2 L beaker with 500 mL sterile water to soak the used Douncers

## 3 Isolation of Nuclei via Dounce Homogenization and Density Centrifugation

- 1) place 20-50 mg frozen tissue or crushed into pre-chilled 7 mL dounce containing 1 mL cold 1x HB
- 2) Dounce with "A" loose pestle until resistance goes away (~10 strokes)
- 3) Dounce with "B" tight pestle until resistance goes away (~15 strokes)
- 4) Place "A" and "B" into sterile water to soak for cleaning later
- 5) Filter during transfer into FACS tube
- 6) Place Dounce into beaker with sterile water to soak for cleaning later
- 7) Pellet nuclei by spinning 5 min at 4c at 350 xg in a fixed angle centrifuge
- 8) Remove 950 uL of supernatant (50 uL remaining)
- 9) Gently resuspend nuclei in 350 uL 1x HB
- 10) Add 1 volume (400 uL) of 50% Iodixanol solution and pipette mix
- 11) Slowly layer 600 uL of 30% Iodixanol solution under the 25% mixture. Wipe side of pipette tip with kimwipe to avoid mixing layers.
- 12) Layer 600 uL 40% Iodixanol solution under the 30% mixture. Wipe side of pipette tip with kimwipe to avoid mixing layers.
- 13) In a pre-chilled swinging bucket centrifuge, spin for 20 min at 4c at 3000 xg with the brake off. Set acceleration level to 1 and deceleration at 0. (centrifuge time=23 min, time to stop=13 min)
- 14) Slowly extract top layers in increments of 200 uL down to 200-300 uL of nuclei band between 30% and 40% interface.
- 15) Take 200 uL of the nuclei band and put into 1.5 mL LoBind tube
- 16) Dilute nuclei by adding 200 uL of wash buffer and mix by pipetting
- 17) Count nuclei using trypan blue
- 18) Centrifuge at 500 xg for 5 min at 4c

- 
- 19) Remove supernatant without disrupting pellet
  - 20) resuspend in x uL of chilled Nuclei Buffer (depending on what isolated nuclei are used for) to achieve target concentration based on count.