



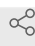
Oct 09, 2022

PCR


Ana Belem García González¹, Jair Alexis Gardea Sáenz¹,
Irán Alessandra Chaparro Rodríguez¹, Georgina Diego¹

¹Tecnologico de Monterrey Campus Chihuahua

1 Works for me

 Share

dx.doi.org/10.17504/protocols.io.bp2l694n1lqe/v1

 Elena Elizabeth Mercado Flores

ABSTRACT

50 µL volume PCR protocol

DOI

dx.doi.org/10.17504/protocols.io.bp2l694n1lqe/v1

PROTOCOL CITATION

Ana Belem García González, Jair Alexis Gardea Sáenz, Irán Alessandra Chaparro Rodríguez, Georgina Diego 2022. PCR. **protocols.io**
<https://dx.doi.org/10.17504/protocols.io.bp2l694n1lqe/v1>



LICENSE

————— This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](https://creativecommons.org/licenses/by/4.0/), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

CREATED

Sep 08, 2022

LAST MODIFIED

Oct 09, 2022

PROTOCOL INTEGER ID

69748

RECOVERY IDT SEQUENCES

- 1 Centrifuge the tubes with the sequences for one minute at 3,000 g x m in the centrifuge

- 2 Centrifuge for 5 seconds in the microcentrifuge
- 3 Add nuclease-free water needed to have a final concentration of 10 ng/uL.
- 4 Vortex, followed by incubating for 20 minutes at 50°C.
- 5 Mix briefly by vortex and centrifuge for 8 seconds in mini centrifuge.

PCR

- 6 Add nuclease-free water to the first forward and reverse until a concentration of 100 uM is reached, centrifuged quickly in microcentrifuge.
- 7 Make dilutions with 10 uL of primers and 90 uL of nuclease-free water to reach a concentration of 10 uM
- 8 In 0.2 mL Eppendorf tubes make the following mixture.
- 9

| A | B |
|-------------------------|----------|
| Nuclease Free Water | 17 uL |
| Buffer 10X PCR | 5 uL |
| 50 mM MgCl ₂ | 1.5 uL |
| 10 mM dNTP's | 1 uL |
| 10 uM Forward primer | 2.5 uL |
| 10 uM Reverse primer | 2.5 uL |
| DNA (25 mg/uL) | 20 uL |
| Taq Platinum Polymerase | 0.5 uL |
| Final Volume | 50 uL |

10 Perform the next cycle in thermocycler

| A | B | C | D |
|-------------|--------------------|------------|----------|
| 1 Cycle | Initial activation | 2 minutes | 94°C |
| 35 Cycles | Denaturation | 30 seconds | 94°C |
| 35 Cycles | Alignment | 30 seconds | 50°C |
| 35 Cycles | Extension | 2 minutes | 72°C |
| 1 Cycle | Final Extension | 5 minutes | 72°C |
| Maintenance | Maintenance | - | 4°C |