




JUN 22, 2023

## Purification and analysis of SKP1-FBXO7 complexes

Frank

Adolf<sup>1</sup>, Brenda A. Schulman<sup>1</sup>

<sup>1</sup>Max-Planck Institute of Biochemistry, Research Department Molecular Machines and Signaling

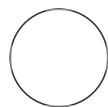

Felix Kraus

OPEN  ACCESS

### DOI:

[dx.doi.org/10.17504/protocols.io.6qpvr4z4ogmk/v1](https://dx.doi.org/10.17504/protocols.io.6qpvr4z4ogmk/v1)

**Protocol Citation:** Frank Adolf, Brenda A. Schulman 2023. Purification and analysis of SKP1-FBXO7 complexes. **protocols.io** <https://dx.doi.org/10.17504/protocols.io.6qpvr4z4ogmk/v1>

**License:** This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](https://creativecommons.org/licenses/by/4.0/), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

**Protocol status:** Working  
We use this protocol and it's working

**Created:** Mar 28, 2023

**Last Modified:** Jun 22, 2023

**PROTOCOL integer ID:**  
79610

**Keywords:** ASAPCRN

### ABSTRACT

Protocol for the biochemical purification and analysis of SKP1-FBXO7 complexes

### GUIDELINES

Please wear appropriate PE while performing the experiment.

Please familiarise yourself with the laboratory safety rules and guidelines and follow these while performing the experiment.

### MATERIALS

Buffer A for size exclusion chromatography (SEC):

- 25mM HEPES pH7.5 (KOH)
- 150 mM NaCl
- 1mM DTT

TB medium

**Preparatory note**

- 1 All constructs were prepared utilizing standard molecular biological techniques and verified by sanger sequencing.

The cDNAs coding for HsFbxo7-129-398 and Skp1 preceding a second RBS were cloned into pGEX4T1 as previously described for other substrate receptor/Skp1 complexes (Schulman et al. 2000) (pGEX4T1-TEV-HsFbxo7-129-398/HsSkp1).

The cDNA coding for HsPI31-1-151 (with N-terminal TEV cleavable His8-tag) was cloned into pRSF1b (pRSF1b-His8-TEV-HsPI31-1-151).

For co-expression of Fbxo7/Skp1/PI31 complexes, both plasmids were co-transformed into *E. coli* BL21 Rosetta (DE3).

## Molecular biological methods and protein expression

- 2 Grow *E. coli* cultures in Terrific Broth (TB) medium at 37°C.  
At OD600 of 0.8, induce expression with 0.5mM IPTG.  
Continue growing the *E. coli* culture for 16h at 18°C. 1d
- 3 Purify FBX07/Skp1/PI31 complexes by sequential standard GST- and His-affinity chromatography.  
Cleave affinity tags by incubation with TEV protease at 4°C for 16h. 1d
- 4 Further purify complexes by preparative size exclusion chromatography (SEC) in buffer A (25mM HEPES pH7.5 (KOH), 150 mM NaCl, 1mM DTT; for details please see next section) on a Superdex 200 Increase 10/300 GL column.
- 5 Pool fractions of interest, aliquoted and snap freeze in liquid N<sub>2</sub>.  
Store fractions at -80°C until further usage.
- 6 HsCul1-1-410 was expressed as GST-fusion protein and purified as described previously (Hopf et al., 2022).

- 7 Analytical SEC was carried out on an ÄKTApure system (GE Healthcare) equipped with a Superdex 200 Increase 10/300 GL column (Cytiva), in buffer A (25mM HEPES pH7.5 (KOH), 150 mM NaCl, 1mM DTT).
- 8 Preincubate samples at 37°C for 10 min before loading.  
Samples:  
HsCul1-1-410,  
HsFbxo7-129-398/HsSkp1/HsPI31-1-151  
HsFbxo7-129-398/HsSkp1/HsPI31-1-151 + HsCul1-410
- 9 Apply samples (100 µl at a concentration of 45 µM) on column.  
Set flow rate to 1 ml/min.
- 10 Record UV absorbance at 280 nm and collect fractions of 200 µl volume.
- 11 Analyze fractions by SDS-PAGE.

30m